

Periotomes

Periotomes are used in a procedure that facilitates the removal of teeth with minimal damage to the surrounding alveolar bone. This procedure may be used whenever an extraction is indicated, but is particularly important when considering the placement of a dental implant.

Exclusive Immunity Steel™ for the sharpest and longest lasting cutting edges

Thin, clinically-designed blades sever periodontal ligament easily

Blades are permanently fixed to handle and will not rotate under pressure

Satin Steel™ handle features maximum diameter for increased comfort and decreased hand fatigue

PT1X

PT5X

PT6X


PT3X

Periotomes

Atraumatic Extractions



Instrumentation and
Procedural Products

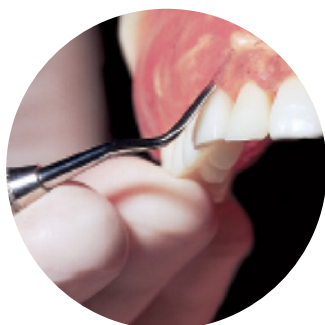

The periotome is used to sever the periodontal ligament (PDL) from the root surface of the tooth. The blade must be angled approximately 20 degrees out from the root to ensure the tip of the periotome is within the crest of the alveolar bone and to prevent the blade from slipping.




As shown here, the periotome is pushed several millimeters into the PDL space and is repeated circumferentially until reaching approximately 3/4 of the distance to the apex.

At this point, the tooth is attached only by the most apical portion of the PDL. The tooth may then be extracted without thrusting the extraction forceps deep into the PDL space, preventing damage to the alveolar bone.



The periotome blade is pushed into the depth of the sulcus to sever the PDL circumferentially. Repeat this step in stages, to ensure all attachments have been severed.

Most Popular Black Line Periotomes

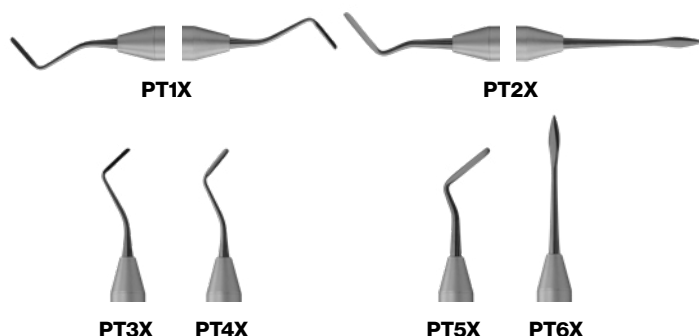

SKU	Description
PT1X	Posterior Periotome, Double-Ended, Black Line, 1.5mm/1.5mm
PT2X	Anterior Periotome, Double-Ended, Black Line, 1.9mm/3.2mm
PT3X	Posterior Periotome, Right, Single-Ended, Black Line, 1.5mm
PT4X	Posterior Periotome, Left, Single-Ended, Black Line, 1.5mm
PT5X	Anterior Periotome, Angled, Single-Ended, Black Line, 1.5mm
PT6X	Anterior Periotome, Straight, Single-Ended, Black Line, 3.2mm



Scan for more information about
Atraumatic Extractions

For more information, visit **HuFriedyGroup.com**