

HU-FRIEDY STEP PLIERS

Hu-Friedy Step Pliers allow clinicians to easily create z-bends (bayonet bends) in one, simple step by using a gentle squeeze of the pliers handle.


Hu-Friedy Pliers are crafted with Immunity Steel® alloy, a high quality surgical grade stainless steel which combines an optimal blend of carbon and chromium ensuring strength and corrosion resistance.

DOUBLE-SIDED OFFSET
BEAKS PERMIT BENDS ON
WIRES UP TO .022" X .025"
(.53 MM X .64 MM)

Step Pliers, 1/2 mm | **678-340**

Step Pliers, 3/4 mm | **678-341**

Step Pliers, 1 mm | **678-342**

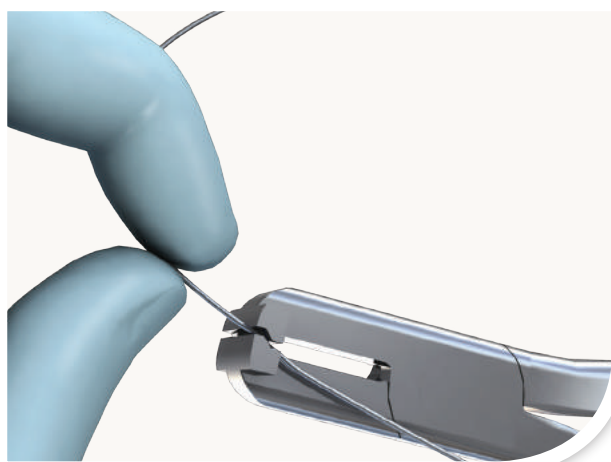

Step Pliers help to create a right-hand or left-hand step-up or step-down in wire at 1/2 mm, 3/4 mm or 1 mm intervals.



Most commonly, z-bends are utilized in cases where there is a need to extrude or intrude a specific tooth, or teeth. The bends can also change inclination and/or rotation, advancing a patient's orthodontic treatment.

Visit us online at **Hu-Friedy.com**

Call 1-800-Hu-Friedy or contact your authorized Hu-Friedy representative for more information.
©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFO-415/0515

How the best perform

Hu-Friedy