

Sterilization Container System

Provides Efficient Sterilization of Cassettes and Instruments



Infection Prevention and Instrument Management

HuFriedyGroup Sterilization Containers are a reusable, rigid system used for the packaging of cassettes and instruments during sterilization. An alternative to the wrap and pouch materials commonly used to package cassettes and instruments, Sterilization Containers provide significant economic and environmental benefits along with enhanced safety.

- Allows sterilization of enclosed contents and maintains sterility during transport and storage, until opened.
- Provides fully enclosed and locked transit of loose instruments and cassettes.
- Offers economic and environmental benefits by providing reusable packaging, low cost of use, a long service life, and minimal use of disposable materials.
- Uniquely sized to fit many Hu-Friedy cassettes, including the Container Cassette that features an integrated handle for easy insertion and removal.
- Includes a full line of Container accessories including filters, indicator cards, and tamper evident seals.




Points of performance

Two external display windows provide space for a process indicator card and container identification

Tamper evident seals on latches provide visual evidence of sterility breach

Latches on the lid act as handles for easy, contamination free lid removal

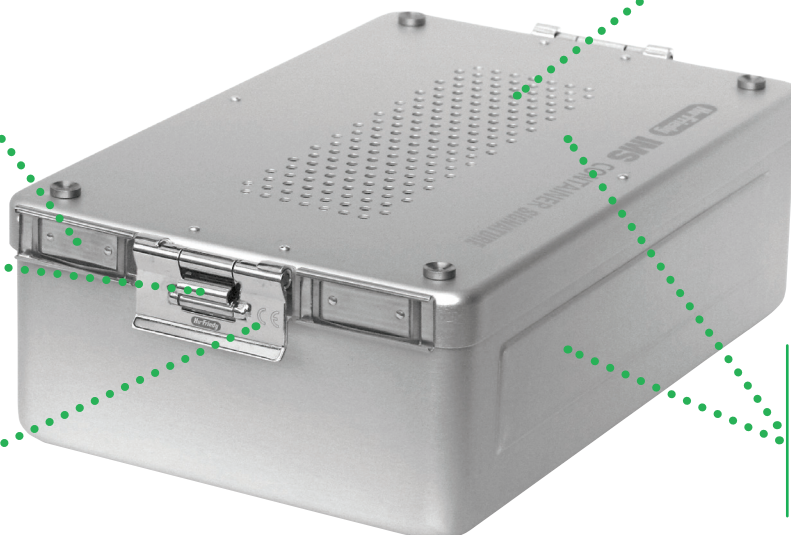

The filter retention plates audibly click to lock into place providing easy filter changes

Lid and bottom are made of anodized aluminum making them lightweight with excellent drying properties. The aluminum alloy with an anodized oxide surface helps resist corrosion

Sterilization Container System

Provides Efficient Sterilization of Cassettes and Instruments



Infection Prevention and Instrument Management



Container Cassette for 12 Instruments

IMECC128

Blue
268 x 176 x 34 mm
10-1/2 x 7 x 1-1/3"

Container Cassette with handle on lid and filter outlets.
Fits 12 instruments, with an accessory area.



IMDINCO0M

312 x 190 x 50 mm
2-1/4 x 7-1/2 x 2"

IMDINCO1M

312 x 190 x 65 mm
12-1/4 x 7-1/2 x 2-1/2"

Top and bottom with filter.
Fits one container cassette.



IMDINCO2M

312 x 190 x 100 mm
12-1/4 x 7-1/2 x 4"

Top and bottom with filter.
Fits two container cassettes.



IMDINCO3M

312 x 190 x 135 mm
12-1/4 x 7-1/2 x 5-1/3"

Top and bottom with filter.
Fits three container cassettes.

CONTAINER ACCESSORIES

IMCO-INDNA

Container Indicator Card
250 pieces, single use

IMCOM-SEAL

Tamper Evident Seal for Latch
100 pieces, single use

IMCO-FIL3M

Disposable Paper Filter
100 pieces, blank, single use



Tamper Evident Seal for Latch

Container Indicator Card

Visit us online at HuFriedyGroup.com/SterilizationContainers

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.
©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved. IMS-372/2605

 **HuFriedyGroup**
The Best In Practice