



SURGICAL **SOLUTIONS**

PROCEDURAL KIT LISTS

How the best perform

Hu-Friedy



HU-FRIEDY SURGICAL KIT

ATRAUMATIC EXTRACTION KIT



PART CODE	DESCRIPTION
EPTSMS	Luxating Hybrid Elevator, Small Straight
EPTSMC	Luxating Hybrid, Small Curved
EPTSSP	Luxating Hybrid Elevator, Straight Spade
FAF150	Apical Forceps 150
FAF151	Apical Forceps 151
CM2/4S6	2/4 Molt DE Curette, #6 Handle
NHCW	Crile-Wood Needle Holder
S16SC	16 G-F Scissors Super Cut
10-130-05	5 Scalpel Handle, Round Straight
P96	9 Molt Periosteal Elevator
IM1002	6-Pack Signature Series® Instrument Clip
IM4161	16 Instrument Signature Series® Cassette, Gray
IMATRAEXT	Complete Atraumatic Extraction Kit

ATRAUMATIC EXTRACTION KIT

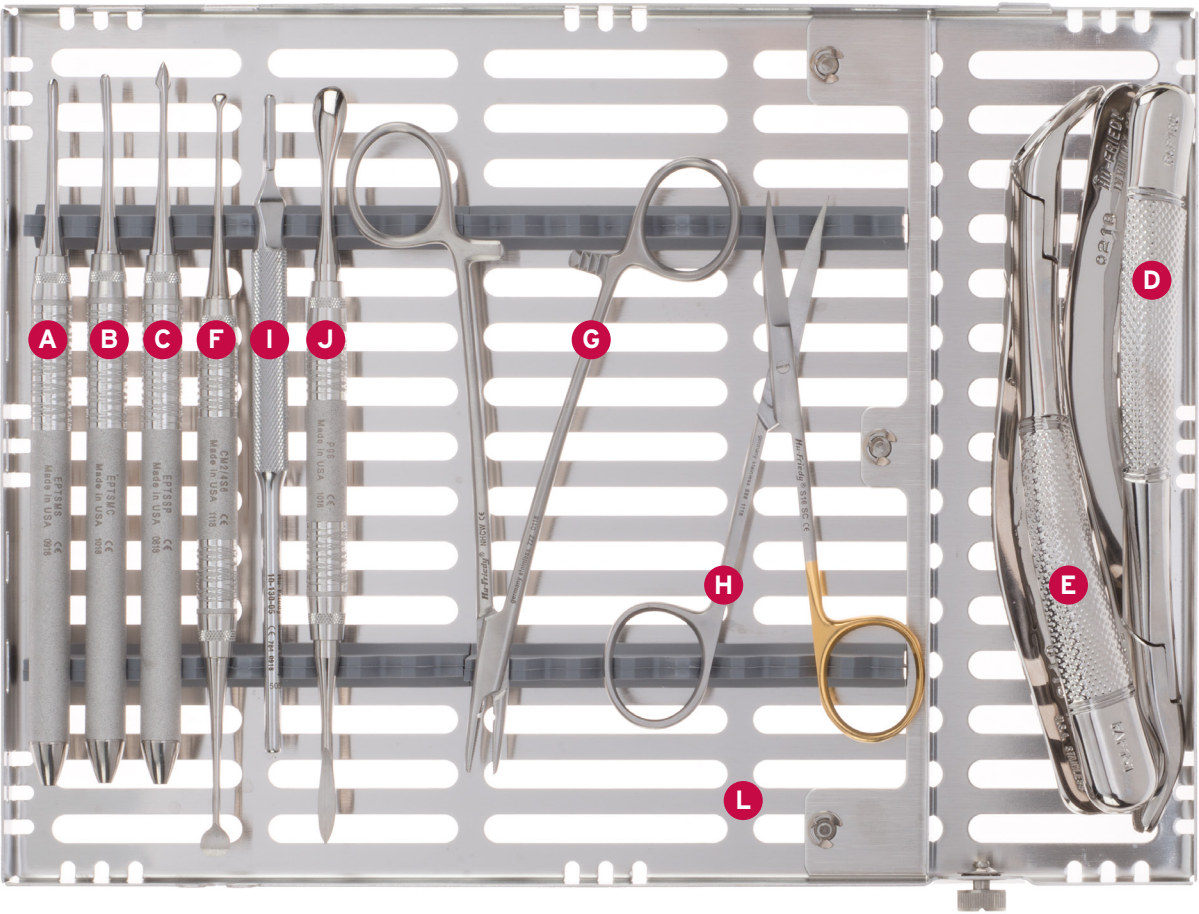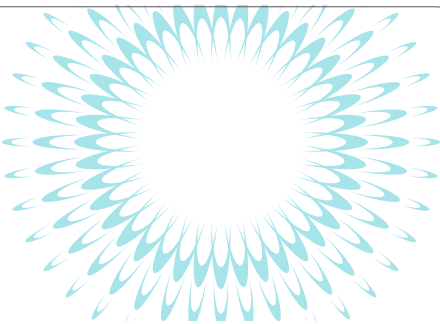

ATRAUMATIC EXTRACTION KIT



- A EPTSMS**
Luxating Hybrid Elevator, Small Straight
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- B EPTSMC**
Luxating Hybrid, Small Curved
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSSP**
Luxating Hybrid Elevator, Straight Spade
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D FAF150**
Apical Forceps 150
 Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- E FAF151**
Apical Forceps 151
 Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- F CM2/4S6**
2/4 Molt DE Curette, #6 Handle
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- G NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- H S16SC**
16 Goldman-Fox Scissors, Super Cut
 Used for trimming tissue or cutting sutures. Featuring tapered blade tips, tips that are not sharp on outside and are used as tonotomy scissors. 12.5cm with 1 serrated blade



- I 10-130-05**
5 Scalpel Handle, Round Straight
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- J P96**
9 Molt Periosteal Elevator
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm

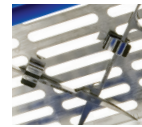

- K IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- L IM4161**
16 Instrument Signature Series® Cassette, Gray
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.


HU-FRIEDY SURGICAL KIT

ATRAUMATIC EXTRACTION KIT, BLACK LINE




PART CODE	DESCRIPTION
EPTSMSX	Luxating Hybrid Elevator, Small Straight, Black Line
EPTSMCX	Luxating Hybrid, Small Curved, Black Line
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
FAF35XS	Extracting Forceps, 35N
FAF36XS	Extracting Forceps, 36N
CM2/4X	2/4 Molt DE Curette, Black Line
NHCW	Crile-Wood Needle Holder
S16SCX	Goldman-Fox, Curv, SC Scissors Black Line
10-130-05E	Scalpel Handle #5
P9X	9 Molt Periosteal Elevator, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMN4200	20 Instrument Infinity Series™ Cassette, Red
IMNATREXTX	Atraumatic Extraction Kit, Black Line

ATRAUMATIC EXTRACTION KIT, BLACK LINE

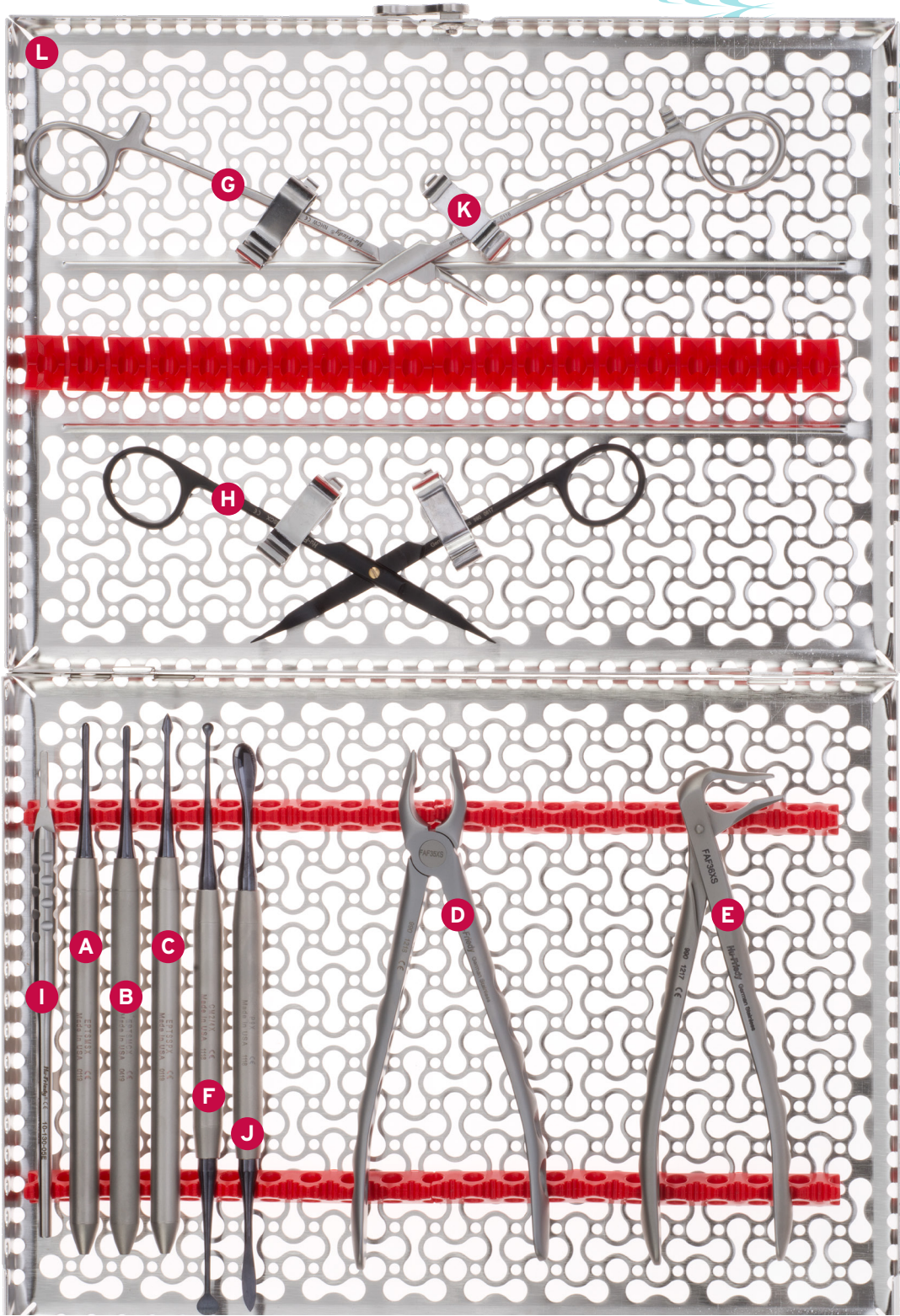

ATRAUMATIC EXTRACTION KIT, BLACK LINE

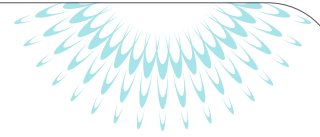

- A EPTSMSX**
Luxating Hybrid Elevator, Small Straight, Black Line
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- B EPTSMCX**
Luxating Hybrid, Small Curved, Black Line
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSSPX**
Luxating Hybrid Elevator, Straight Spade, Black Line
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D FAF35X**
Extracting Forceps, 35N
 Atraumair Forceps used to extract upper canines & premolars.



- E FAF36X**
Extracting Forceps, 36N
 Atraumair Forceps used to extract lower incisors, canines and premolars.



- F CM2/4X**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- G NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- H S16SCX**
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- I 10-130-05E**
5 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.





- J P9X**
9 Molt Periosteal Elevator, Black Line
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm

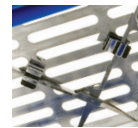

- K IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- L IMN4200**
20 Instrument Infinity Series™ Cassette, Red
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.

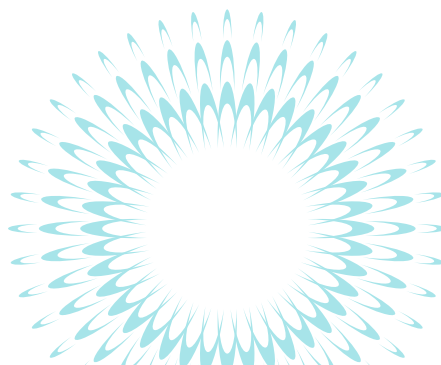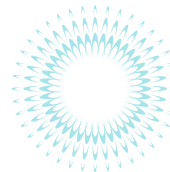



HU-FRIEDY SURGICAL KIT

BASIC PERIO KIT



PART CODE	DESCRIPTION
IM4121	12 Instrument Signature Series® Cassette, Gray
IM1002	6-Pack Signature Series® Instrument Clip
MIR5	5 FS CS Mouth Mirror
MH6	6 Cone Socket Mirror Handle
10-130-03	3 Scalpel Handle
KO12KPO3AR	1/2 Modified Allen Orban Knife Round
XP23/QOW6	23/QOW Expro, Wms Mkgs, #6 Handle
SG11/126	11/12 Gracey Curette, #6 Handle
SG13/146	13/14 Gracey Curette, #6 Handle
PPR36	Prichard Periosteal, #6 Handle
CM2/4S6	2/4 Molt DE Curette, #6 Handle
S13K/TG6	13/TG DE Curette/Chisel, #6 Handle
NHCW	Crile-Wood Needle Holder
S16	16 Goldman-Fox Scissors
IMBASPER	Complete Basic Perio Kit

BASIC PERIO KIT

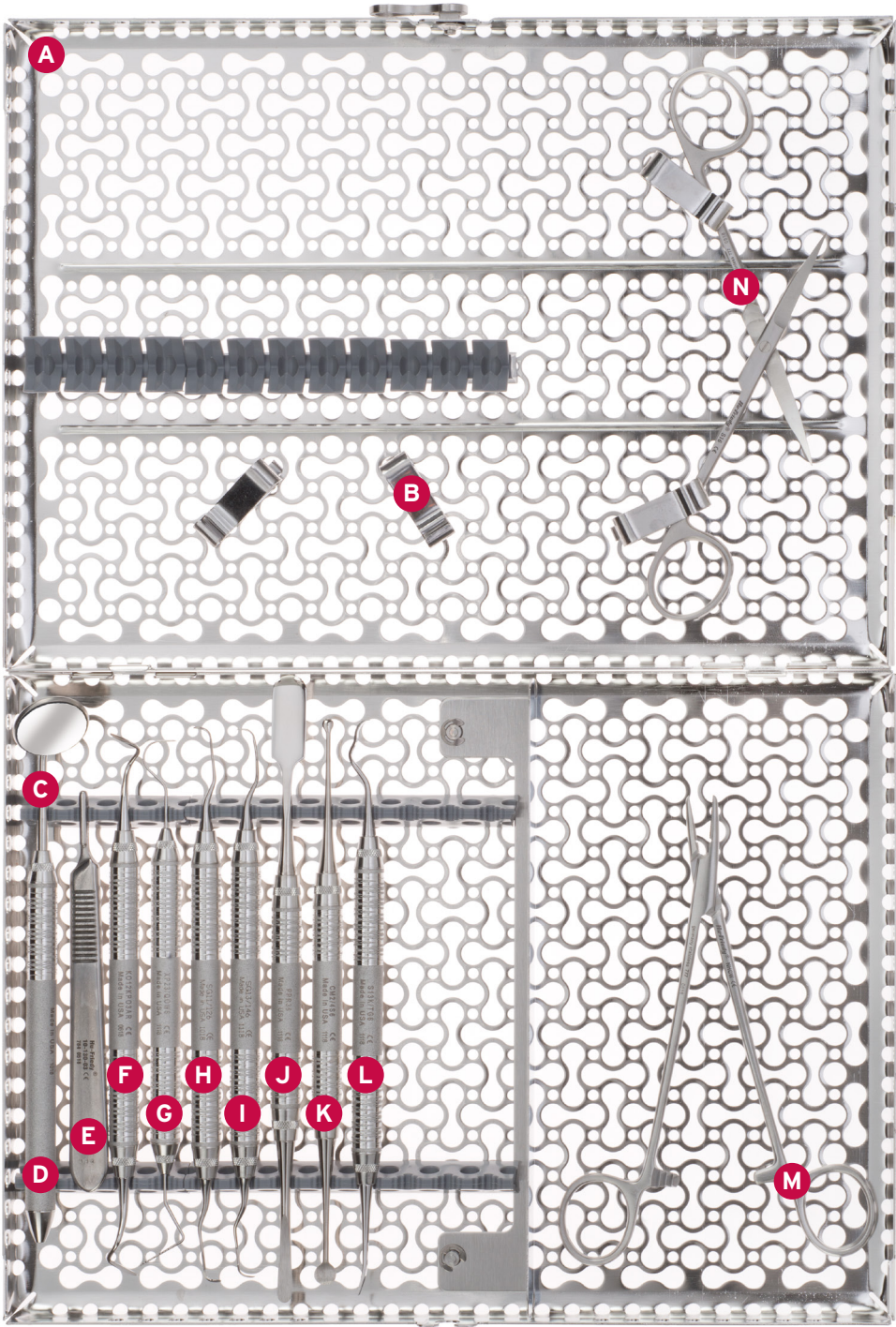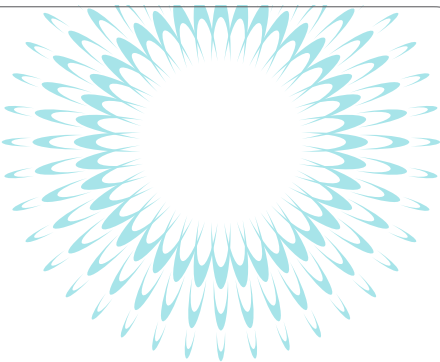

* pictured in IMN8121 - 12 Instrument Infinity Series™ Cassette, Gray

BASIC PERIO KIT



A IM4121
12 Instrument Signature Series® Cassette, Gray

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



C MIR5
5 FS CS Mouth Mirror

Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



D MH6
6 Cone Socket Mirror Handle

Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



E 10-130-03
3 Scalpel Handle

Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



F KO12KPO3AR
1/2 Modified Allen Orban Knife Round

30% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



G XP23/QOW6
23/QOW Expro, Wms Mkgs, #6 Handle

A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



H SG11/126
11/12 Gracey Curette, #6 Handle

Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



I SG13/146
13/14 Gracey Curette, #6 Handle

Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



J PPR36
Prichard Periosteal, #6 Handle

For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



K CM2/4S6
2/4 Molt DE Curette, #6 Handle

For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



L S13K/TG6
13/TG DE Curette/Chisel, #6 Handle

Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with 3.3mm and 3.2mm working ends.



M NHCW
Crile-Wood Needle Holder

Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



N S16
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



HU-FRIEDY SURGICAL KIT

BASIC PERIO KIT, BLACK LINE



PART CODE	DESCRIPTION
IM1002	6-Pack Signature Series® Instrument Clip
MIR5	5 FS CS Mouth Mirror
MH6	6 Cone Socket Mirror Handle
10-130-03	3 Scalpel Handle
KO12KP03RX	1/2 Allen Orban Knife Round, Black Line
XP23/QOW6	23/QOW Expro, Wms Mkgs, #6 Handle
SG11/126	11/12 Gracey Curette, #6 Handle
SG13/146	13/14 Gracey Curette, #6 Handle
PPR3X	PR-3 Prichard Periosteal, Black Line
CM2/4X	2/4 Molt DE Curette, Black Line
S13K/13KLX	13/TG DE Curette/Chisel, Black Line
NHCW	Crile-Wood Needle Holder
IMN4121	12 Instrument Infinity Series™ Cassette, Gray
S16SCX	16 Goldman-Fox Scissors, Black Line
IMNBASPERX	Complete Basic Perio Kit, Black Line

BASIC PERIO KIT, BLACK LINE

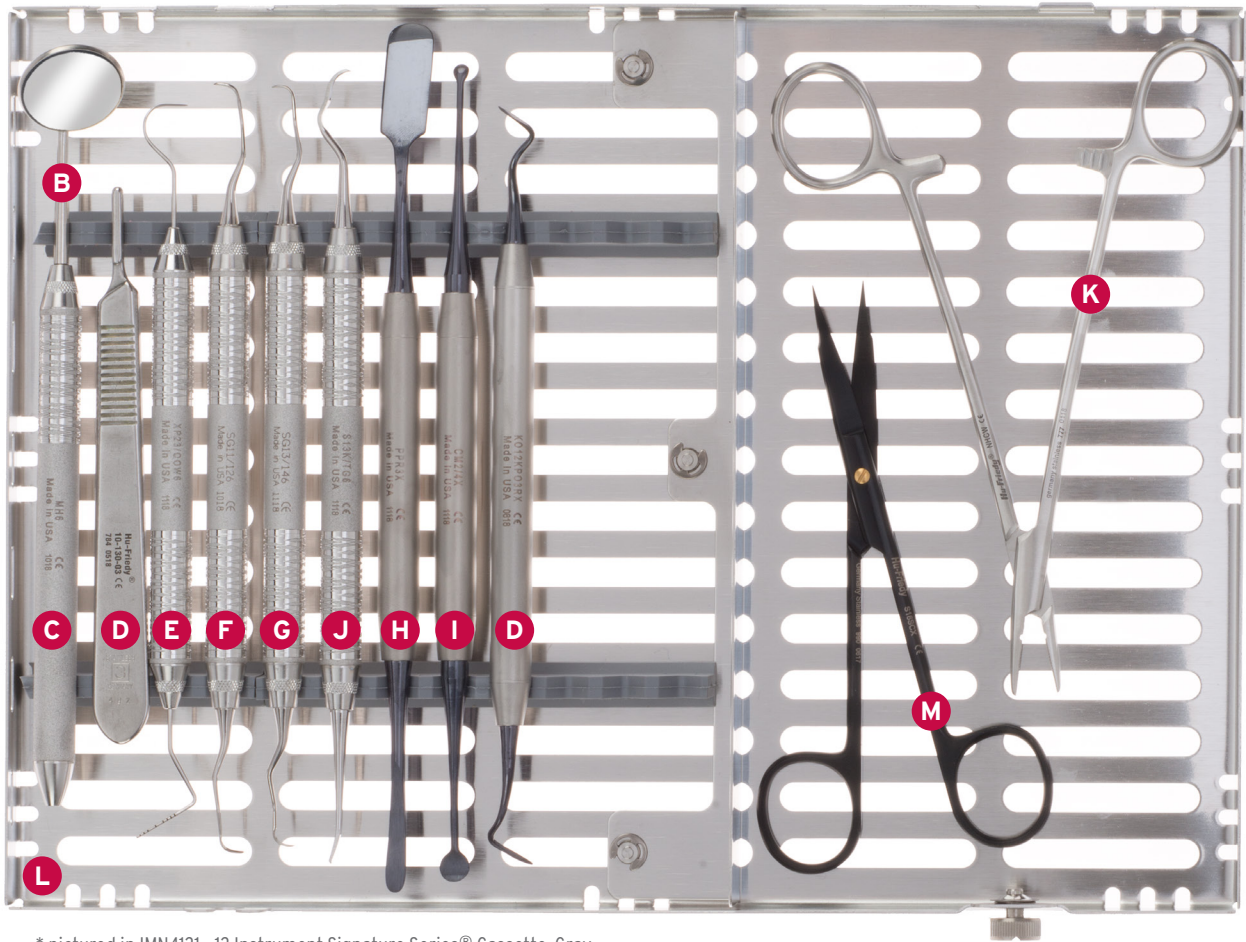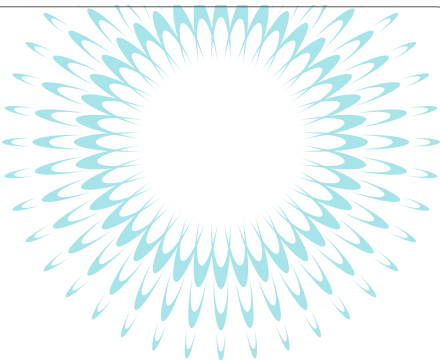

* pictured in IMN4121 - 12 Instrument Signature Series® Cassette, Gray

BASIC PERIO KIT, BLACK LINE



- A IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- B MIR5**
5 FS CS Mouth Mirror
 Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- C MH6**
6 Cone Socket Mirror Handle
 Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



- D 10-130-03**
3 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- E KO12KPO3RX**
1/2 Allen Orban Knife Round, Black Line
 20% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



- F XP23/QOW6**
23/QOW Expro, Wms Mkgs, #6 Handle
 A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- G SG11/126**
11/12 Gracey Curette, #6 Handle
 Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



- H SG13/146**
13/14 Gracey Curette, #6 Handle
 Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



- I PPR3X**
PR-3 Prichard Periosteal, Black Line
 Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



- J CM2/4X**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- K S13K/13KLX**
13/TG DE Curette/Chisel, Black Line
 Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with curved working ends 3.2mm and 4.9mm wide.



- L NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



- M IMN4121**
12 Instrument Infinity Series™ Cassette, Gray
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- N S16SCX**
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



HU-FRIEDY SURGICAL KIT

BONE GRAFTING KIT



PART CODE	DESCRIPTION
IM4128	12 Instrument Signature Series® Cassette, Blue
10-130-05	5 Scalpel Handle, Round Straight
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
SP20	20 Corn Suture Pliers
PPBUSER6	Buser Periosteal, #6 Handle
P96	9 Molt Periosteal Elevator
PPBUSE1/36	PJO/1 & PJO/3 Combined 13KL, #6 Handle
PPBUSE2/36	PJO/2 & PJO/3 Combined 13KL, #6 Handle
CL866	86 Lucas DE Surgical Curette, #6 Handle
SPR1/26	PR1/2 Prichard Curette, #6 Handle
MCUPE	Measuring cup 25ml with Marking
S16	16 Goldman-Fox Scissors
NHHB	Hegar-Baumgartner Needle Holder
PMPI	Membrane Placement Instrument
IM1002	6-Pack Signature Series® Instrument Clip
IMBGKIT	Complete Bone Grafting Kit

BONE GRAFTING KIT

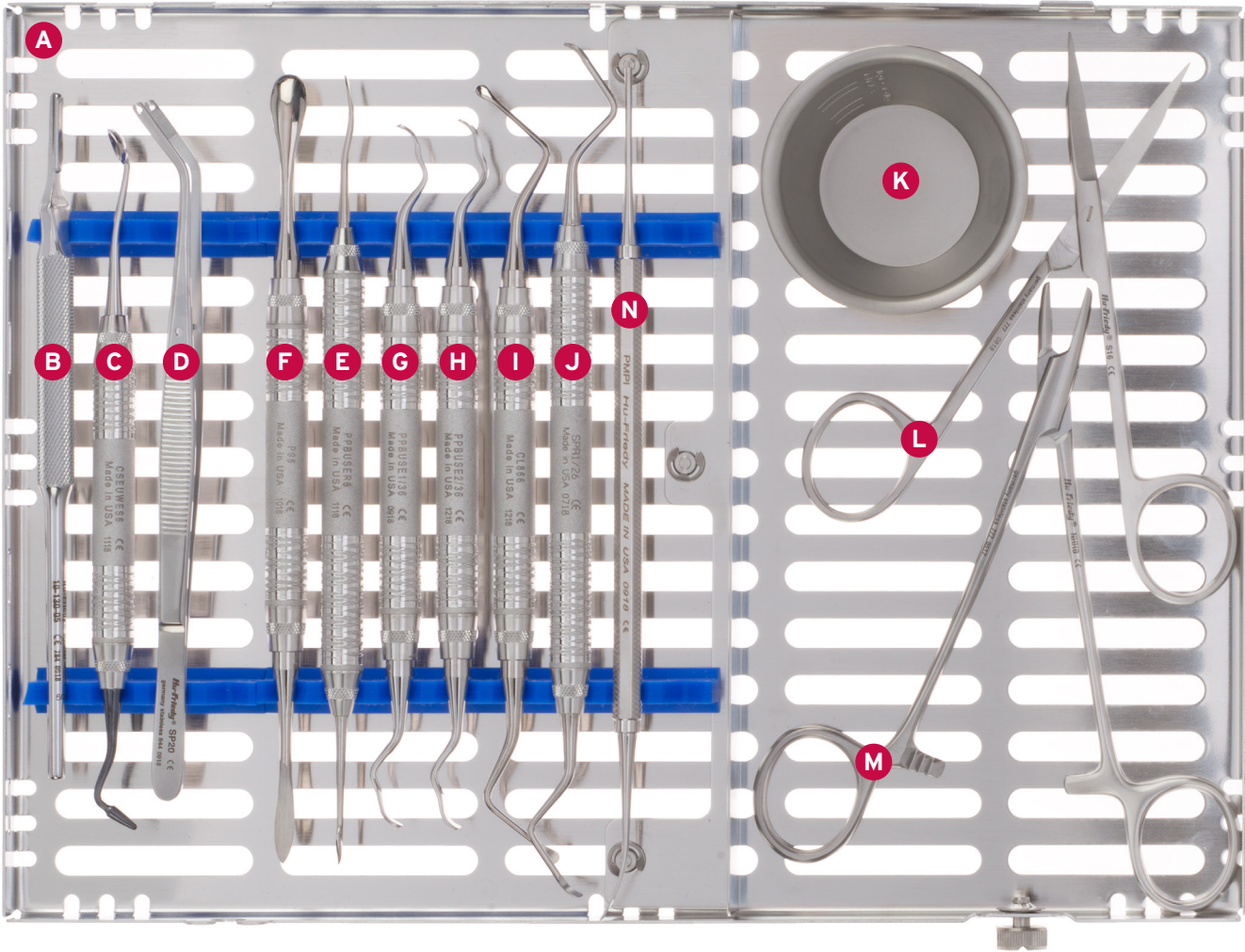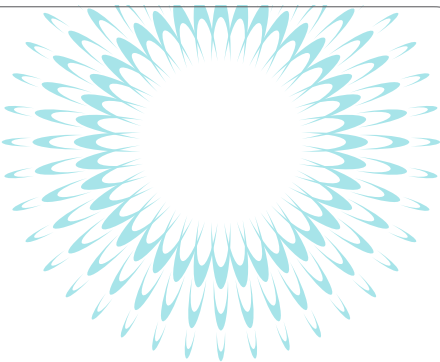

BONE GRAFTING KIT



A IM4128
12 Instrument Signature Series® Cassette, Blue

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-05
5 Scalpel Handle, Round Straight

Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



C CSEUWES6
Euwe Bone Scoop/Plugger, #6 Handle

Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



D SP20
20 Corn Suture Pliers

Ideal for use with guided tissue regeneration procedures. 15cm



E PPBUSER
Buser Periosteal, #6 Handle

Designed for reflecting and retracting the muco-periosteum after incisions of the gingival tissue. 3.8mm/4.3mm



F P96
9 Molt Periosteal Elevator

Designed for reflecting and retracting the muco-periosteum after incisions of the gingival tissue. 4.3mm/8.4mm



G PPBUSE1/36
PJO/1 & PJO/3 Combined 13KL, #6 Handle

Used to remove autologous bone for transplantation material with 4mm and 6mm wide working ends.



H PPBUSE2/36
PJO/2 & PJO/3 Combined 13KL, #6 Handle

Used to to remove autologous bone for transplantation material with 5mm and 6mm wide working ends.



I CL866
86 Lucas DE Surgical Curette, #6 Handle

Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



J SPR1/26
PR1/2 Prichard Curette, #6 Handle

Larger and heavier curette for the removal of granulation tissue and tenacious subgingival deposits measuring 7.5" long with 1.5mm wide working ends.



K MCUPE
Measuring Cup 25ml with Marking



L S16
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.




M NHHB
Hegar-Baumgartner Needle Holder

Used to grasp and guide the needle when suturing. 14 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



N PMPI
Membrane Placement Instrument

Micropoint tip for precise placement and manipulation of membranes.

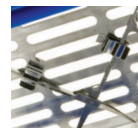

O IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



HU-FRIEDY SURGICAL KIT

BONE GRAFTING KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4128	12 Instrument Infinity Series™ Cassette, Blue
10-130-05E	5 Scalpel Handle
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
SPM20	Micro Suture Pliers
PPBUSERX	Buser Periosteal Elevator, Black Line
P9X	9 Molt Periosteal Elevator, Black Line
PPBUSER13X	PJO/1 & PJO/3 Combined, 13KL, Black Line
PPBUSER23X	PJO/2 & PJO/3 Combined, 13KL, Black Line
CL86X	86 Lucas DE Surgical Curette, Black Line
SPR1/26	PR1/2 Prichard Curette, #6 Handle
MCUPE	Measuring Cup 25ml with Markings
S16SCX	16 Goldman-Fox Scissors, Black Line
NHCW	Crile-Wood Needle Holder
PMPIX	Membrane Placement Instrument, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMNBGKITX	Complete Bone Grafting Kit, Black Line

BONE GRAFTING KIT, BLACK LINE

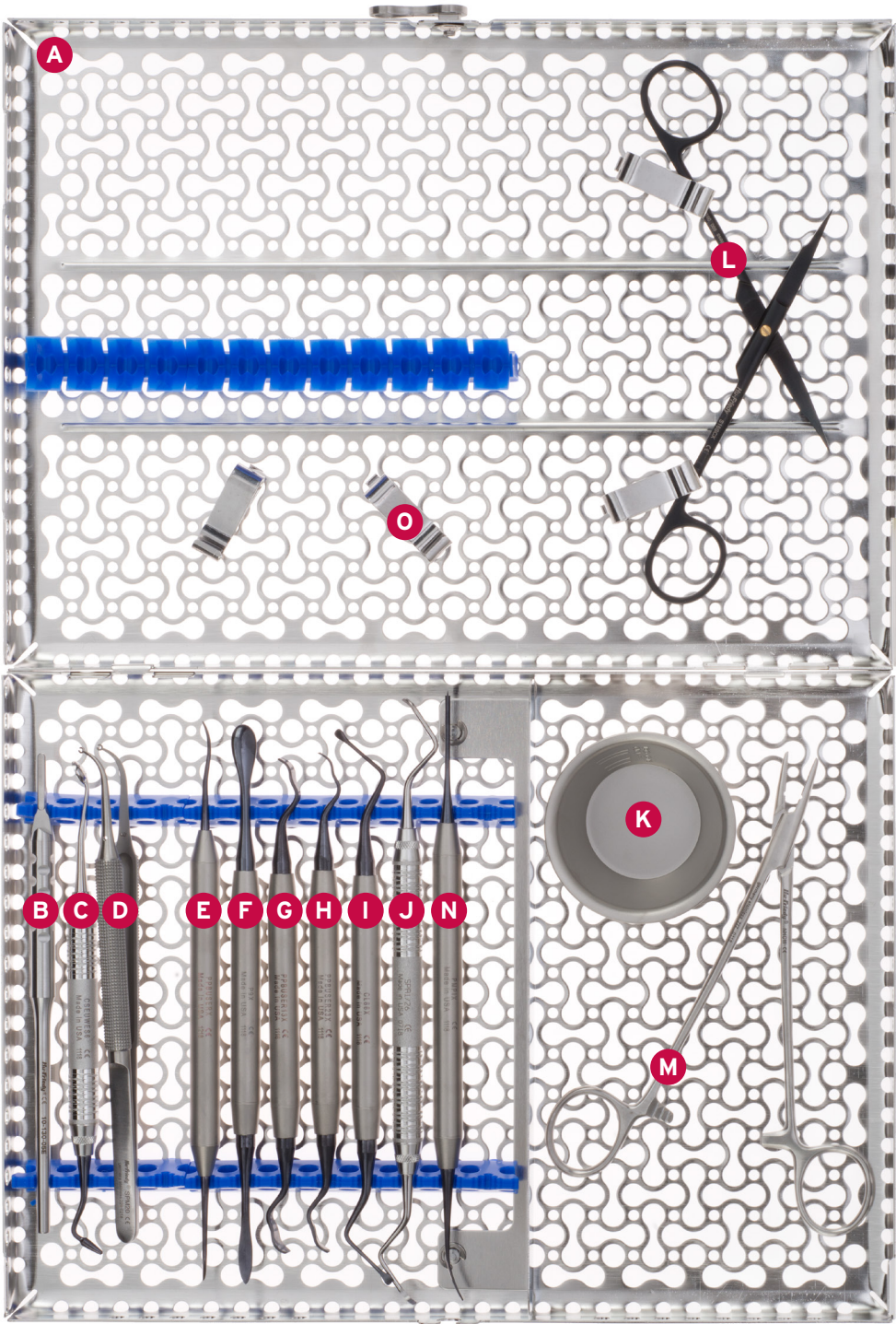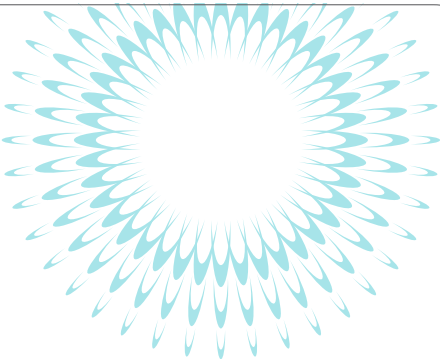

BONE GRAFTING KIT, BLACK LINE



- A IM4128**
12 Instrument Infinity Series™ Cassette, Blue
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.

10-130-05E

- B 5 Scalpel Handle**
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



CSEUWES6

- C Euwe Bone Scoop/Plugger, #6 Handle**
 Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



SPM20

- D Micro Suture Pliers**
 Ideal for use with guided tissue regeneration procedures. 15cm



PPBUSERX

- E Buser Periosteal Elevator, Black Line**
 1/3 Buser Modified Periodontal Chisel 4mm/6mm.



- F P9X**
9 Molt Periosteal Elevator, Black Line
 Featuring a wide sharp blade, non-reflective surface and improved contrast between tissue and elevator. Used for reflecting and retracting the mucoperiosteum.



- G PPBUSER13X**
PJO/1 & PJO/3 Combined 13KL, Black Line
 1/3 Buser Modified Periodontal Chisel 4mm/6mm.



- H PPBUSER23X**
PJO/2 & PJO/3 Combined 13KL, Black Line
 2/3 Buser Modified Periodontal Chisel 5mm/6mm.



- I CL86X**
86 Lucas DE Surgical Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



- J SPR1/26**
PR1/2 Prichard Curette, #6 Handle
 Larger and heavier curette for the removal of granulation tissue and tenacious subgingival deposits measuring 7.5" long with 1.5mm wide working ends.



- K MCUPE**
Measuring Cup 25ml with Marking



- L S16SCX**
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- M NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



- N PMPIX**
Membrane Placement Instrument, Black Line
 Micropoint tip for precise placement and manipulation of membranes.



- O IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



HU-FRIEDY SURGICAL KIT

CRESTAL SINUS LIFT KIT



PART CODE	DESCRIPTION
IM4162	16 Instrument Signature Series® Cassette, White
10-130-03	3 Scalpel Handle
P96	9 Molt Periosteal Elevator
CM2/4S6	2/4 Molt DE Curette, #6 Handle
TP43	43 Adson-Brown Tissue Pliers
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
PCP65006	20mm Probe 8-10-12-14-16-18-20, #6 Handle
OSTPAL1	1 Palti Osteotome 2.7mm
OSTPAL2	2 Palti Osteotome 3.2mm
OSTPAL3	3 Palti Osteotome 3.7mm
OSTPAL4	4 Palti Osteotome 4.2mm
PFIWDS1	1 Woodson DE Plastic Filling Instrument
PPR3S6	Prichard Periosteal, #6 Handle
IMCRESTKIT	Complete Crestal Sinus Lift Kit

CRESTAL SINUS LIFT KIT

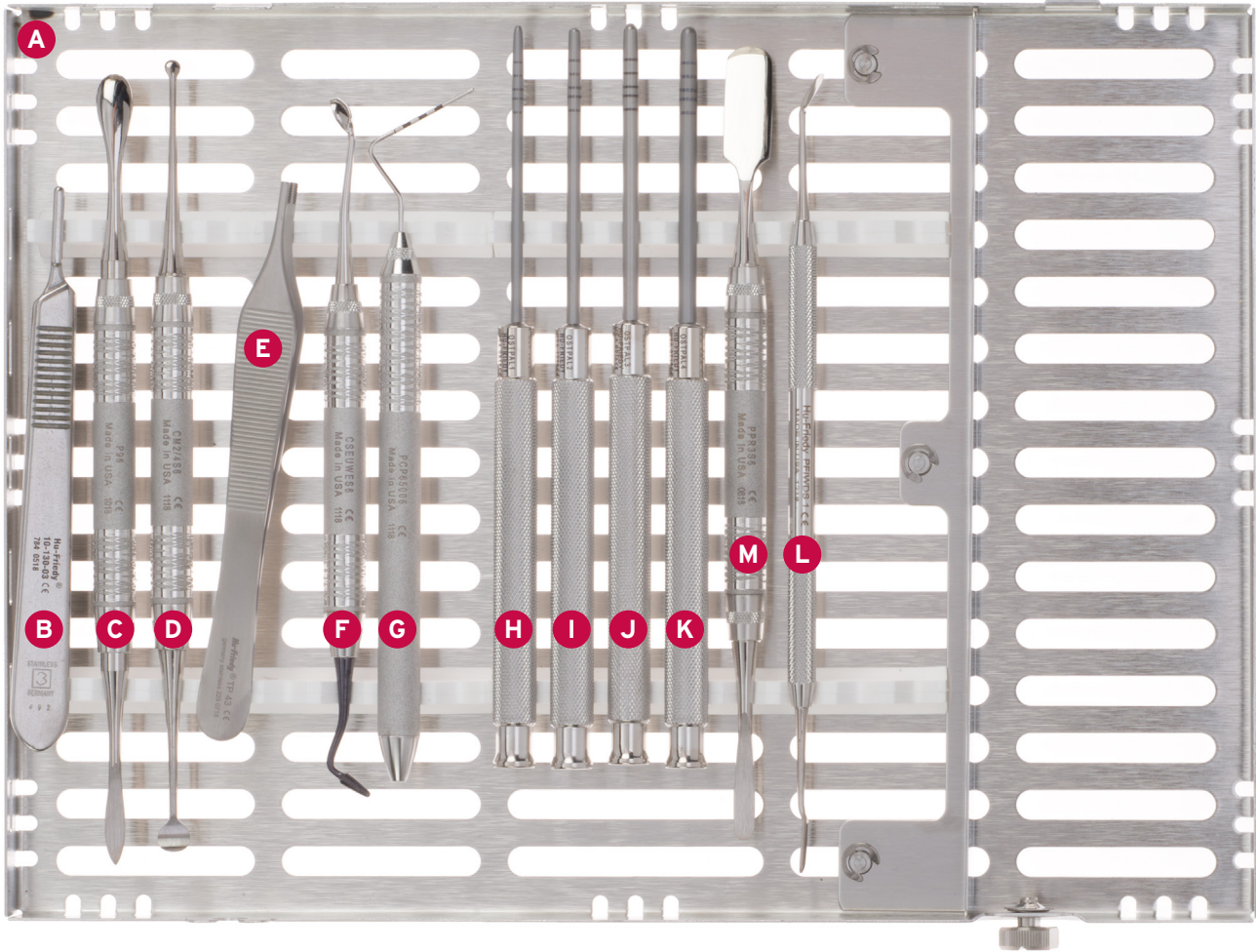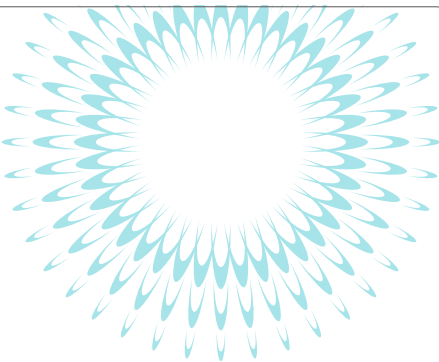

CRESTAL SINUS LIFT KIT



A IMN4162
16 Instrument Signature Series® Cassette, White

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-03
3 Scalpel Handle

Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



C P96
9 Molt Periosteal Elevator

Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D CM2/4S6
2/4 Molt DE Curette, #6 Handle

Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



E TP43
43 Adson-Brown Tissue Pliers

Used to grasp and stabilize soft tissue during suturing. 12mm



F CSEUWES6
Euwe Bone Scoop/Plugger, #6 Handle

Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



G PCP65006
20mm Probe 8-10-12-14-16-18-20, #6 Handle

A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



H OSTPAL1
1 Palti Osteotome 2.7mm

Used for preparing bone. 7-10-13-15-18mm



I OSTPAL2
2 Palti Osteotome 3.2mm

Used for preparing bone. 7-10-13-15-18mm



J OSTPAL3
3 Palti Osteotome 3.7mm

Used for preparing bone. 7-10-13-15-18mm



K OSTPAL4
4 Palti Osteotome 4.2mm

Used for preparing bone. 7-10-13-15-18mm



L PFIWDS1
1 Woodson Plastic Filling Instrument

Designed to reflect and retract the mucoperiosteum after incisions of the gingival tissue with 2.2mm and 3.8mm width.



M PPR3S6
Prichard Periosteal, #6 Handle

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



HU-FRIEDY SURGICAL KIT

CRESTAL SINUS LIFT KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4162	16 Instrument Infinity Series™ Cassette, White
10-130-03	3 Scalpel Handle
P9X	9 Molt Periosteal Elevator, Black Line
CM2/4X	2/4 Molt DE Curette, Black Line
TP5043	Adson-Brown Perma Sharp® Tissue Pliers
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
PCP65006	20mm Probe 8-10-12-14-16-18-20, #6 Handle
OSTPAL1	1 Palti Osteotome 2.7mm
OSTPAL2	2 Palti Osteotome 3.2mm
OSTPAL3	3 Palti Osteotome 3.7mm
OSTPAL4	4 Palti Osteotome 4.2mm
PFIWDS16	1 Woodson Plastic Filling, #6 Handle
IM1002	6-Pack Signature Series® Instrument Clip
NH5038	Crile-Wood Perma Sharp® Needle Holder
IMNCREKITX	Complete Crestal Sinus Lift Kit, Black Line

CRESTAL SINUS LIFT KIT, BLACK LINE

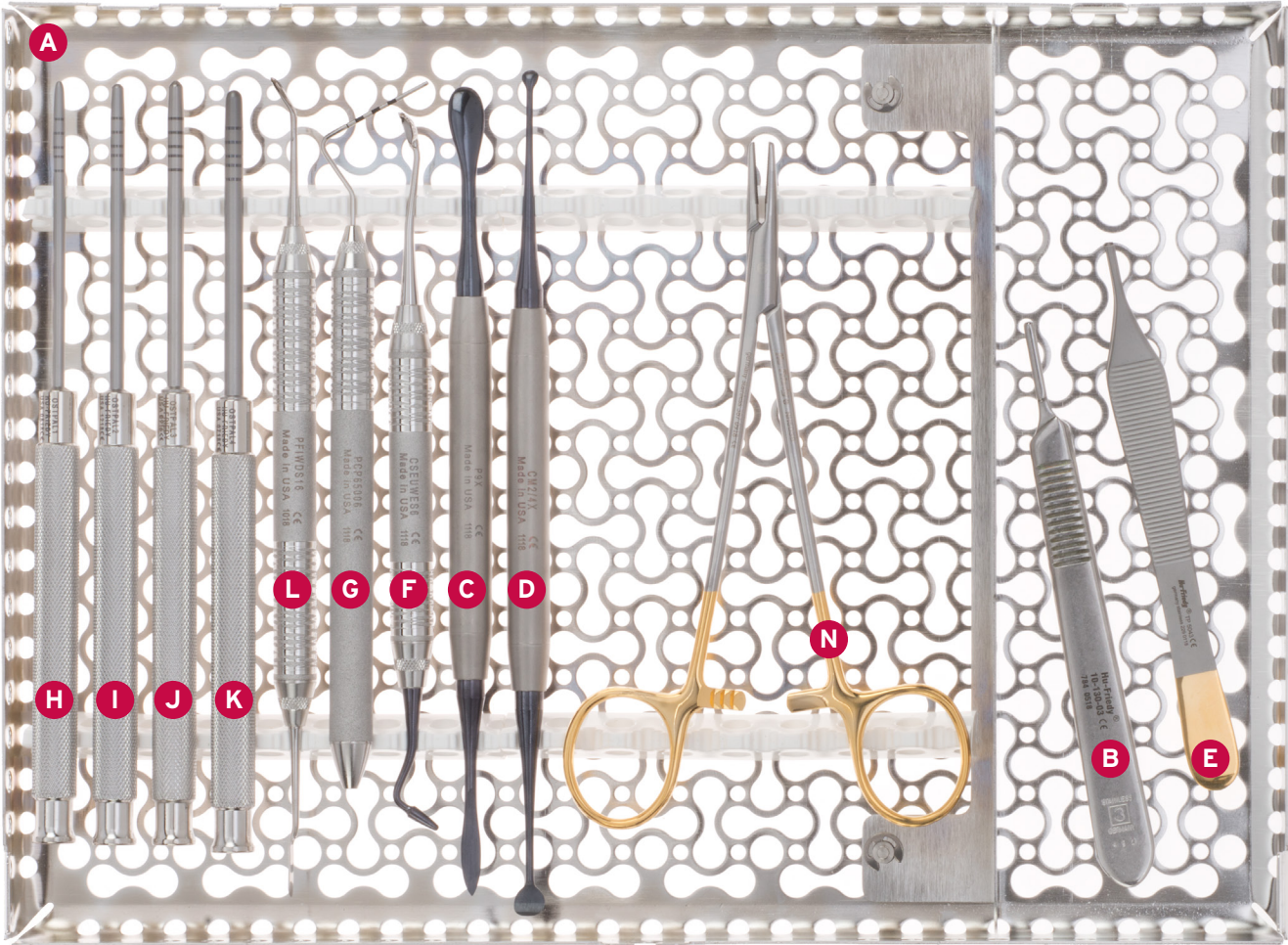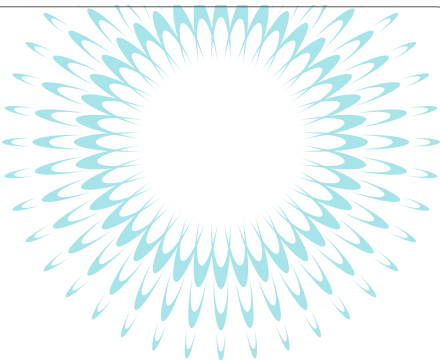

CRESTAL SINUS LIFT KIT, BLACK LINE



- A IMN4162**
16 Instrument Infinity Series™ Cassette, White
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- B 10-130-03**
3 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- C P9X**
9 Molt Periosteal Elevator, Black Line
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



- D CM2/4X**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- E TP5043**
Adson-Brown Perma Sharp® Tissue Pliers
 Used to grasp and stabilize soft tissue during suturing. 12cm



- F CSEUWES6**
Euwe Bone Scoop/Plugger, #6 Handle
 Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



- G PCP65006**
20mm Probe 8-10-12-14-16-18-20, #6 Handle
 A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- H OSTPAL1**
1 Palti Osteotome 2.7mm
 Used for preparing bone. 7-10-13-15-18mm



- I OSTPAL2**
2 Palti Osteotome 3.2mm
 Used for preparing bone. 7-10-13-15-18mm



- J OSTPAL3**
3 Palti Osteotome 3.7mm
 Used for preparing bone. 7-10-13-15-18mm



- K OSTPAL4**
4 Palti Osteotome 4.2mm
 Used for preparing bone. 7-10-13-15-18mm



- L PFIWDS16**
1 Woodson Plastic Filling, #6 Handle
 Designed to reflect and retract the mucoperiosteum after incisions of the gingival tissue with 2.2mm and 3.8mm width.



- M IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- N NH5038**
Crile-Wood Perma Sharp® Needle Holder
 Used to grasp and guide the needle when suturing. Featuring 3-0 to 5-0 carbide beaks, finger ring style.



HU-FRIEDY SURGICAL KIT

EXTRACTION KIT



PART CODE	DESCRIPTION
IM4163E2C4	16 Instrument Signature Series® Cassette, Orange
10-130-05	5 Scalpel Handle, Round Straight
P96	9 Molt Periosteal Elevator
NHMH	Mayo-Hegar Needle Holder
S16	16 Goldman-Fox Scissors
CL866	86 Lucas DE Surgical Curette, #6 Handle
E301	301 Apexo Elevator
E34S	34S Seldin Straight Elevator
BF2X6	2X Miller-Colburn Bone File
EHB13/146	13/14 Heidbrink DE Root Pick
RF	Friedman Rongeurs
TP43	43 Adson-Brown Tissue Pliers
EPTSSP	Luxating Hybrid Elevator, Straight Spade
IMEXTKIT	Complete Extraction Kit

EXTRACTION KIT

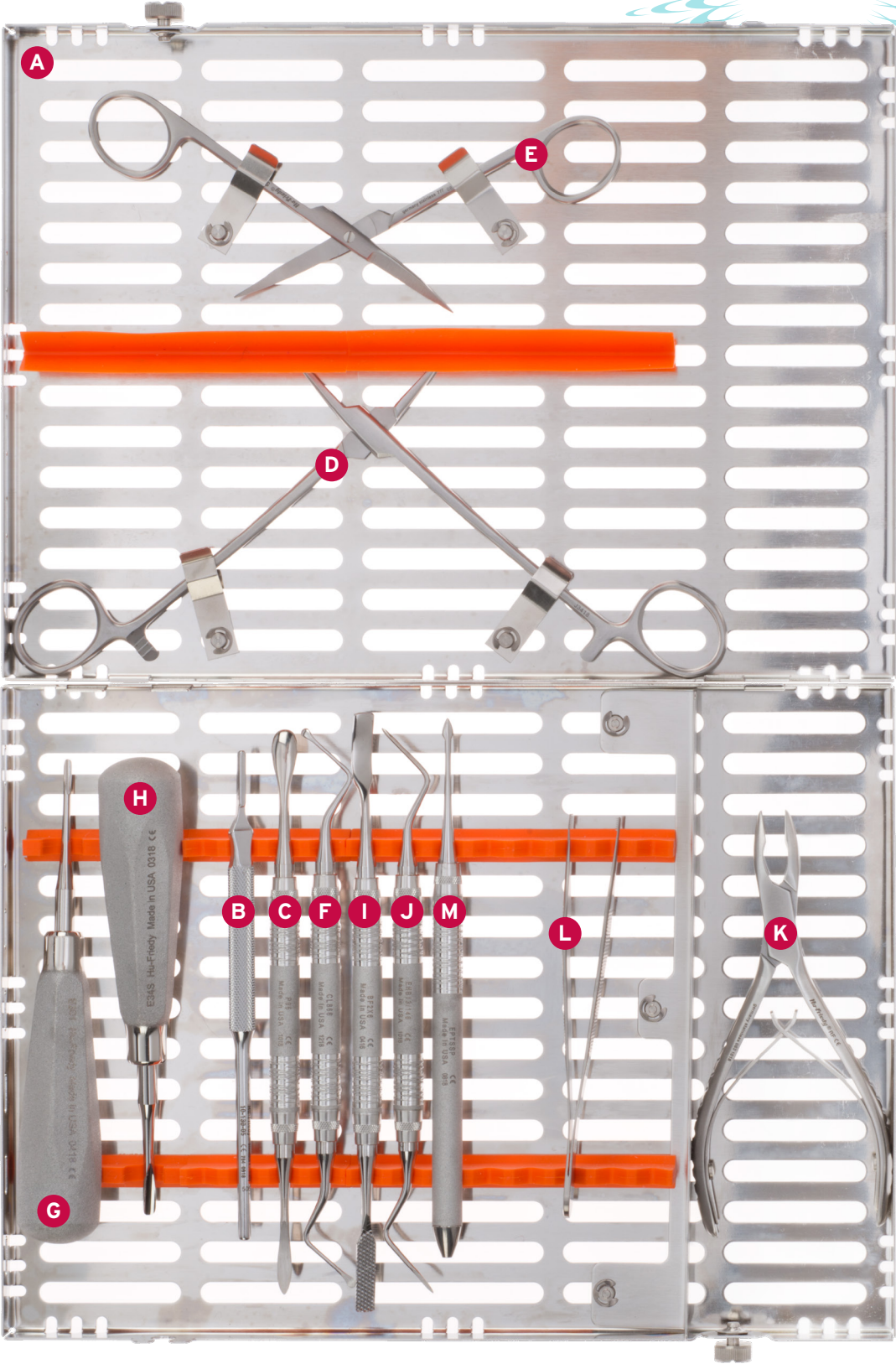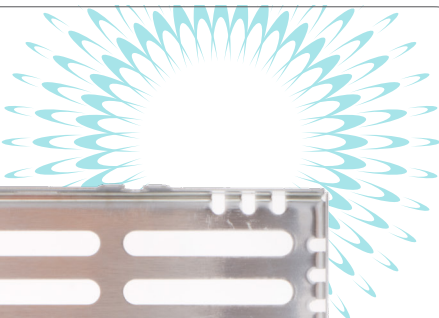

EXTRACTION KIT



A IM4163E2C4
16 Instrument Signature Series® Cassette, Orange

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-05
5 Scalpel Handle, Round Straight

Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



C P96
9 Molt Periosteal Elevator

Designed for reflecting and retracting the muco-periosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D NHMH
Mayo-Hegar Needle Holder

Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



E S16
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



F CL866
86 Lucas DE Surgical Curette, #6 Handle

Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



G E301
301 Apexo Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction.



H E34S
34S Seldin Straight Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction. 4.0mm



I BF2X6
2X Miller-Colburn Bone File

Used for final smoothing and shaping of bone with 6.2mm wide working ends.



J EHB13/146
13/14 Heidbrink DE Root Pick

Used for the removal of root tips in difficult-to-reach areas. Designed for delicate root teasing procedures. Also available in a "heavy" pattern for more substantial root tip removal procedures.



K RF
Friedman Rongeurs

Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



L TP43
43 Adson-Brown Tissue Pliers

Used to grasp and stabilize soft tissue during suturing. 12mm



M EPTSSP
Luxating Hybrid Elevator, Straight Spade

A unique cross between a periosteal and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



HU-FRIEDY SURGICAL KIT

EXTRACTION KIT, BLACK LINE



PART CODE	DESCRIPTION
IMNOS3	Oral Surgery Infinity Series™ Cassette
10-130-05E	Scalpel Handle #5
P9X	9 Molt Periosteal Elevator, Black Line
NHMH	Mayo-Hegar Needle Holder
S16SCX	16 Goldman-Fox Scissors, Black Line
CL86X	86 Lucas DE Surgical Curette, Black Line
E301	301 Apexo Elevator
E21	Mini Cryer Left Elevator
E22	Mini Cryer Right Elevator
BF2X6	2X Miller-Colburn Bone File
EHB13/14X	13/14 Heidbrink DE Root Pick, Black Line
RF	Friedman Rongeurs
TP43	43 Adson-Brown Tissue Pliers
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMNEXTKITX	Complete Extraction Kit, Black Line

EXTRACTION KIT, BLACK LINE

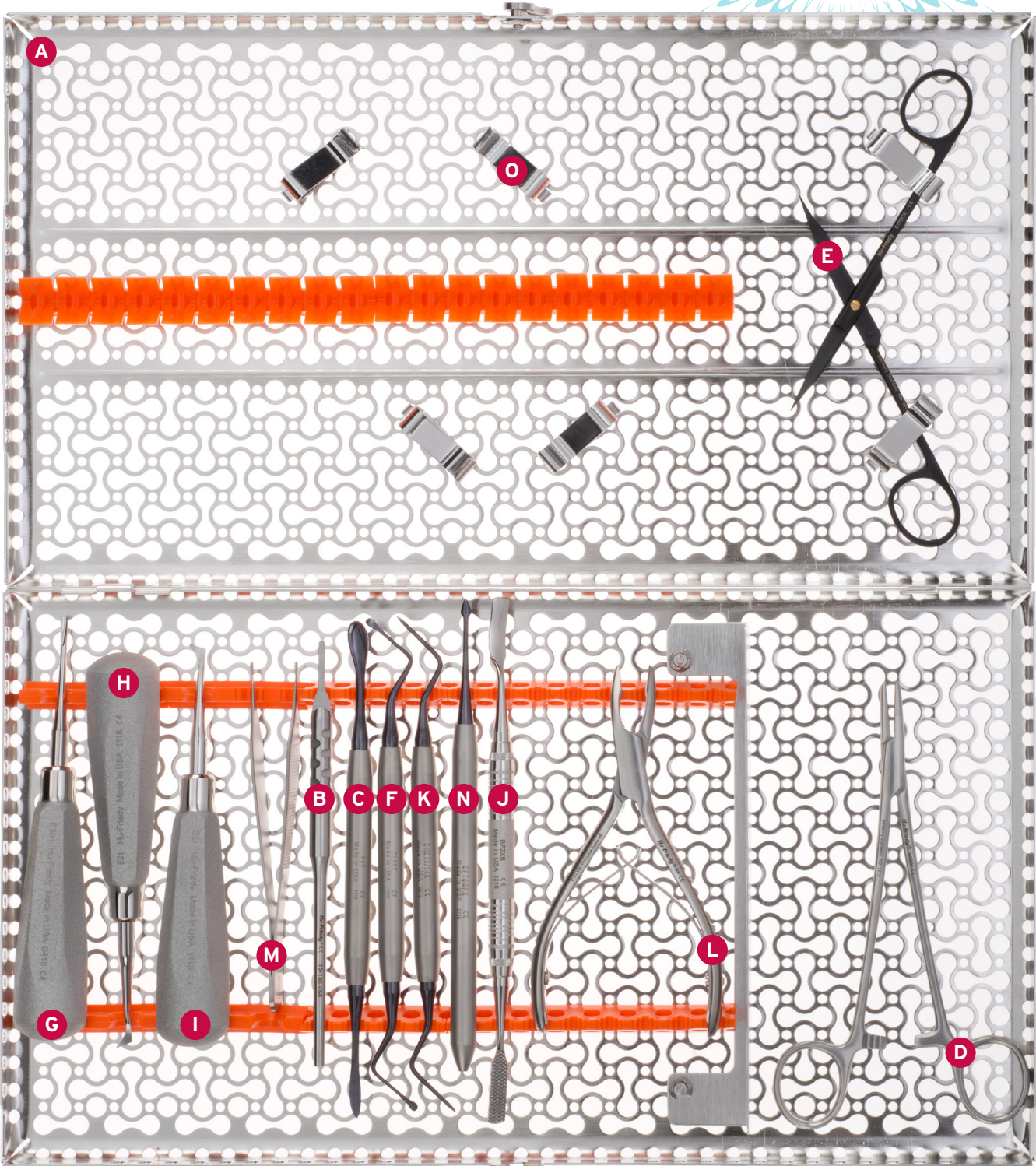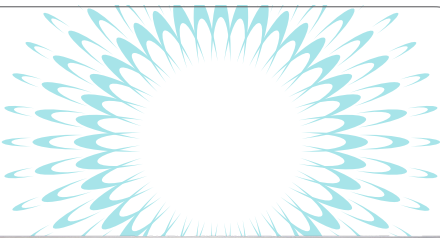

EXTRACTION KIT, BLACK LINE



A IMNOS3
Oral Surgery Infinity Series™ Cassette, Orange

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-05E
Scalpel Handle #5

Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



C P9X
9 Molt Periosteal Elevator, Black Line

Designed for reflecting and retracting the muco-periosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D NHMH
Mayo-Hegar Needle Holder

Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



E S16SCX
16 Goldman-Fox Scissors, Black Line

Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



F CL86X
86 Lucas DE Surgical Curette, Black Line

Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



G E301
PJO/2 & PJO/3 Combined 13KL, #6 Handle

Used to loosen the tooth from the periodontal ligament and ease extraction.



H E21
Mini Cryer Left Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



I E22
Mini Cryer Right Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



J BF2X6
2X Miller-Colburn Bone File

Used for final smoothing and shaping of bone with 6.2mm wide working ends.



K EHB13/14X
13/14 Heidbrink DE Root Pick, Black Line

For the removal of root tips in difficult-to-reach areas.



L RF
Friedman Rongeurs

Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



M TP43
43 Adson-Brown Tissue Pliers

Used to grasp and stabilize soft tissue during suturing. 12mm



N EPTSSPX
Luxating Hybrid Elevator, Straight Spade, Black Line

A unique cross between a periosteal and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



O IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.


HU-FRIEDY SURGICAL KIT

KEYS SOCKET GRAFTING KIT



PART CODE	DESCRIPTION
IMN4164	16 Instrument Infinity Series™ Cassette, Purple
XP23/86	23/CP-8 Expro
H2	2 Kelly Hemostat, Curved
PLGGR1X	Grisdale Bone Condenser 8-10-13-15-18, Black Line
CM2/4SX	2/4 Molt DE Curette, Black Line
10-130-05	5 Scalpel Handle, Round, Straight
PPR3X	PR-3 Prichard Periosteal, Black Line
NH5037	Halsey Perma Sharp® Needle Holder, Serrated
S5055	Metzenbaum Perma Sharp® Scissors, Curved/Blunt
TP5090	DeBakey Perma Sharp® Tissue Pliers, 6"
DP1	1 Dressing Pliers, Perry
PMP1X	Membrane Placement Instrument, Black Line
MCUPE	Measuring Cup 25 mL w/ Marking
SYBG	Syringe Bone Grafting 7mm
CL856	85 Lucas DE Surgical Curette
PJO/16	S13K/TG Jovanovic 4mm Chisel
IMNKEYSGKT	Complete Keys Socket Grafting Kit

KEYS SOCKET GRAFTING KIT

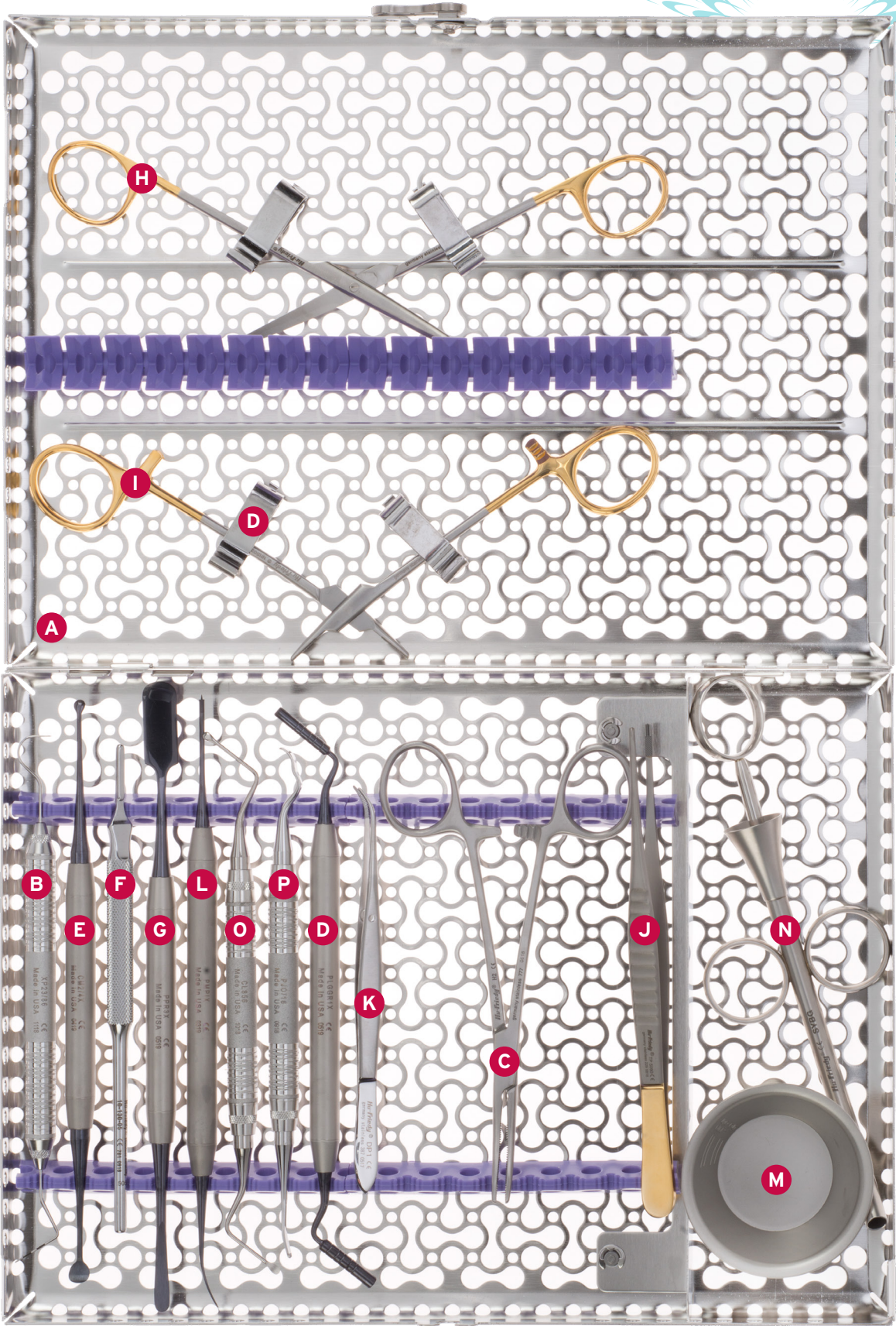

KEYS SOCKET GRAFTING KIT



- A IMN4164**
16 Instrument Infinity Series™ Cassette, Purple
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- B XP23/86**
23/CP-8 Expro
 Double-ended instrument with an explorer on one end and a probe (color-coded) on the other. ergonomic, Satin Steel™ handle.



- C H2**
2 Kelly Hemostat, Curved
 Used to clamp off blood vessels, remove small root tips and grasp loose objects. Featuring kelly curved beaks.



- D PLGGR1X**
Grisdale Bone Condenser 8-10-13-15-18, Black Line
 Used during grafting procedures with 8-11-13-15-18 mm markings. Featuring a non-reflective surface, improved contrast between tissue/scissor.



- E CM2/4SX**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- F 10-130-05**
5 Scalpel Handle, Round, Straight
 Used for cutting gingival tissue and making surgical incisions.



- G PPR3X**
PR-3 Prichard Periosteal, Black Line
 Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



- H NH5037**
Halsey Perma Sharp® Needle Holder, Serrated
 Used to grasp and guide the needle when suturing. Featuring carbide inserts, finger ring style. Use on all types of sutures.



- I S5055**
Metzenbaum Perma Sharp® Scissors, Curved/Blunt
 Used for trimming tissue or cutting sutures. 14.5cm



- J TP5090**
DeBakery Perma Sharp® Tissue Pliers, 6"
 Used to grasp and stabilize soft tissue during suturing.



- K DP1**
#1 Dressing Pliers, Perry
 Used to grip or transfer materials in and out of the oral cavity. Features a serrated tip to provides a positive grip.



- L PMPIX**
Membrane Placement Instrument, Black Line
 Micropoint tip for precise placement and manipulation of membranes.



- M MCUPE**
Measuring Cup 25 mL w/ Marking
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- N SYBG**
Syringe Bone Grafting 7mm
 Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.



- O CL856**
85 Lucas DE Surgical Curette
 Designed for curettage, cyst removal and tooth socket debridement. 2.5mm



- P PJO/16**
S13K/TG Jovanovic 4mm Chisel
 For removing and recontouring bone.



HU-FRIEDY SURGICAL KIT

LATERAL SINUS LIFT KIT



PART CODE	DESCRIPTION
IM4162	16 Instrument Signature Series® Cassette, White
IMP65786	Kramer-Nevins Sinus Lift Instr., Obtuse Angles, #6 Handle
IMP6577SC6	Kramer-Nevins Sinus Lift Instr., Acute Angles, #6 Handle
BPACK	Bone Packer DE 6-8mm
SYBG	Syringe Bone Grafting 7mm
BONT	Bone Tamper 7mm
NHCW	Crile-Wood Needle Holder
S16	16 Goldman-Fox Scissors
CRURBAN	Urban Cheek Retractor
IM1002	6-Pack Signature Series® Instrument Clip
IMLATKIT	Complete Lateral Sinus Lift Kit

LATERAL SINUS LIFT KIT

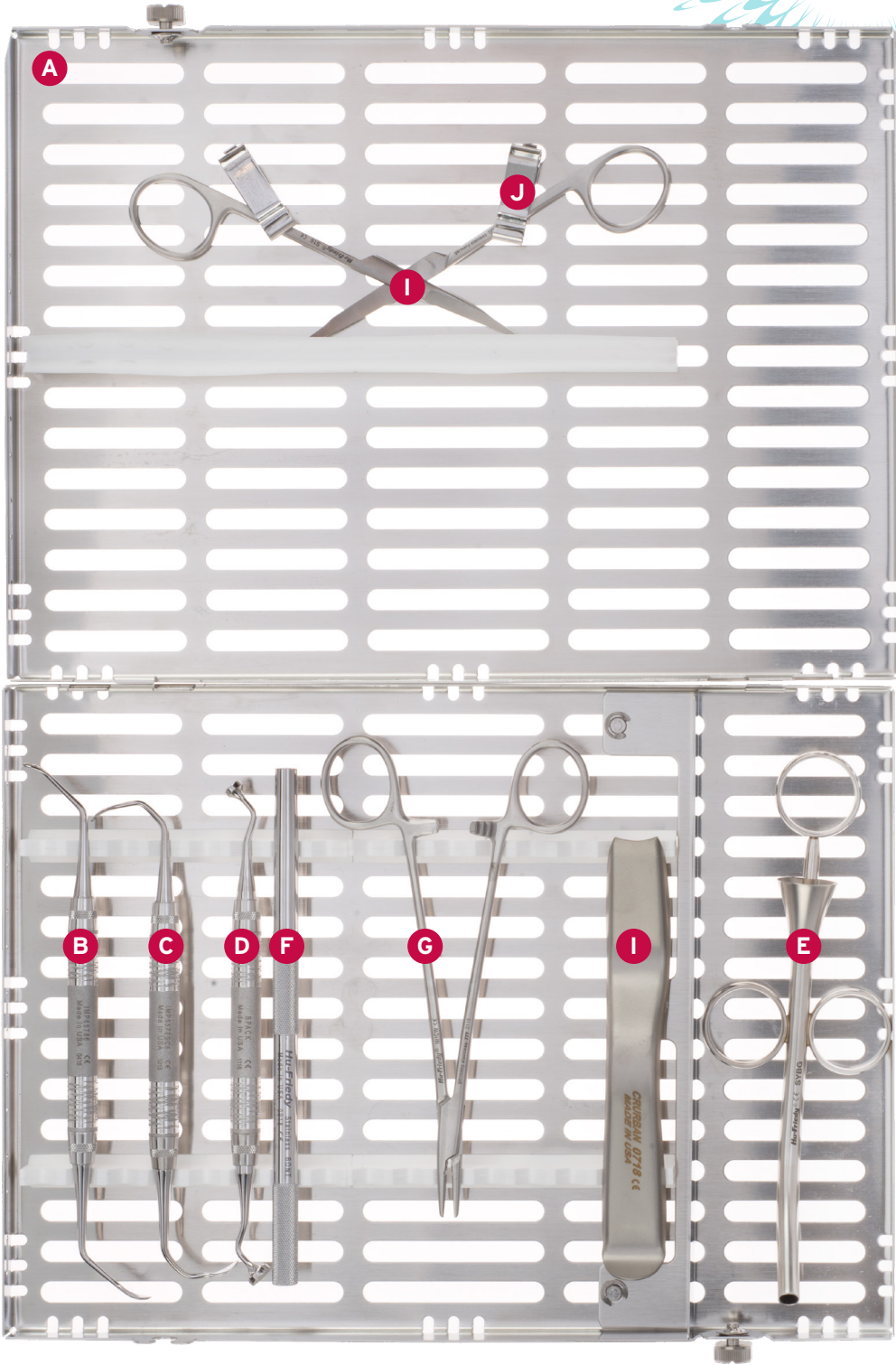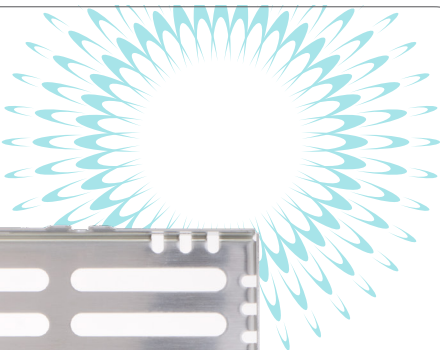

LATERAL SINUS LIFT KIT



A IM4162
16 Instrument Signature Series® Cassette, White

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B IMP65786
Kramer-Nevins Sinus Lift Instr., Obtuse Angles, #6 Handle

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With obtuse angles 4.1mm/7.6mm.



C IMP6577SC6
Kramer-Nevins Sinus Lift Instr., Acute Angles, #6 Handle

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With acute angles 4.1mm/7.6mm.



D BPACK
Bone Packer DE 6-8mm

Used to condense bone during grafting procedures with working ends 6mm and 8mm wide.



E SYBG
Syringe Bone Grafting 7mm

Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.



F BONT
Bone Tamper 7mm

Used to separate the Schneiderian membrane from surrounding bone during lateral sinus lift procedures.



G NHCW
Crile-Wood Needle Holder

Used to separate the Schneiderian membrane from surrounding bone during lateral sinus lift procedures.



H S16
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



I CRURBAN
Urban Cheek Retractor

Used to hold mucoperiosteal flaps, cheeks, lips and tongue away from the surgical area. 20.6mm



J IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



HU-FRIEDY SURGICAL KIT

LATERAL SINUS LIFT KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4162	16 Instrument Infinity Series™ Cassette, White
IMP6578SX	Kramer-Nevins Sinus Lift Instr., 47-Degree, Black Line
IMP6577SCX	Kramer-Nevins Sinus Lift Instr., Strong Curve, Black Line
BPACK	Bone Packer DE 6-8mm
SYBG	Syringe Bone Grafting 7mm
BONT	Bone Tamper 7mm
NHCW	Crile-Wood Needle Holder
S16SCX	16 Goldman-Fox Scissors, Black Line
CRURBAN	Urban Cheek Retractor
IM1002	6-Pack Signature Series® Instrument Clip
IMNLATKITX	Complete Lateral Sinus Lift Kit, Black Line

LATERAL SINUS LIFT KIT, BLACK LINE

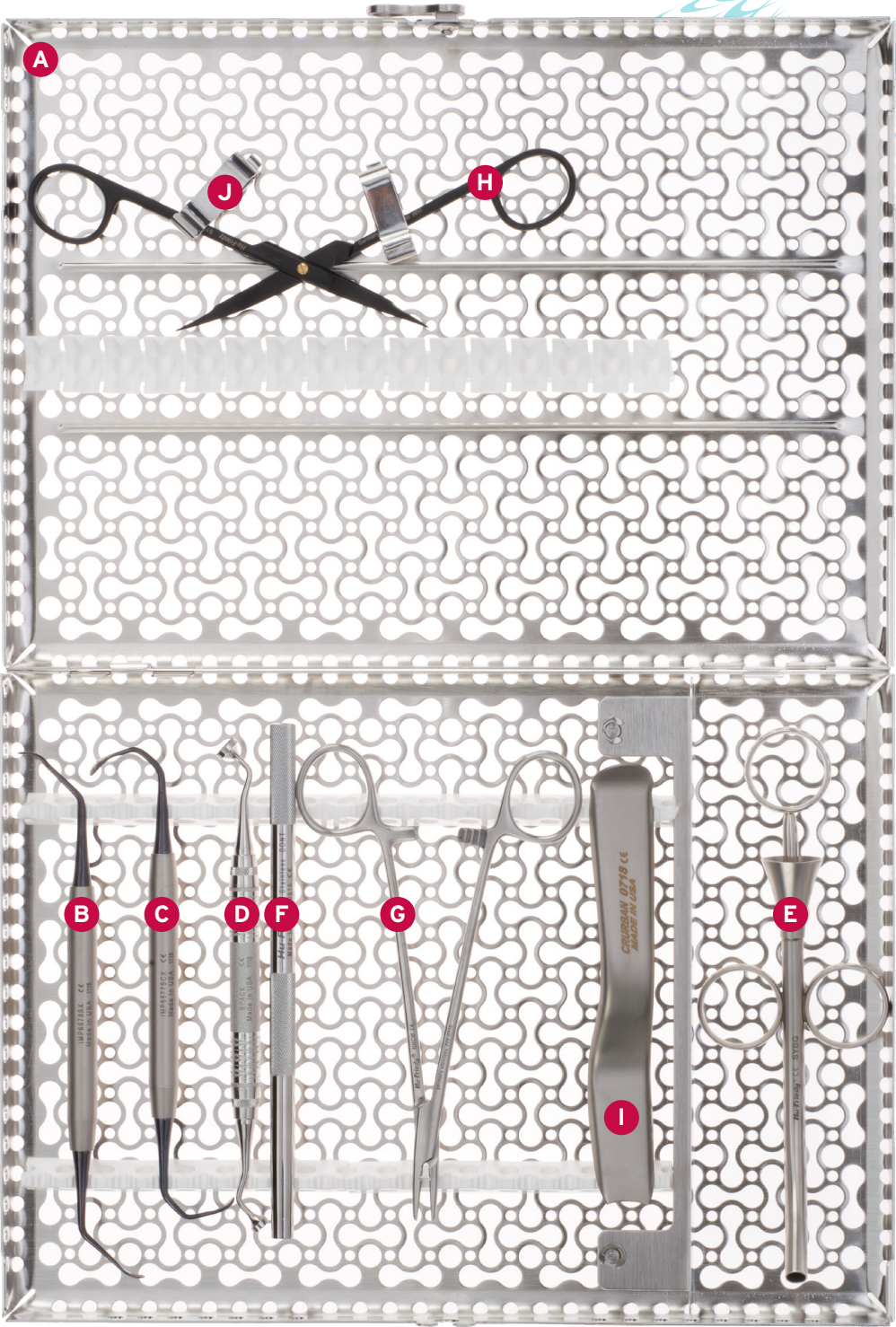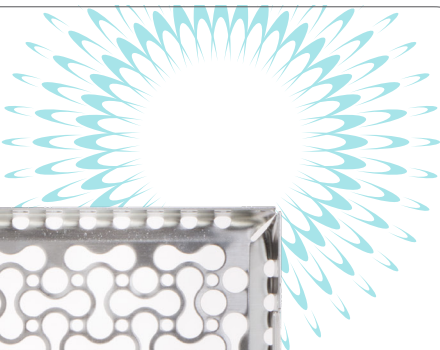

LATERAL SINUS LIFT KIT, BLACK LINE



A IM4162
16 Instrument Infinity Series™ Cassette, White

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B IMP6578SX
Kramer-Nevins Sinus Lift Instr., 47-Degree, Black Line

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With obtuse angles 4.1mm/7.6mm.



C IMP6577SCX
Kramer-Nevins Sinus Lift Instr., Strong Curve, Black Line

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With acute angles 4.1mm/7.6mm.



D BPACK
Bone Packer DE 6-8mm

Used to condense bone during grafting procedures with working ends 6mm and 8mm wide.



E SYBG
Syringe Bone Grafting 7mm

Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.



F BONT
Bone Tamper 7mm

Used to separate the Schneiderian membrane from surrounding bone during lateral sinus lift procedures.



G NHCW
Crile-Wood Needle Holder

Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



H S16SCX
16 Goldman-Fox Scissors, Black Line

Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



I CRURBAN
Urban Cheek Retractor

Used to hold mucoperiosteal flaps, cheeks, lips and tongue away from the surgical area. 20.6mm



J IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



HU-FRIEDY SURGICAL KIT

LUXATING HYBRID KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN5081	8 Instrument Infinity Series™ Cassette, Gray
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
EPTSMSX	Luxating Hybrid Elevator, Small Straight, Black Line
EPTSMCX	Luxating Hybrid - Small Curved, Black Line
EPTSMDX	Luxating Hybrid - Small Distal, Black Line
EPTSMMX	Luxating Hybrid - Small Mesial, Black Line
EPTLGDX	Luxating Hybrid - Large Distal, Black Line
EPTLGMX	Luxating Hybrid - Large Mesial, Black Line
EPTKITX	Complete Luxating Hybrid Kit, Black Line

LUXATING HYBRID KIT, BLACK LINE

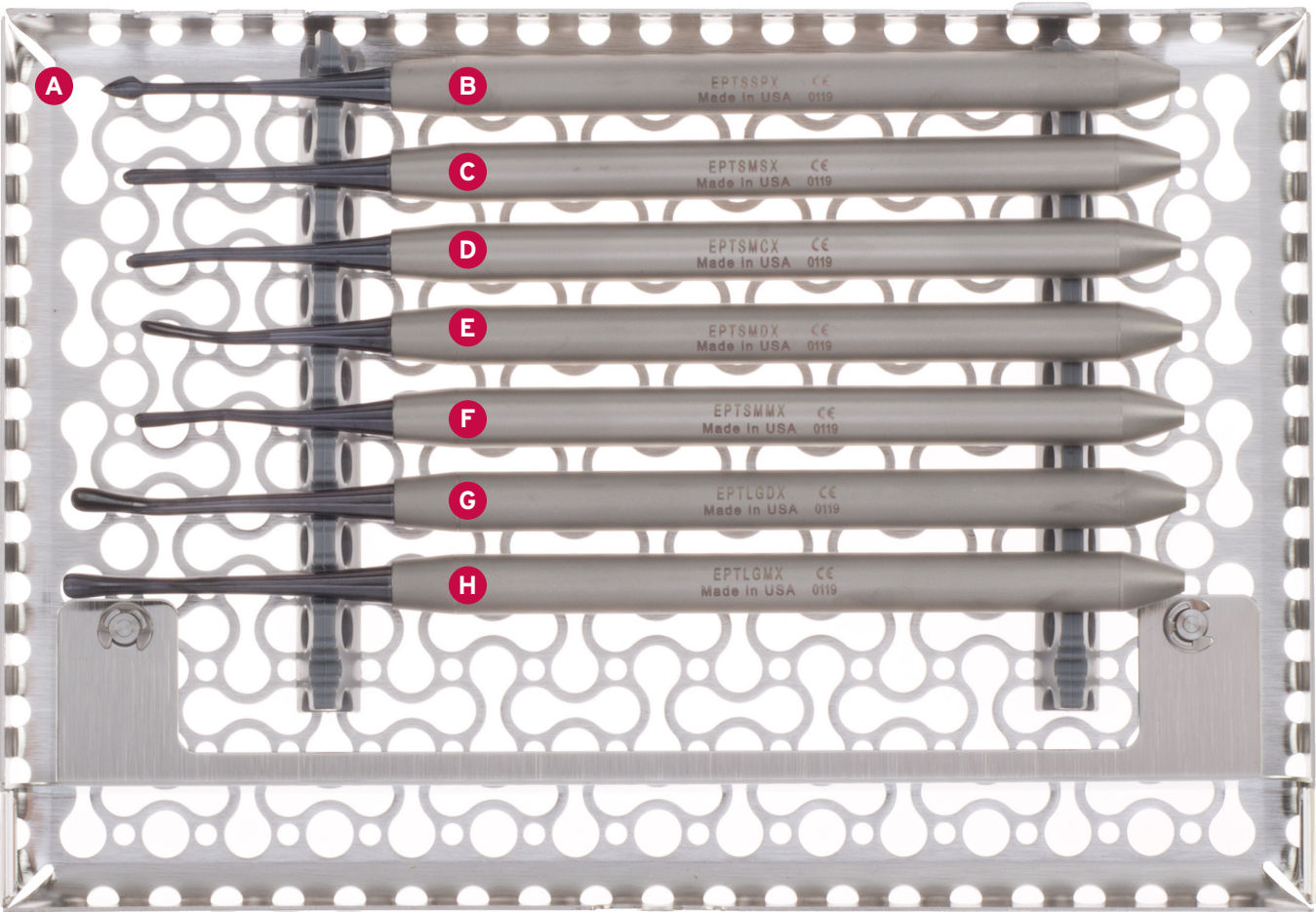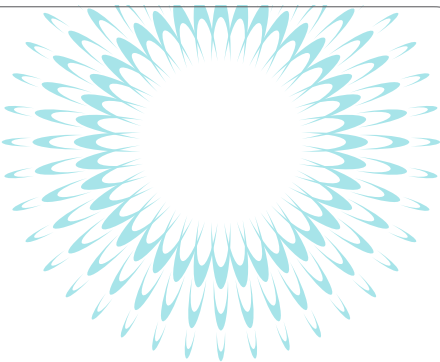

LUXATING HYBRID KIT, BLACK LINE



- A IMN5081**
8 Instrument Infinity Series™ Cassette, Gray
IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- B EPTSSPX**
Luxating Hybrid Elevator, Straight Spade, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSMSX**
Luxating Hybrid Elevator, Small Straight, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D EPTSMCX**
Luxating Hybrid - Small Curved, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- E EPTSMDX**
Luxating Hybrid - Small Distal, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- F EPTSMMX**
Luxating Hybrid - Small Mesial, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- G EPTLGDX**
Luxating Hybrid - Large Distal, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- H EPTLGMX**
Luxating Hybrid - Large Mesial, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



HU-FRIEDY SURGICAL KIT

MICROSURGICAL KIT



PART CODE	DESCRIPTION
10-130-70	7 Mini Scalpel Blade Holder
S13K/13KL6	13K/13KL DE Surgical Chisel
C36/376	36/37 Back-Action DE Chisel
8-905DD	Precision Dressing Forceps, Diamond Dusted
NH5024CORT	Cortellini Micro Castro Needle Holder
PPBUSER6	Buser Periosteal, #6 Handle
SPV	Microsurgical Scissors 18cm
IM8120	12 Instrument Signature Series® Cassette, Red
PCPUNC156	CC SE UNC Probe, 1-15, #6 Handle
MIR5	5 FS CS Mouth Mirror
MH6	Cone Socket Mirror Handle, #6 Handle
IMMICSUR	Complete Microsurgical Kit

MICROSURGICAL KIT

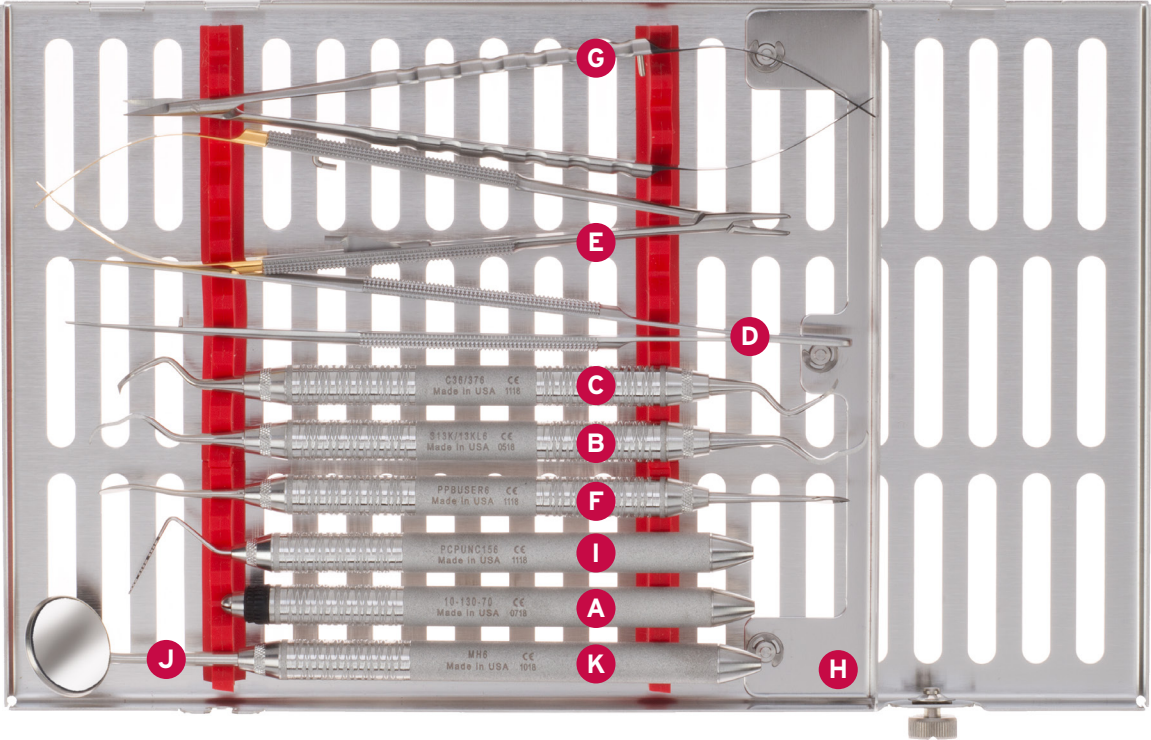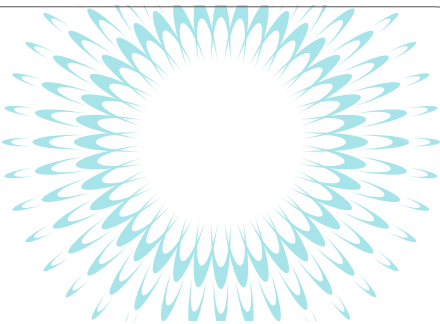

MICROSURGICAL KIT



- A 10-130-70**
7 Mini Scalpel Blade Holder
 Designed to hold micro scalpel blades for minimally invasive procedures and microsurgery



- B S13K/13KL6**
13K/13KL DE Surgical Chisel
 Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with curved working ends 3.2mm and 4.9mm wide.



- C C36/376**
36/37 Back-Action DE Chisel
 Designed for removing bone adjacent to the tooth without causing trauma and is especially useful on the distal of last molars with working ends 3.8mm wide.



- D 8-905DD**
Precision Dressing Forceps, Diamond Dusted
 Used to grasp and stabilize soft tissue, 7" length.



- E NH5024CORT**
Cortellini Micro Castro Needle Holder
 Used to grasp and guide the needle when suturing.



- F PPBUSER6**
Buser Periosteal, Satin #6 Handle
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 3.8mm/4.3mm



- G SPV**
Microsurgical Scissors
 Microsurgical scissors unique design features allow extremely smooth cutting of fine and coarse tissue with reduced trauma. 18cm



- H IM8120**
12 Instrument Signature Series® Cassette, Red
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- I PCPUNC156**
CC SE UNC Probe, 1-15, #6 Handle
 Single-ended, color-coded probe. ergonomic, Satin Steel™ handle. University of North Carolina.



- J MIR5**
5 FS CS Mouth Mirror
 Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- K MH6**
Cone Socket Mirror Handle, #6 Handle
 Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



HU-FRIEDY SURGICAL KIT

MICROSURGICAL KIT, BLACK LINE




PART CODE	DESCRIPTION
10-130-70	7 Mini Scalpel Blade Holder
S13K/13KL9	13K/13KL DE Surgical Chisel, #9 Handle
C36/379	36/37 Back-Action DE Chisel, #9 Handle
8-905DD	Precision Dressing Forceps, Diamond Dusted
NH5024CORT	Cortellini Micro Castro Needle Holder
PPBUSERX	Buser Periosteal Elevator, Black Line
SPVX	Castro Microsurgical Scissors, Black Line
IMN4121	12 Instrument Infinity Series™ Cassette, Gray
PCPUNC156	CC SE UNC Probe, 1-15, #6 Handle
MIR5	5 FS CS Mouth Mirror
MH6	Cone Socket Mirror Handle, #6 Handle
IMNMICSURX	Complete Microsurgical Kit, Black Line

MICROSURGICAL KIT, BLACK LINE

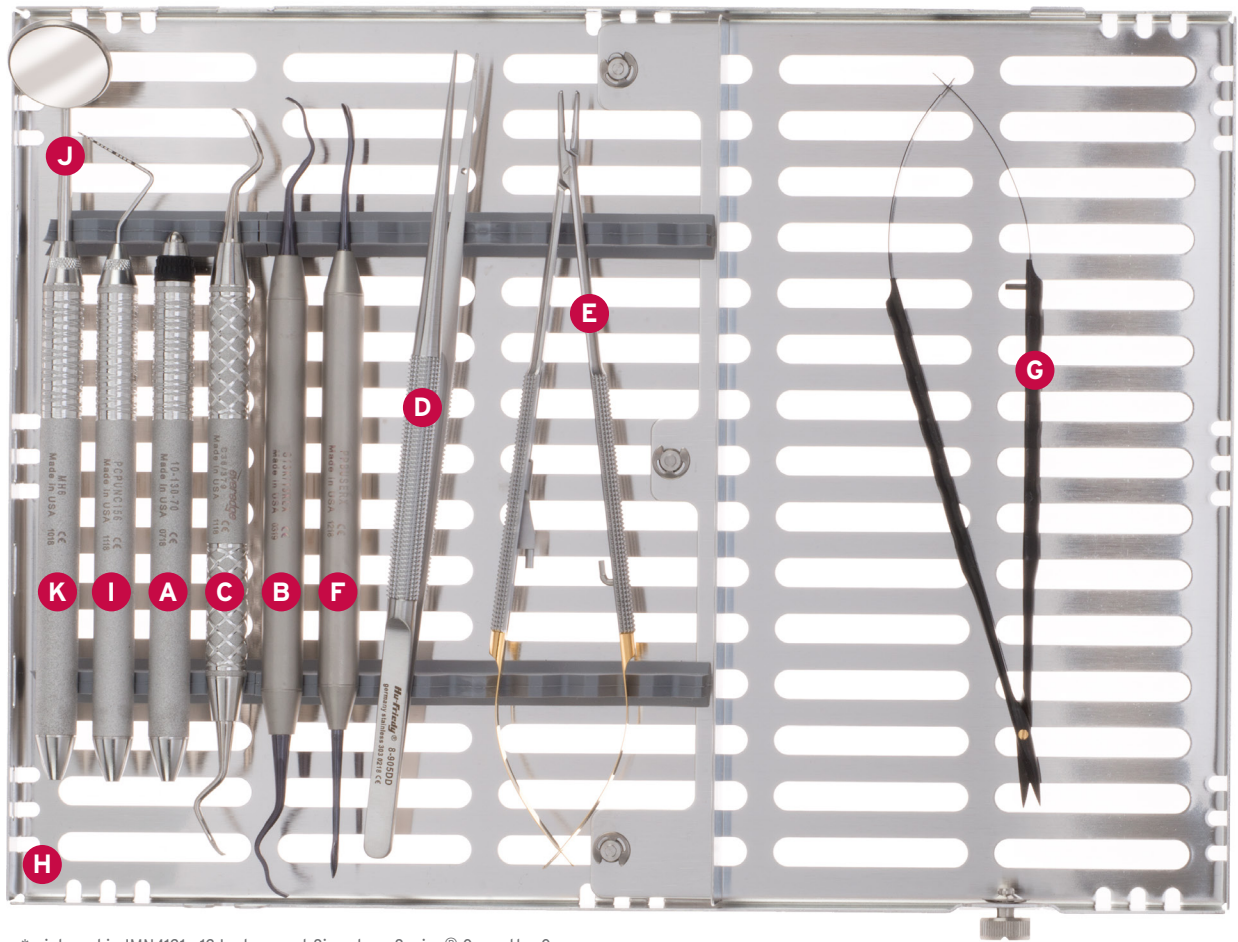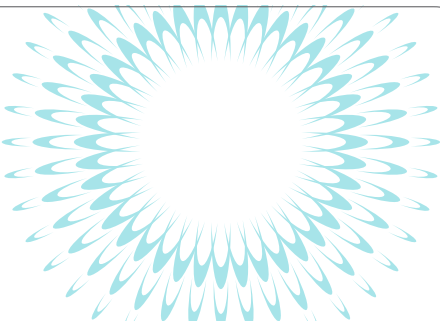

* pictured in IMN4121 - 12 Instrument Signature Series® Cassette, Gray

MICROSURGICAL KIT BLACK LINE

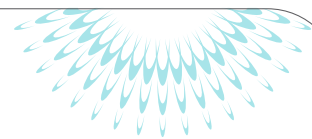

- A 10-130-70**
7 Mini Scalpel Blade Holder
 Designed to hold micro scalpel blades for minimally invasive procedures and microsurgery



- B S13K/13KL6**
13K/13KL DE Surgical Chisel
 Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with curved working ends 3.2mm and 4.9mm wide.



- C C36/376**
36/37 Back-Action DE Chisel
 Designed for removing bone adjacent to the tooth without causing trauma and is especially useful on the distal of last molars with working ends 3.8mm wide.



- D 8-905DD**
Precision Dressing Forceps, Diamond Dusted
 Used to grasp and stabilize soft tissue, 7" length.



- E NH5024CORT**
Cortellini Micro Castro Needle Holder
 Used to grasp and guide the needle when suturing.



- F PPBUSERX**
Buser Periosteal, #6 Handle
 Featuring a sharp round blade and narrow pointed end for interproximal use. For reflecting and retracting the mucoperiosteum.



- G SPVX**
Microsurgical Scissors
 Used for trimming tissue or cutting sutures. 16.5cm



- H IMN4121**
12 Instrument Infinity Series™ Cassette, Gray
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.





- I PCPUNC156**
CC SE UNC Probe, 1-15, #6 Handle
 Single-ended, color-coded probe. ergonomic, Satin Steel™ handle. University of North Carolina.



- J MIR5**
5 FS CS Mouth Mirror
 Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- K MH6**
Cone Socket Mirror Handle, Satin Steel
 Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.

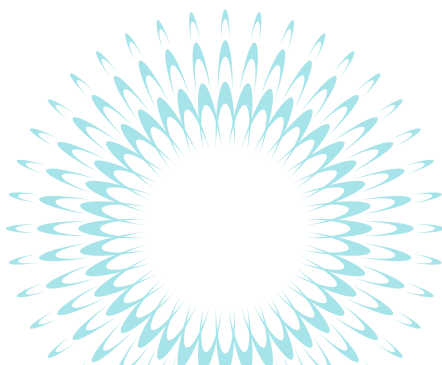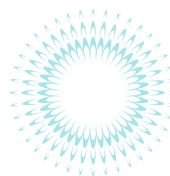



HU-FRIEDY SURGICAL KIT

RESCUE KIT



PART CODE	DESCRIPTION
EPTSSP	Luxating Hybrid Elevator, Straight Spade
EPTSMC	Luxating Hybrid, Small Curved
EPTSMS	Luxating Hybrid Elevator, Small Straight
EHB13/146	13/14 Heidbrink DE Root Pick
CL866	86 Lucas DE Surgical Curette #6 Handle
RMF	Mini Friedman Rongeurs
TRE040M	Trephine Bur 4mm Inn. Diam. Right Angle
TRE060M	Trephine Bur 6mm Inn. Diam. Right Angle
IMS-1372SH	IMS 1/2 Bur Cushion - Short Lid
IM4160	16 Instrument Signature Series® Cassette, Red
IMRESCUE	Complete Rescue Kit

RESCUE KIT

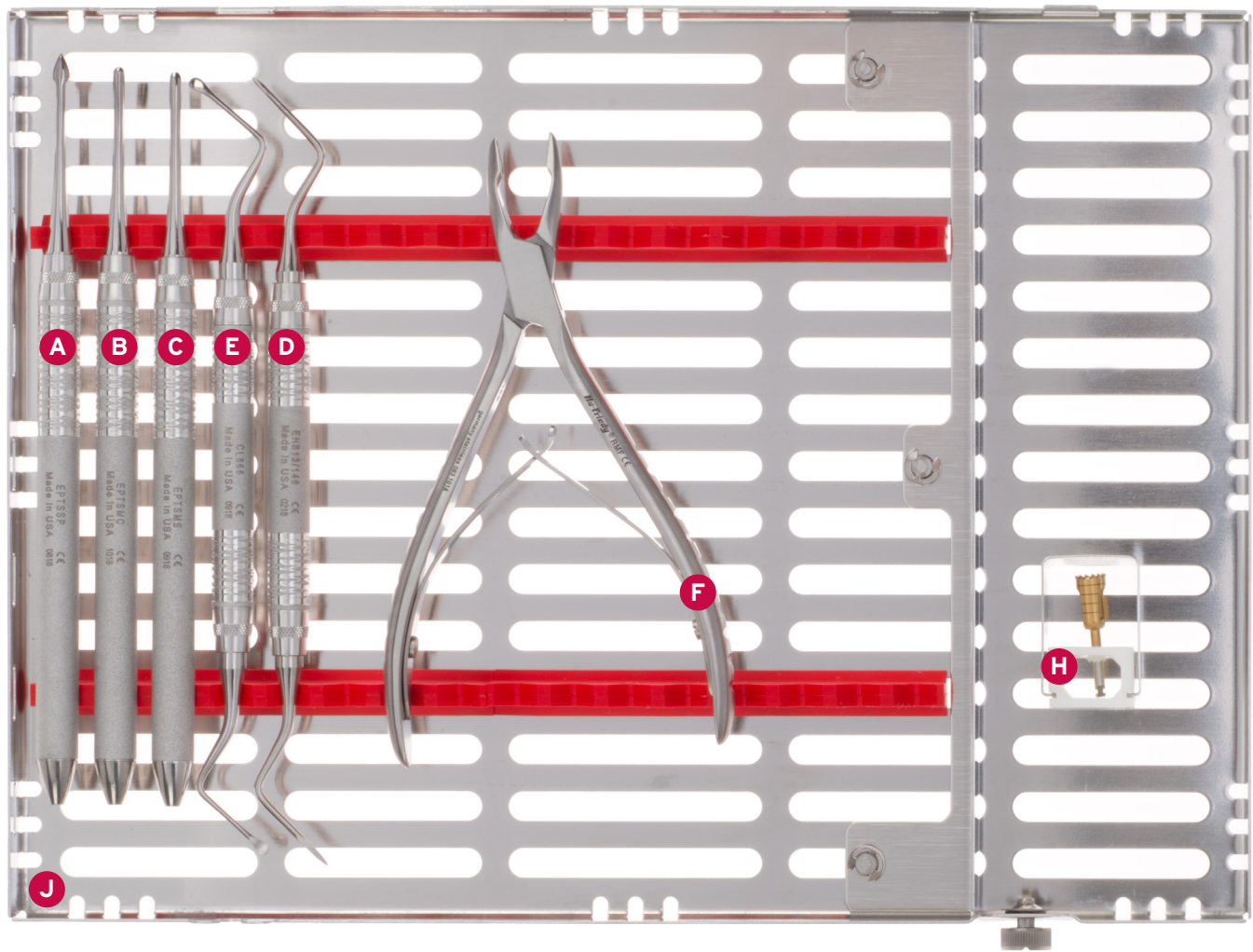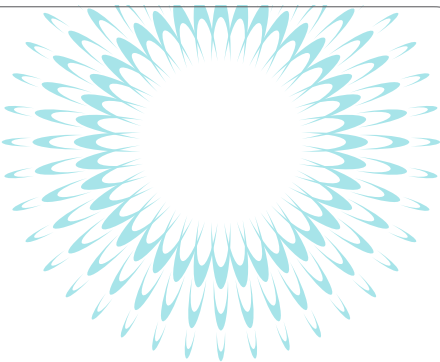

G Trephine Bur 4mm Inn. Diam Right Angle *[not pictured]*

RESCUE KIT



- A EPTSSP**
Luxating Hybrid Elevator, Straight Spade
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- B EPTSMC**
Luxating Hybrid, Small Curved
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSMS**
Luxating Hybrid Elevator, Small Straight
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D EHB13/146**
13/14 Heidbrink DE Root Pick
 Used for the removal of root tips in difficult-to-reach areas. Designed for delicate root teasing procedures.



- E CL866**
86 Lucas DE Surgical Curette #6 Handle
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



- F RMF**
Mini Friedman Rongeurs
 Used for trimming and recontouring alveolar bone and gross tissue removal.



- G TRE040M**
Trephine Bur 4mm Inn. Diam Right Angle
 Used primarily for harvesting bone cores.



- H TRE060M**
Trephine Bur 6mm Inn. Diam. Right Angle
 Used primarily for harvesting bone cores.



- I IMS-1372SH**
IMS 1/2 Bur Cushion - Short Lid
 Designed to hold burs securely in place during cleaning and sterilization. Hollow design allows drainage and proper cleaning.



- J IM4160**
16 Instrument Signature Series® Cassette, Red
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.


HU-FRIEDY SURGICAL KIT

RESCUE KIT, BLACK LINE



PART CODE	DESCRIPTION
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
EPTSMCX	Luxating Hybrid, Small Curved, Black Line
EPTSMSX	Luxating Hybrid Elevator, Small Straight, Black Line
EHB13/14X	13/14 Heidbrink DE Root Pick, Black Line
CL86X	86 Lucas DE Surgical Curette, Black Line
RMF	Mini Friedman Rongeurs
IMN8120	12 Instrument Infinity Series™ Cassette, Red
IMNRESCUEX	Complete Rescue Kit, Black Line

RESCUE KIT, BLACK LINE

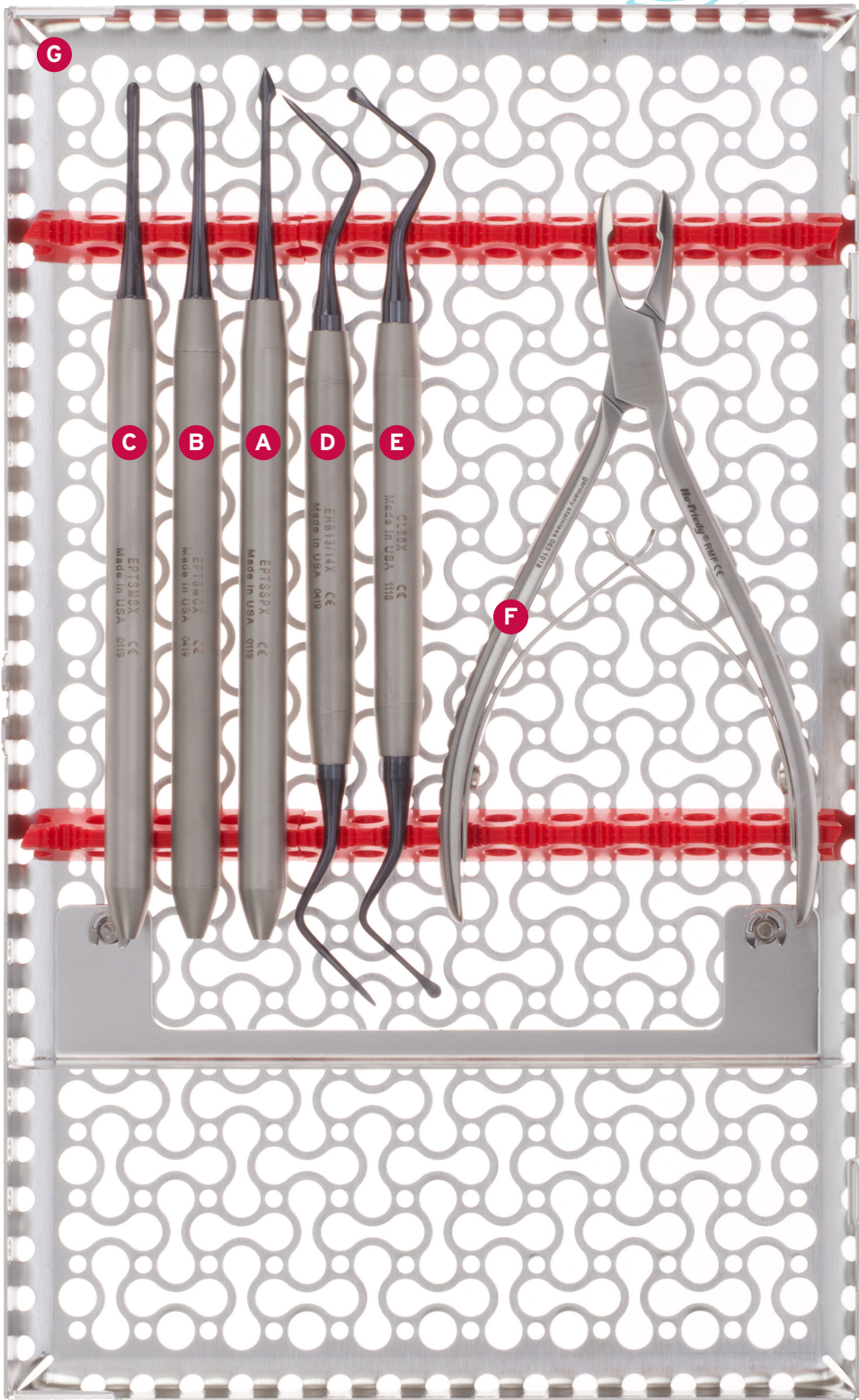

RESCUE KIT, BLACK LINE



A EPTSSPX
Luxating Hybrid Elevator, Straight Spade, Black Line

A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



B EPTSMCX
Luxating Hybrid, Small Curved, Black Line

A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



C EPTSMSX
Luxating Hybrid Elevator, Small Straight, Black Line

A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



D EHB13/14X
13/14 Heidbrink DE Root Pick, Black Line

Used for the removal of root tips in difficult-to-reach areas.




E CL86X
86 Lucas DE Surgical Curette #6 Handle, Black Line

Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



F RMF
Mini Friedman Rongeurs

Used for trimming and recontouring alveolar bone and gross tissue removal. 1.5mm/14.5cm

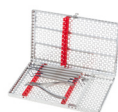

G IMN8120
12 Instrument Infinity Series™ Cassette, Red

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.