



HU-FRIEDY SURGICAL KIT

ATRAUMATIC EXTRACTION KIT



PART CODE	DESCRIPTION
EPTSMS	Luxating Hybrid Elevator, Small Straight
EPTSMC	Luxating Hybrid, Small Curved
EPTSSP	Luxating Hybrid Elevator, Straight Spade
FAF150	Apical Forceps 150
FAF151	Apical Forceps 151
CM2/4S6	2/4 Molt DE Curette, #6 Handle
NHCW	Crile-Wood Needle Holder
S16SC	16 G-F Scissors Super Cut
10-130-05	5 Scalpel Handle, Round Straight
P96	9 Molt Periosteal Elevator
IM1002	6-Pack Signature Series® Instrument Clip
IM4161	16 Instrument Signature Series® Cassette, Gray
IMATRAEXT	Complete Atraumatic Extraction Kit

ATRAUMATIC EXTRACTION KIT

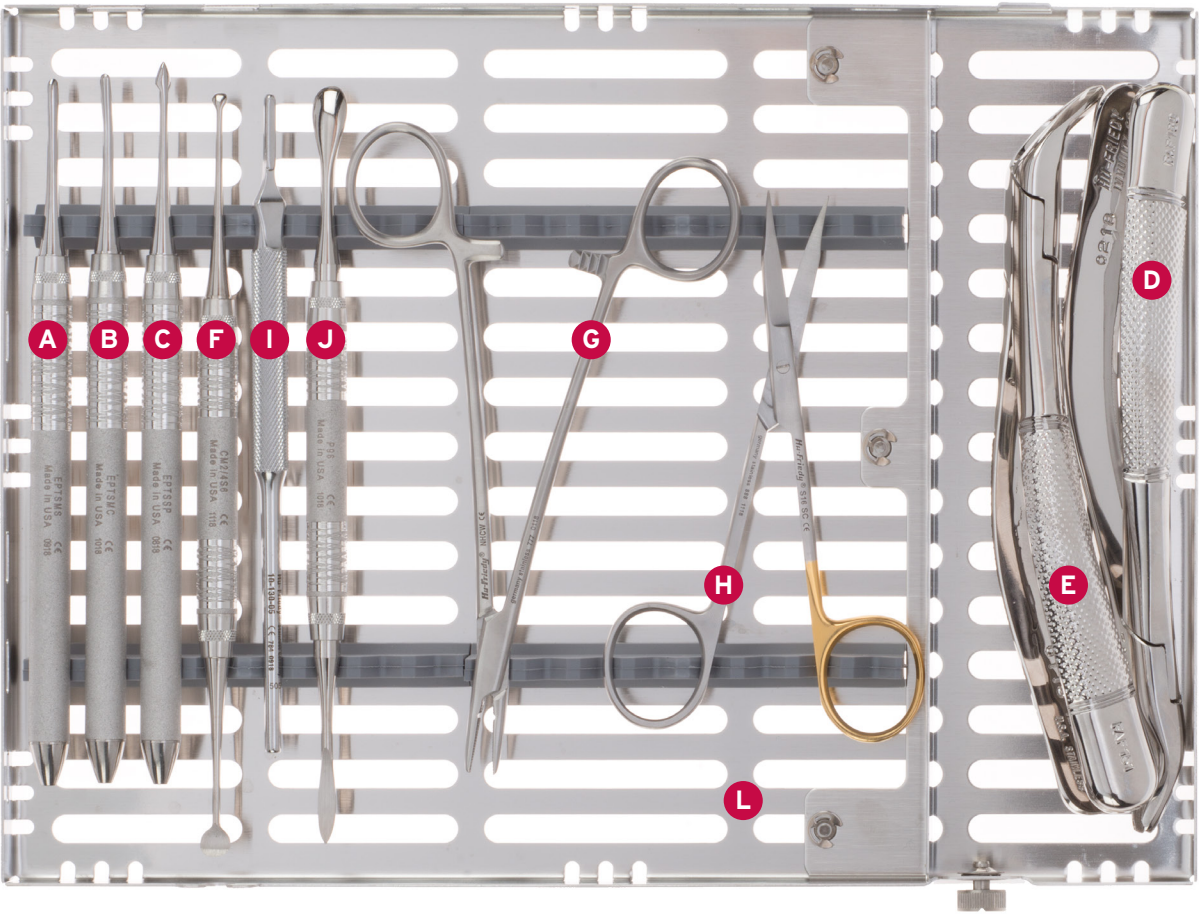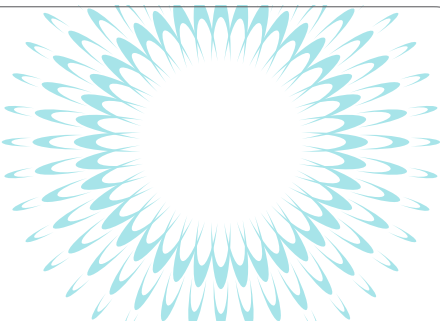

ATRAUMATIC EXTRACTION KIT



- A EPTSMS**
Luxating Hybrid Elevator, Small Straight
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- B EPTSMC**
Luxating Hybrid, Small Curved
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSSP**
Luxating Hybrid Elevator, Straight Spade
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D FAF150**
Apical Forceps 150
Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- E FAF151**
Apical Forceps 151
Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- F CM2/4S6**
2/4 Molt DE Curette, #6 Handle
Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- G NHCW**
Crile-Wood Needle Holder
Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- H S16SC**
16 Goldman-Fox Scissors, Super Cut
Used for trimming tissue or cutting sutures. Featuring tapered blade tips, tips that are not sharp on outside and are used as tonotomy scissors. 12.5cm with 1 serrated blade



- I 10-130-05**
5 Scalpel Handle, Round Straight
Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- J P96**
9 Molt Periosteal Elevator
Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm

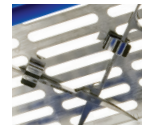

- K IM1002**
6-Pack Signature Series® Instrument Clip
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- L IM4161**
16 Instrument Signature Series® Cassette, Gray
IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.