



HU-FRIEDY SURGICAL KIT

BASIC PERIO KIT



PART CODE	DESCRIPTION
IM4121	12 Instrument Signature Series® Cassette, Gray
IM1002	6-Pack Signature Series® Instrument Clip
MIR5	5 FS CS Mouth Mirror
MH6	6 Cone Socket Mirror Handle
10-130-03	3 Scalpel Handle
KO12KPO3AR	1/2 Modified Allen Orban Knife Round
XP23/QOW6	23/QOW Expro, Wms Mkgs, #6 Handle
SG11/126	11/12 Gracey Curette, #6 Handle
SG13/146	13/14 Gracey Curette, #6 Handle
PPR36	Prichard Periosteal, #6 Handle
CM2/4S6	2/4 Molt DE Curette, #6 Handle
S13K/TG6	13/TG DE Curette/Chisel, #6 Handle
NHCW	Crile-Wood Needle Holder
S16	16 Goldman-Fox Scissors
IMBASPER	Complete Basic Perio Kit

BASIC PERIO KIT

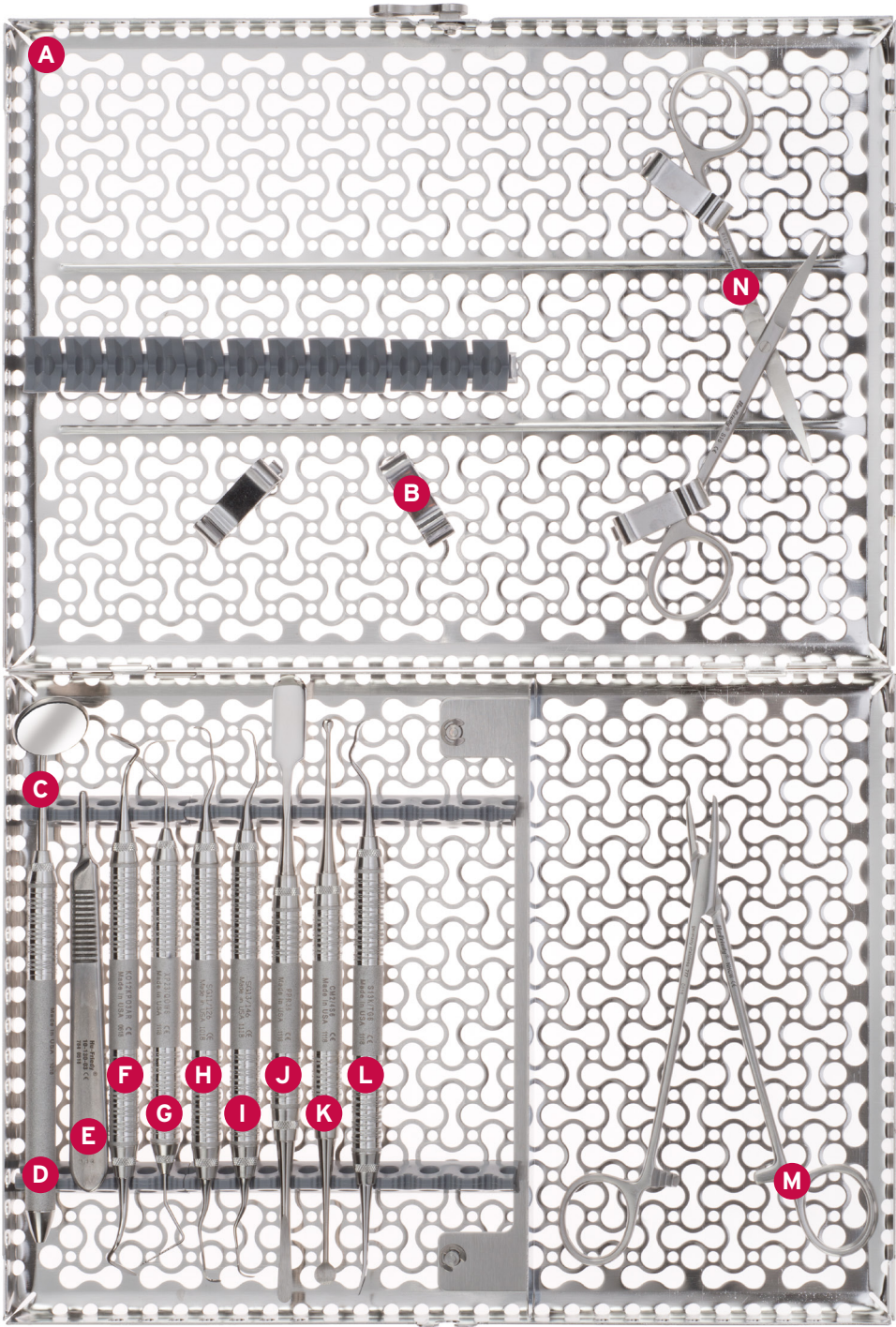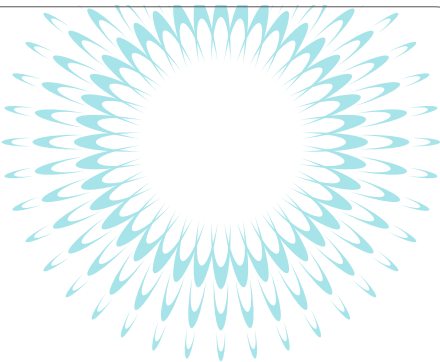

* pictured in IMN8121 - 12 Instrument Infinity Series™ Cassette, Gray

BASIC PERIO KIT



A IM4121
12 Instrument Signature Series® Cassette, Gray

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B IM1002
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



C MIR5
5 FS CS Mouth Mirror

Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



D MH6
6 Cone Socket Mirror Handle

Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



E 10-130-03
3 Scalpel Handle

Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



F KO12KPO3AR
1/2 Modified Allen Orban Knife Round

30% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



G XP23/QOW6
23/QOW Expro, Wms Mkgs, #6 Handle

A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



H SG11/126
11/12 Gracey Curette, #6 Handle

Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



I SG13/146
13/14 Gracey Curette, #6 Handle

Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



J PPR36
Prichard Periosteal, #6 Handle

For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



K CM2/4S6
2/4 Molt DE Curette, #6 Handle

For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



L S13K/TG6
13/TG DE Curette/Chisel, #6 Handle

Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with 3.3mm and 3.2mm working ends.



M NHCW
Crile-Wood Needle Holder

Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



N S16
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.