


HU-FRIEDY SURGICAL KIT

ATRAUMATIC EXTRACTION KIT, BLACK LINE




PART CODE	DESCRIPTION
EPTSMSX	Luxating Hybrid Elevator, Small Straight, Black Line
EPTSMCX	Luxating Hybrid, Small Curved, Black Line
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
FAF35XS	Extracting Forceps, 35N
FAF36XS	Extracting Forceps, 36N
CM2/4X	2/4 Molt DE Curette, Black Line
NHCW	Crile-Wood Needle Holder
S16SCX	Goldman-Fox, Curv, SC Scissors Black Line
10-130-05E	Scalpel Handle #5
P9X	9 Molt Periosteal Elevator, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMN4200	20 Instrument Infinity Series™ Cassette, Red
IMNATREXTX	Atraumatic Extraction Kit, Black Line

ATRAUMATIC EXTRACTION KIT, BLACK LINE

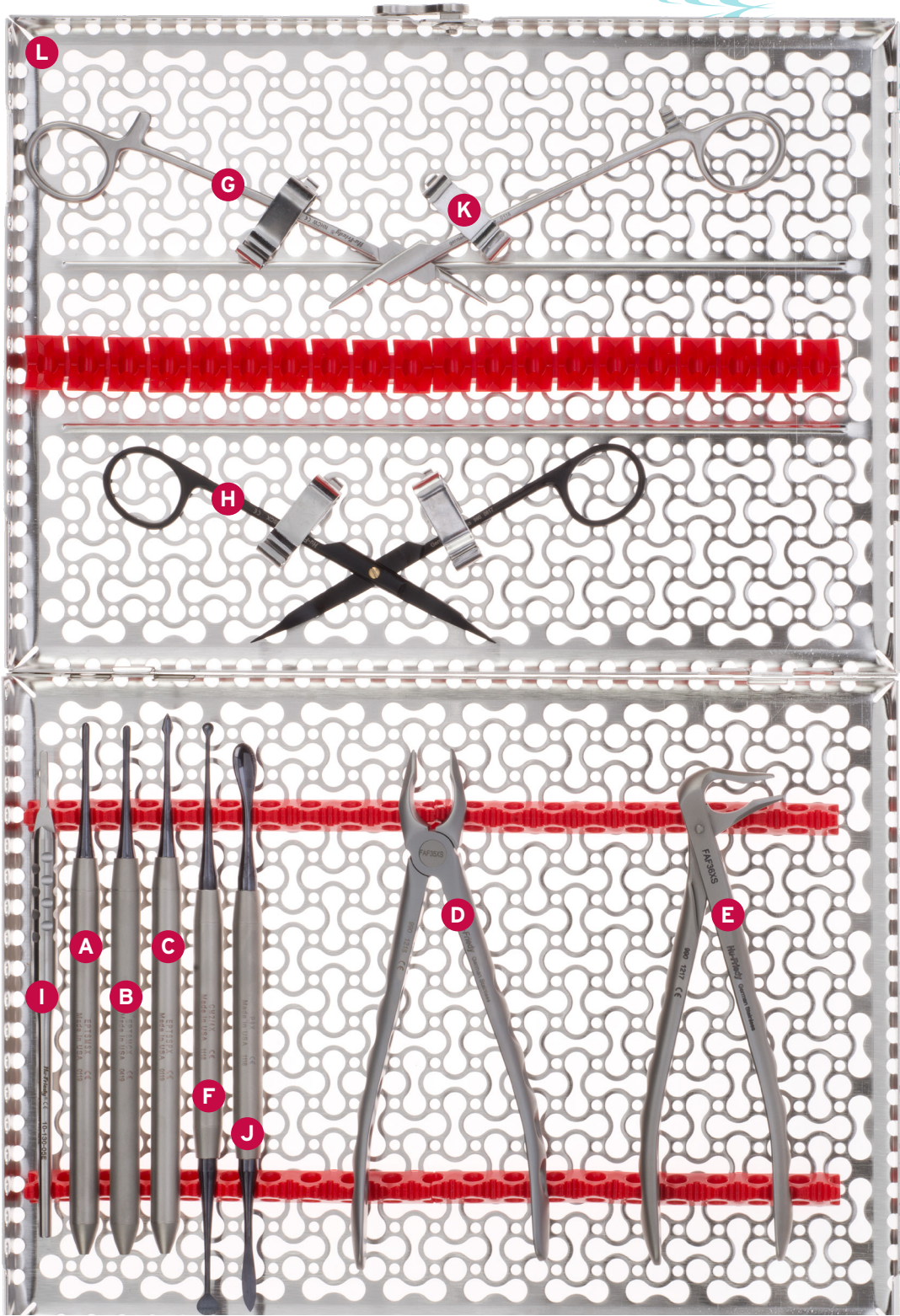

ATRAUMATIC EXTRACTION KIT, BLACK LINE

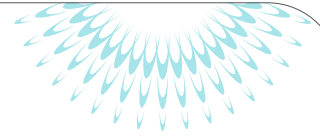

- A EPTSMSX**
Luxating Hybrid Elevator, Small Straight, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- B EPTSMCX**
Luxating Hybrid, Small Curved, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- C EPTSSPX**
Luxating Hybrid Elevator, Straight Spade, Black Line
A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- D FAF35X**
Extracting Forceps, 35N
Atraumair Forceps used to extract upper canines & premolars.



- E FAF36X**
Extracting Forceps, 36N
Atraumair Forceps used to extract lower incisors, canines and premolars.



- F CM2/4X**
2/4 Molt DE Curette, Black Line
Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- G NHCW**
Crile-Wood Needle Holder
Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- H S16SCX**
16 Goldman-Fox Scissors, Black Line
Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- I 10-130-05E**
5 Scalpel Handle
Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- J P9X**
9 Molt Periosteal Elevator, Black Line
Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm

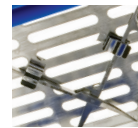

- K IM1002**
6-Pack Signature Series® Instrument Clip
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- L IMN4200**
20 Instrument Infinity Series™ Cassette, Red
IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.

