



HU-FRIEDY SURGICAL KIT

BASIC PERIO KIT, BLACK LINE



PART CODE	DESCRIPTION
IM1002	6-Pack Signature Series® Instrument Clip
MIR5	5 FS CS Mouth Mirror
MH6	6 Cone Socket Mirror Handle
10-130-03	3 Scalpel Handle
KO12KP03RX	1/2 Allen Orban Knife Round, Black Line
XP23/QOW6	23/QOW Expro, Wms Mkgs, #6 Handle
SG11/126	11/12 Gracey Curette, #6 Handle
SG13/146	13/14 Gracey Curette, #6 Handle
PPR3X	PR-3 Prichard Periosteal, Black Line
CM2/4X	2/4 Molt DE Curette, Black Line
S13K/13KLX	13/TG DE Curette/Chisel, Black Line
NHCW	Crile-Wood Needle Holder
IMN4121	12 Instrument Infinity Series™ Cassette, Gray
S16SCX	16 Goldman-Fox Scissors, Black Line
IMNBASPERX	Complete Basic Perio Kit, Black Line

BASIC PERIO KIT, BLACK LINE

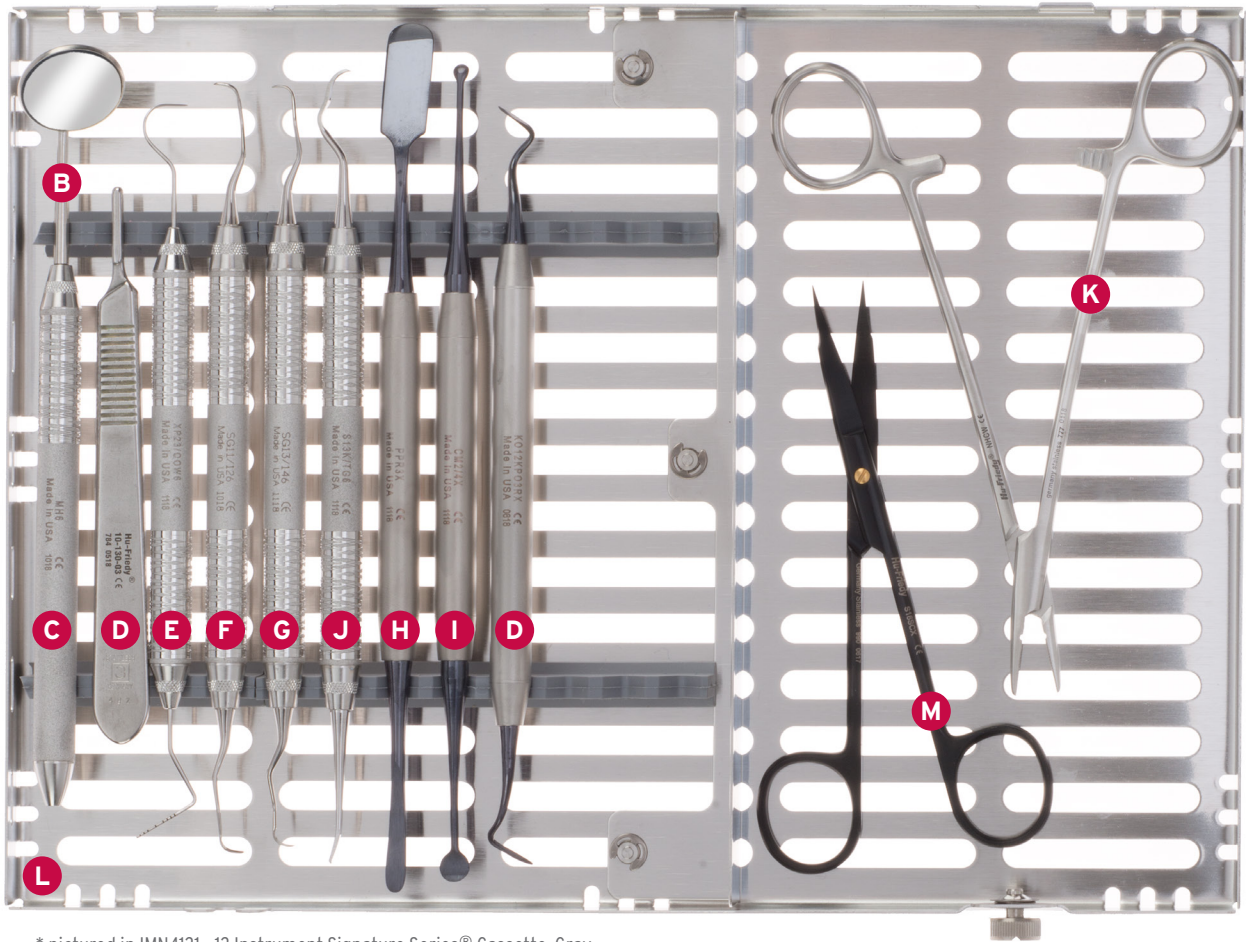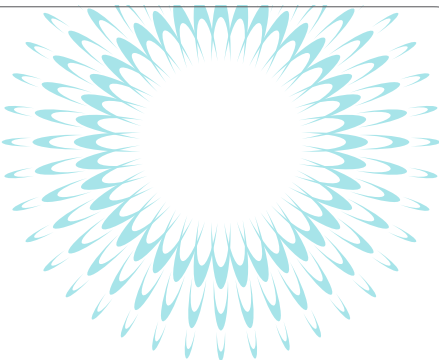

* pictured in IMN4121 - 12 Instrument Signature Series® Cassette, Gray

BASIC PERIO KIT, BLACK LINE



- A IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- B MIR5**
5 FS CS Mouth Mirror
 Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- C MH6**
6 Cone Socket Mirror Handle
 Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



- D 10-130-03**
3 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- E KO12KPO3RX**
1/2 Allen Orban Knife Round, Black Line
 20% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



- F XP23/QOW6**
23/QOW Expro, Wms Mkgs, #6 Handle
 A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- G SG11/126**
11/12 Gracey Curette, #6 Handle
 Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



- H SG13/146**
13/14 Gracey Curette, #6 Handle
 Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



- I PPR3X**
PR-3 Prichard Periosteal, Black Line
 Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



- J CM2/4X**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- K S13K/13KLX**
13/TG DE Curette/Chisel, Black Line
 Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with curved working ends 3.2mm and 4.9mm wide.



- L NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



- M IMN4121**
12 Instrument Infinity Series™ Cassette, Gray
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- N S16SCX**
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.