



HU-FRIEDY SURGICAL KIT

CRESTAL SINUS LIFT KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4162	16 Instrument Infinity Series™ Cassette, White
10-130-03	3 Scalpel Handle
P9X	9 Molt Periosteal Elevator, Black Line
CM2/4X	2/4 Molt DE Curette, Black Line
TP5043	Adson-Brown Perma Sharp® Tissue Pliers
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
PCP65006	20mm Probe 8-10-12-14-16-18-20, #6 Handle
OSTPAL1	1 Palti Osteotome 2.7mm
OSTPAL2	2 Palti Osteotome 3.2mm
OSTPAL3	3 Palti Osteotome 3.7mm
OSTPAL4	4 Palti Osteotome 4.2mm
PFIWDS16	1 Woodson Plastic Filling, #6 Handle
IM1002	6-Pack Signature Series® Instrument Clip
NH5038	Crile-Wood Perma Sharp® Needle Holder
IMNCREKITX	Complete Crestal Sinus Lift Kit, Black Line

CRESTAL SINUS LIFT KIT, BLACK LINE

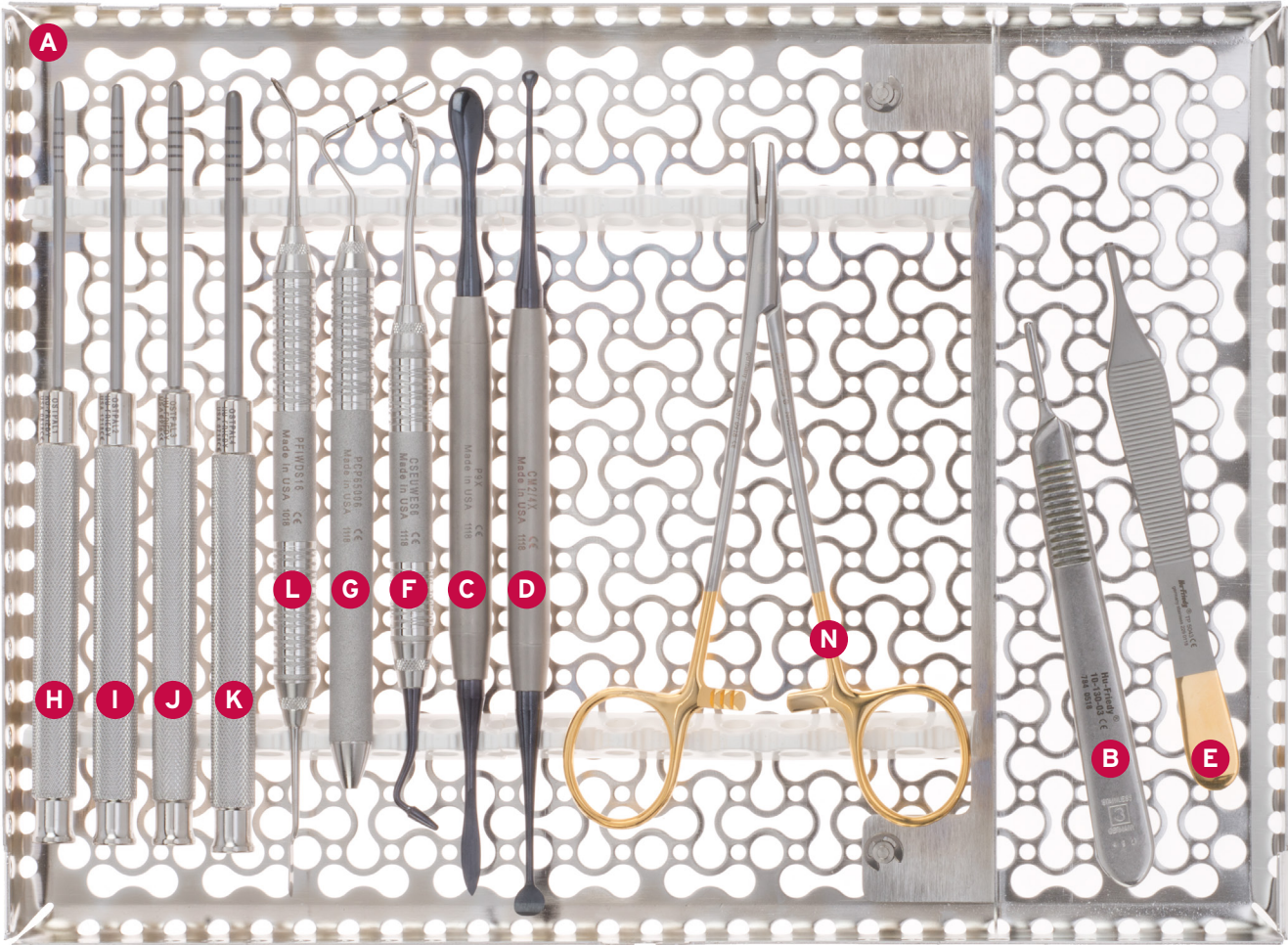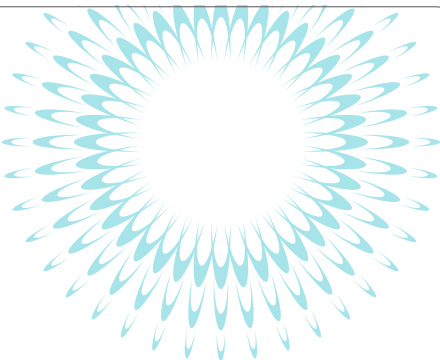

CRESTAL SINUS LIFT KIT, BLACK LINE



- A IMN4162**
16 Instrument Infinity Series™ Cassette, White
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- B 10-130-03**
3 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- C P9X**
9 Molt Periosteal Elevator, Black Line
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



- D CM2/4X**
2/4 Molt DE Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- E TP5043**
Adson-Brown Perma Sharp® Tissue Pliers
 Used to grasp and stabilize soft tissue during suturing. 12cm



- F CSEUWES6**
Euwe Bone Scoop/Plugger, #6 Handle
 Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



- G PCP65006**
20mm Probe 8-10-12-14-16-18-20, #6 Handle
 A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- H OSTPAL1**
1 Palti Osteotome 2.7mm
 Used for preparing bone. 7-10-13-15-18mm



- I OSTPAL2**
2 Palti Osteotome 3.2mm
 Used for preparing bone. 7-10-13-15-18mm



- J OSTPAL3**
3 Palti Osteotome 3.7mm
 Used for preparing bone. 7-10-13-15-18mm



- K OSTPAL4**
4 Palti Osteotome 4.2mm
 Used for preparing bone. 7-10-13-15-18mm



- L PFIWDS16**
1 Woodson Plastic Filling, #6 Handle
 Designed to reflect and retract the mucoperiosteum after incisions of the gingival tissue with 2.2mm and 3.8mm width.



- M IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- N NH5038**
Crile-Wood Perma Sharp® Needle Holder
 Used to grasp and guide the needle when suturing. Featuring 3-0 to 5-0 carbide beaks, finger ring style.