

PINK MATTERS® X HUFRIEDYGROUP SURGICAL KIT

Minimally Invasive Periodontal and Implant Surgery Solution




Pink Matters® was initiated in Europe by a group of clinicians in order to improve both the patient and practitioner's quality of life. "The Pink Matters Group¹ believes in the synergy of skills. Because of this, our selective range of products are designed by practitioners and developed in collaboration with the most recognized and innovative manufacturers and researchers. The Pink Matters core value is based on ethics and high quality standards through science, practice and knowledge.

Discover more on
Pink Matters Group website
www.pink-matters.com

“ Minimally invasive periodontal and implant surgery has become the new standard to improve healing processes and patient quality of life. With the development of many innovative surgical procedures the need for specifically designed instruments has become crucial. Pink Matters®, a group of dedicated clinicians, partnered with HuFriedyGroup, the world leader of high quality surgical instruments, to develop a universal surgical kit, including five new instruments. It allows clinicians to perform all modern periodontal and implant surgery procedures. ”

Pink Matters Group

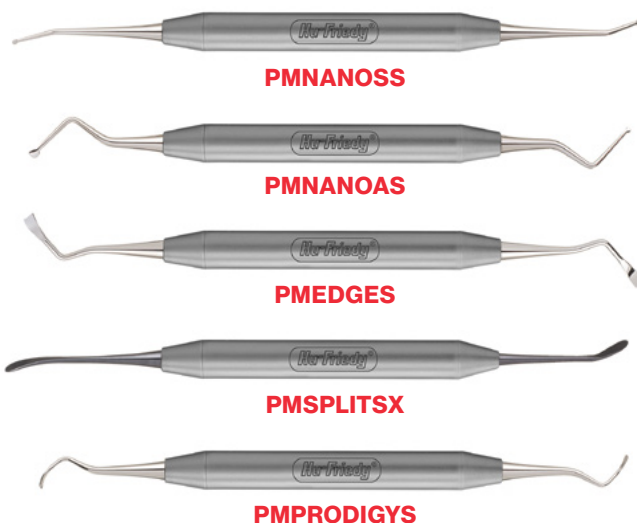

VERSATILITY AND EFFICIENCY: PINK MATTERS KEYS FOR SUCCESS

The Pink Matters® X HuFriedyGroup Surgical Kit is the result of the partnership between the Pink Matters group and HuFriedyGroup. The kit includes five highly versatile instruments designed for minimally invasive periodontal and implant surgery.

In the kit, there are five instruments: two periodontal knives (which can also be used as micro-elevators), a chisel, and an innovative periodontal curette. Instruments benefit from EverEdge™ 2.0 Technology which utilizes optimized heat control and enhanced finishing processes in order to create sharp edges needed for effective procedures and continued use. The periosteal is offered with Black Line coating. This dark, hard coating provides enriched contrast at the surgical site and enhanced lubricity.

THE PINK MATTERS® INSTRUMENTS




WATCH OUR
TUTORIAL VIDEOS
FOR THE PINK MATTERS®
INSTRUMENTS



PMNANOSS

PM Nanos S (Straight), EverEdge™ 2.0 Perio Knife, Short

Used for supra periosteal dissection apically and laterally during tunneling and soft tissue surgical procedures. This instrument can also be used to elevate thin papilla during partial or full thickness flap surgery and is convenient for papillary preservation technique.

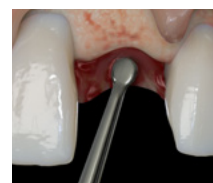

PMNANOAS

PM Nanos A (Angulated), EverEdge™ 2.0 Perio Knife, Short

Used for supra periosteal dissection apically and laterally during tunneling and soft tissue surgical procedures. This instrument can also help in the creation of a partial or full thickness envelope flap during connective tissue grafting procedures in perio plastic surgery and implant surgery. Nano A is also useful when elevating a full thickness flap in the lingual and palatal area during flap surgery.



PMEDGES

PM Edge, EverEdge™ 2.0 Crestal Bone Chisel, Short

Designed to elevate full thickness flap in the posterior area in preparation for implant placement. The Edge can be also very useful early on in the implant placement procedure when reflecting the thicker buccal soft tissue. This instrument allows for bone recontouring during crown lengthening and is very efficient at bone harvesting in areas with delicate anatomy, such as the nasal spine.



PMSPLITSX

PM Split, Black Line Periosteal, Short


Two sides thoughtfully developed: one side is sharp to safely elevate full thickness flaps in perio and implant surgery. The other side is blunt to reduce the surgical risk when separating elastic fibers during buccal and lingual flap advancement.



PMPRODIGYS

PM Prodigy, EverEdge™ 2.0 Surgical Curette, Short

Developed to improve granulation tissue removal during flap surgery. The angulation of the curette improves the access to narrow periodontal and peri-implant defects, particularly in the posterior area.

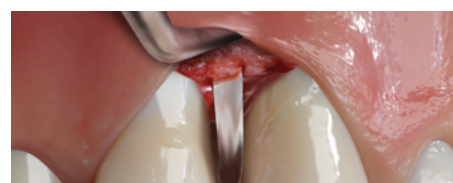

ADDITIONAL KIT INSTRUMENTS



S204SCCHE2

S204S
Matte finish



PCPUNC15CH

Probe UNC
#UNC15 CH
Hdl Qulix, 1-15



10-130-5EM

European Round
Scalpel Handle
matte finish



MIR5DSHDB

Double-sided Mirror
Black Line

MH6B

Mirror Handle
Black Line

REDUCES GLARE UP TO 80%²

See your patients – not your instruments. Glare interferes with visibility. Reducing this disruption to your visibility allows you to see more detail and offer more accurate patient care. HD Black Line Mirrors are engineered to optimize clinical outcomes by delivering superior visibility throughout any dental procedure.


SPV

Curved Microsurgical
Scissor



S18SC

Iris Curved Super-Cut™
Scissor

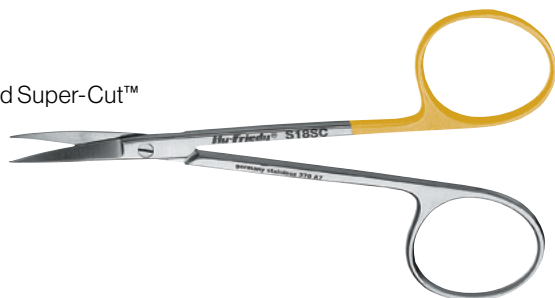

CRMX

Minnesota Retractor,
Black Line




FEATURES AND BENEFITS OF BLACK LINE:

- Performance engineered coating for a harder, smoother surface for optimal edge retention and enhanced lubricity
- Distinct black finish for enriched contrast and visual acuity at the surgical site and underlying tissue
- Reduced light reflection afforded by a matte finished handle and black working ends
- Unique, smooth, large diameter lightweight handle for maximum comfort, reduced hand fatigue and increased control
- Handcrafted with Immunity Steel™ Stainless Steel Alloy for optimal strength and corrosion resistance



TPSLCOSM

Multi Tooth
Straight
Stainless Steel
Micro Tissue Plier

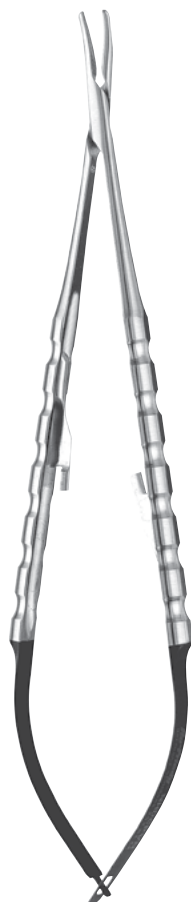

NHDCPVN

Micro Castroviejo
Curved Diamond
Dusted Needle
Holder



TP5090

DeBakey
Perma Sharp™
Tissue Plier



SHDPV

Microsurgical
Scalpel Handle

WHY YOU SHOULD USE AN IMS™ CASSETTE SYSTEM:

By standardizing your instruments with IMS™ System, additional time and revenue is generated while extending the useful life of your HuFriedyGroup Instruments.

Dental professionals state that by implementing IMS™ System, they save 5 to 10 minutes per procedure allowing the staff to spend more time on value-added activities.³

In addition, safer transportation of soiled instruments reduces risk of injury and exposure to blood borne pathogens.



IMN4167

Large Cassette with Accessory Area, Black


IM1000

2-Pack Hinged Instrument Clips



Learn more about
Pink Matters® Surgical Kit

ESSENTIAL SURGICAL KIT

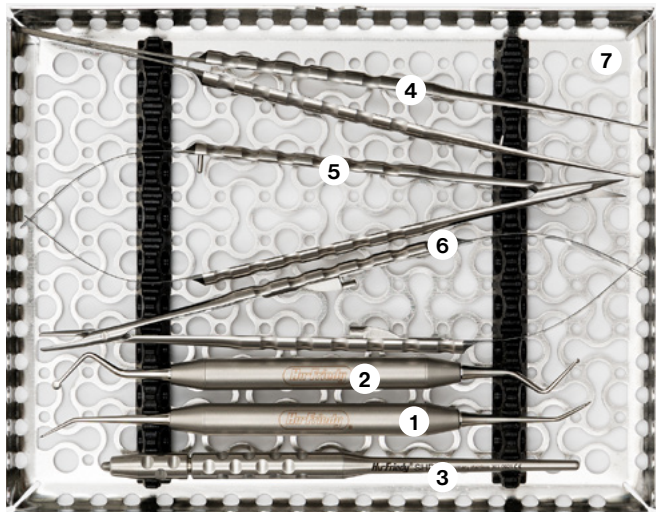

IMN5107

IMS™ Infinity Series™ 10 Instrument Cassette, Black

Dimensions: 152 x 184 x 34 mm



PMNANOSS

PM Nanos S (Straight),
EverEdge™ 2.0
Perio Knife, Short



PMNANOAS

PM Nanos A
(Angulated),
EverEdge™ 2.0
Perio Knife, Short



SHDPV

Microsurgical
Scalpel Handle




TPSLCOSM

Multi Tooth Straight
Stainless Steel
Micro Tissue Plier



SPV

Castro
Microsurgical

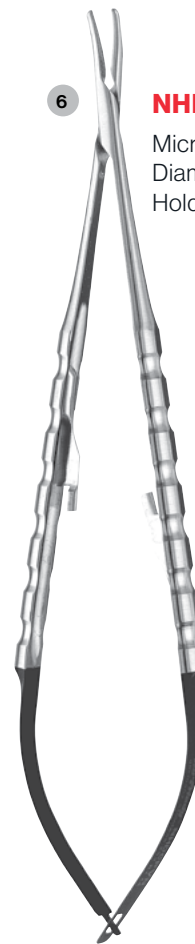

NHDCPVN

Micro Castroviejo Curved
Diamond Dusted Needle
Holder

1) Pink Matters® is not owned by or affiliated with HuFriedyGroup.

2) When comparing the HuFriedyGroup DLC coated mirror head and handle to the HuFriedyGroup non-coated stainless-steel mirror head and handle.

3) When compared to single instrument decontamination. Based on market survey results. Data on file. Available upon request.