
Self-Seal Nylon Sterilization Pouches

REF NSP-400 / NSP-405 / NSP-410 / NSP-420 / NSP-430

INSTRUCTIONS FOR USE

1. Choose the appropriate pouch size for the items to be sterilized.
2. Carefully place instruments into the pouch and remove any excess air.
3. Utilize internal and external dry heat chemical indicators and biological indicators in accordance with facility guidelines and local standards. External dry heat indicators are supplied with this product for your convenience. When used apply one indicator to the outside of the pouch, avoiding placement at the seal. Ensure indicator is facing up during sterilization. Indicators color transition from green to dark when exposed to dry heat sterilization.
4. Remove the paper from the adhesive strip, fold evenly along perforation and secure firmly. Ensure the seal is flat and free of wrinkles or creases.
5. Ensure the pouch contents are a minimum of 1" (2.5 cm) away from all sides of the pouch.
6. Follow directions for sterilization according to the manufacturer of the sterilizer.
7. Allow pouches to cool in the sterilizer 5-10 minutes after cycle completion.



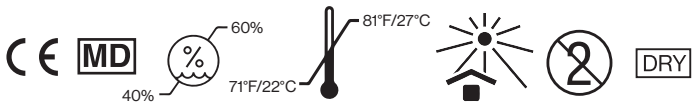

EN ISO 11140-1 Type 1

CROSSTEX

Crosstex International, Inc.
6789 W. Henrietta Road, Rush, NY 14543 USA
585.359.0130 | 631.582.6777 | hufriedygroup.com

EU REP MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany

UK REP MDSS-UK RP Ltd.
Parkway House, Palatine Rd,
Northenden, Wythenshawe,
Manchester M22 4DB
United Kingdom

CH REP MDSS CH GmbH
Laurenzenvorstadt 61
5000 Aarau
Switzerland

Made in USA

NSP-IFU / REV 1 - 04/2026

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted.
©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved.

CROSSTEX