

Surgical Kit

Keys Socket Grafting Kit (IMNKEYSGKT)



1 IMN4164
16 Instrument Infinity Series™ Cassette, Purple

IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



2 XP23/86
23/CP-8 Expro

Double-ended instrument with an explorer on one end and a probe (color-coded) on the other. ergonomic, Satin Steel™ Handle.



3 H2
2 Kelly Hemostat, Curved

Used to clamp off blood vessels, remove small root tips and grasp loose objects. Featuring kelly curved beaks.



4 PLGGR1X
Grisdale Bone Condenser 8-10-13-15-18, Black Line

Used during grafting procedures with 8-11-13-15-18 mm markings. Featuring a non-reflective surface, improved contrast between tissue/scissor.



5 CM2/4SX
2/4 Molt DE Curette, Black Line

Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



6 10-130-05
5 Scalpel Handle, Round, Straight

Used for cutting gingival tissue and making surgical incisions.



7 PPR3X
PR-3 Prichard Periosteal, Black Line

Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



8 NH5037
Halsey Perma Sharp™ Needle Holder, Serrated

Used to grasp and guide the needle when suturing. Featuring carbide inserts, finger ring style. Use on all types of sutures.



9 S5055
Metzenbaum Perma Sharp™ Scissors, Curved/Blunt

Used for trimming tissue or cutting sutures. 14.5cm



10 TP5090
DeBakery Perma Sharp™ Tissue Pliers, 6"

Used to grasp and stabilize soft tissue during suturing.



11 DP1
#1 Dressing Pliers, Perry

Used to grip or transfer materials in and out of the oral cavity. Features a serrated tip to provides a positive grip.



12 PMPIX
Membrane Placement Instrument, Black Line

Micropoint tip for precise placement and manipulation of membranes.



13 MCUPE
Measuring Cup 25 mL w/ Marking

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



14 SYBG
Syringe Bone Grafting 7mm

Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.




15 CL856
85 Lucas DE Surgical Curette

Designed for curettage, cyst removal and tooth socket debridement. 2.5mm



16 PJO/16
S13K/TG Jovanovic 4mm Chisel


For removing and recontouring bone.

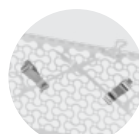

17 IM1002
6-Pack Hinged Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

Keys Socket Grafting Kit (IMNKEYSGKT)

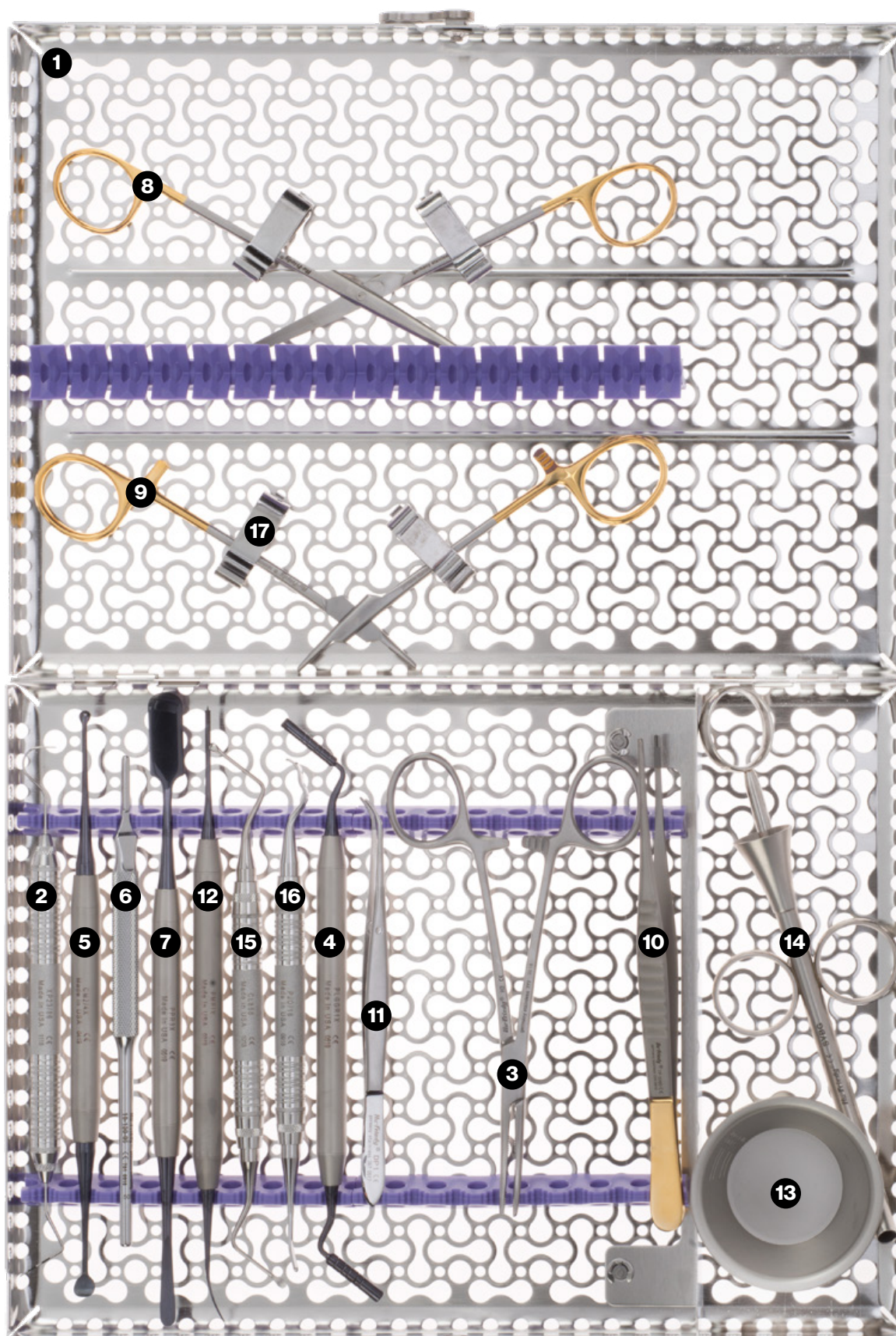

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.

©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-347/1023