

Surgical Kit

Basic Perio Kit, Black Line (IMNBASPERX)



- 1** **IM4121**
12 Instrument Infinity Series™
Cassette, Gray
IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2** **MIR5**
5 FS CS Mouth Mirror
Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- 3** **MH6**
6 Cone Socket Mirror Handle
Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



- 4** **10-130-03**
3 Scalpel Handle
Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- 5** **XP23/QOW6**
23/QOW Expro, Wms Mkgs, #6 Handle
A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- 6** **SG11/126**
11/12 Gracey Curette, #6 Handle
Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



- 7** **SG13/146**
13/14 Gracey Curette, #6 Handle
Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



- 8** **S13K/13KLX**
13/TG DE Curette/Chisel, Black Line
Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with curved working ends 3.2mm and 4.9mm wide.



- 9** **PPR3X**
PR-3 Prichard Periosteal, Black Line
Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



- 10** **CM2/4X**
2/4 Molt DE Curette, Black Line
Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.




- 11** **KO12KPO3RX**
1/2 Allen Orban Knife Round, Black Line
20% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



- 12** **NHCW**
Crile-Wood Needle Holder
Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- 13** **S16SCX**
16 Goldman-Fox Scissors, Black Line
Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.

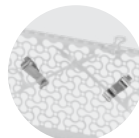

- IM1002**
6-Pack Hinged Instrument Clip
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

Basic Perio Kit, Black Line (IMNBASPERX)

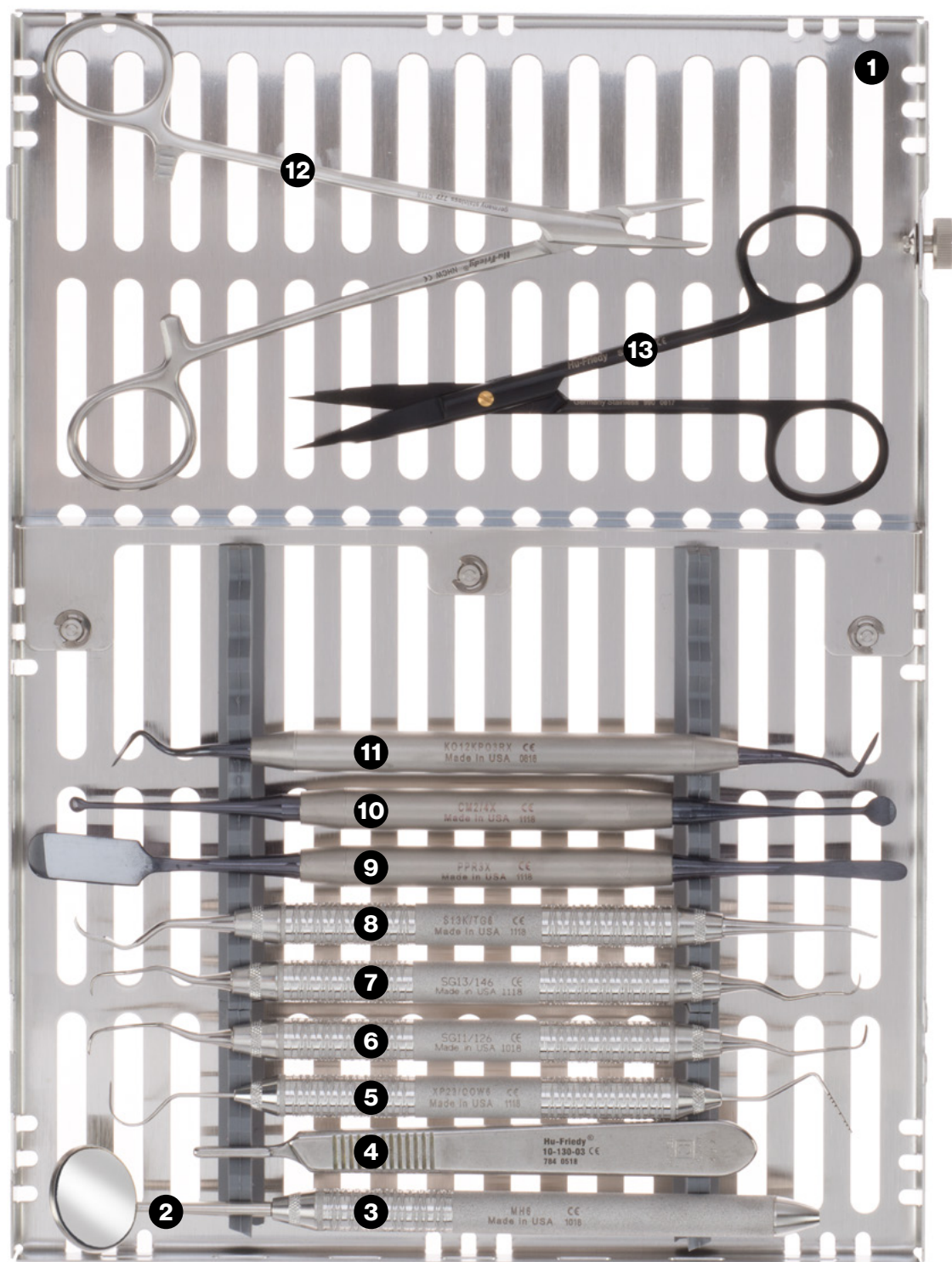

Pictured in IMN4121 - 12 Instrument Signature Series™ Cassette, Gray

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.
©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-343/1023