

Surgical Kit

Extraction Kit (IMEXTKIT)



- 1 IM4163E2C4**
16 Instrument Signature Series™
Cassette, Orange
 IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2 10-130-05**
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- 3 P96**
9 Molt Periosteal Elevator
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



- 4 NHMH**
Mayo-Hegar Needle Holder
 Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



- 5 S16**
16 Goldman-Fox Scissors
 Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.




- 6 CL866**
86 Lucas DE Surgical Curette, #6 Handle
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



- 7 E301**
301 Apexo Elevator
 Used to loosen the tooth from the periodontal ligament and ease extraction.



- 8 E34S**
34S Seldin Straight Elevator
 Used to loosen the tooth from the periodontal ligament and ease extraction. 4.0mm

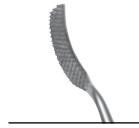

- 9 BF2X6**
2X Miller-Colburn Bone File
 Used for final smoothing and shaping of bone with 6.2mm wide working ends.



- 10 EHB13/146**
13/14 Heidbrink DE Root Pick
 Used for the removal of root tips in difficult-to-reach areas. Designed for delicate root teasing procedures. Also available in a "heavy" pattern for more substantial root tip removal procedures.



- 11 RF**
Friedman Rongeurs
 Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



- 12 TP43**
43 Adson-Brown Tissue Pliers
 Used to grasp and stabilize soft tissue during suturing. 12mm




- 13 EPTSSP**
Luxating Hybrid Elevator, Straight Spade
 A unique cross between a periotome and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- 14 IM1002**
6-Pack Hinged Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

Extraction Kit (IMEXTKIT)

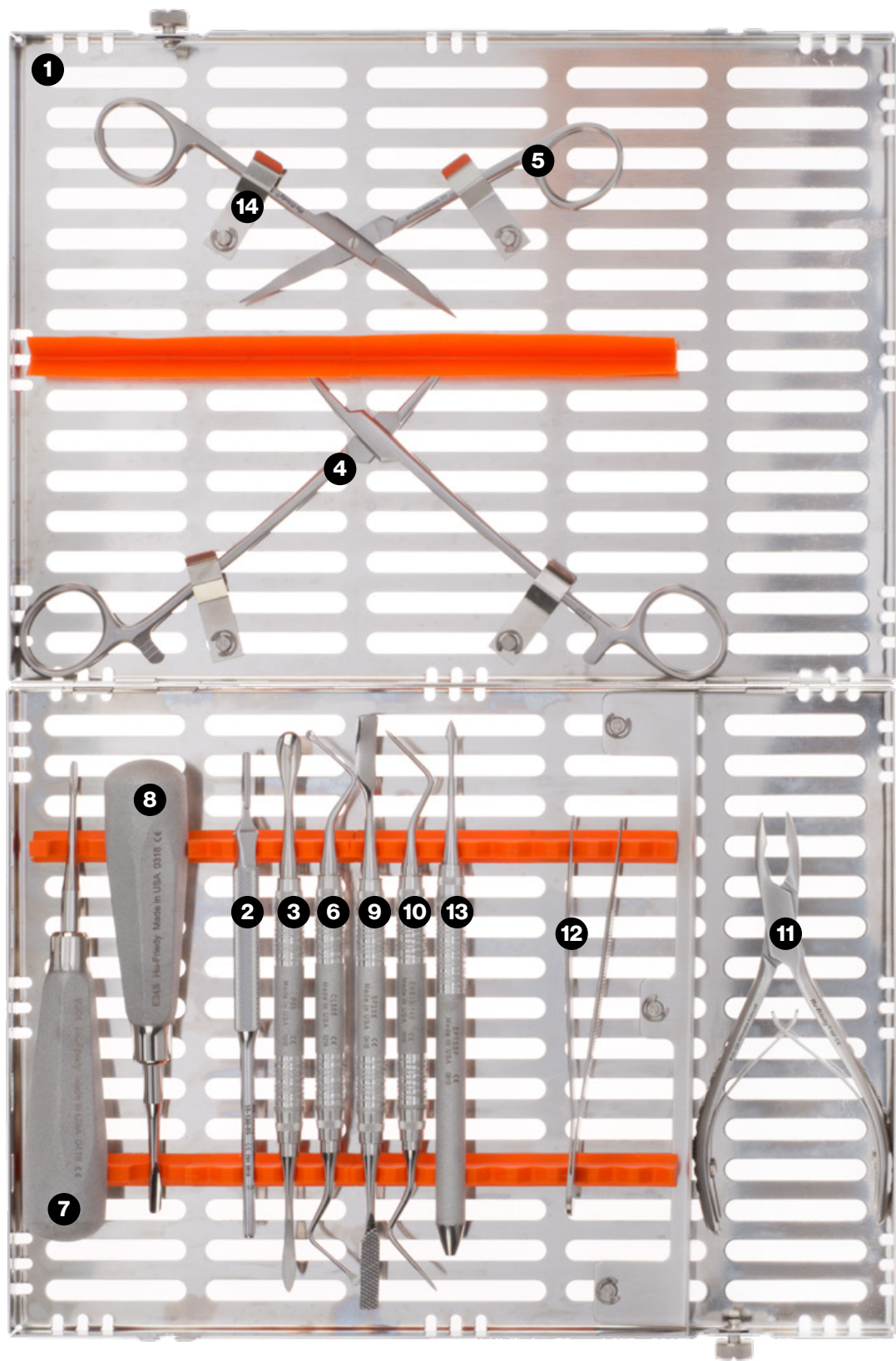

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.
©2026 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-339 /1023