

# Surgical Kit

## Basic Perio Kit (IMBASPER)



- 1** **IM8121**  
**12 Instrument Signature Series™ Cassette, Gray\***  
IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2** **MIR5**  
**5 FS CS Mouth Mirror**  
Single-sided Cone Socket mirror. 15/16 in. (24mm) diameter.



- 3** **MH6**  
**6 Cone Socket Mirror Handle**  
Used with Single and Double-Sided Cone Socket Mirrors. The blunt end of the mirror handle is ideal for testing sensitivity to percussion and assessing tooth mobility.



- 4** **10-130-03**  
**3 Scalpel Handle**  
Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



- 5** **KO12KPO3AR**  
**1/2 Modified Allen Orban Knife Round**  
30% thinner than the standard Orban, features a rounded angled blade and contra-angled shank for posterior use. Also useful for excising interproximal tissue.



- 6** **XP23/QOW6**  
**23/QOW Expro, Wms Mkgs, #6 Handle**  
A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



- 7** **SG11/126**  
**11/12 Gracey Curette, #6 Handle**  
Angulated to reach mesial surfaces of posterior teeth. Double-ended, easy-grip handle.



- 8** **SG13/146**  
**13/14 Gracey Curette, #6 Handle**  
Angulated to reach distal surfaces of posterior teeth. Double-ended, easy-grip handle.



- 9** **PPR36**  
**Prichard Periosteal, #6 Handle**  
For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



- 10** **CM2/4S6**  
**2/4 Molt DE Curette, #6 Handle**  
For reflecting and retracting the mucoperiosteum. Featuring two blades, one large flat blade and one narrow sharp blade on round tip. 4.7mm/11.0mm



- 11** **S13K/TG6**  
**13/TG DE Curette/Chisel, #6 Handle**  
Used to recontour periodontal bone or to remove remaining tooth fibers during periodontal surgery with 3.3mm and 3.2mm working ends.



- 12** **NHCW**  
**Crile-Wood Needle Holder**  
Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- 13** **S16**  
**16 Goldman-Fox Scissors**  
Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



- 14** **IM1002**  
**6-Pack Hinged Instrument Clip**  
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

\*Pictured in IMN8121 - 12 Instrument Infinity Series™ Cassette, Gray

## Basic Perio Kit (IMBASPER)

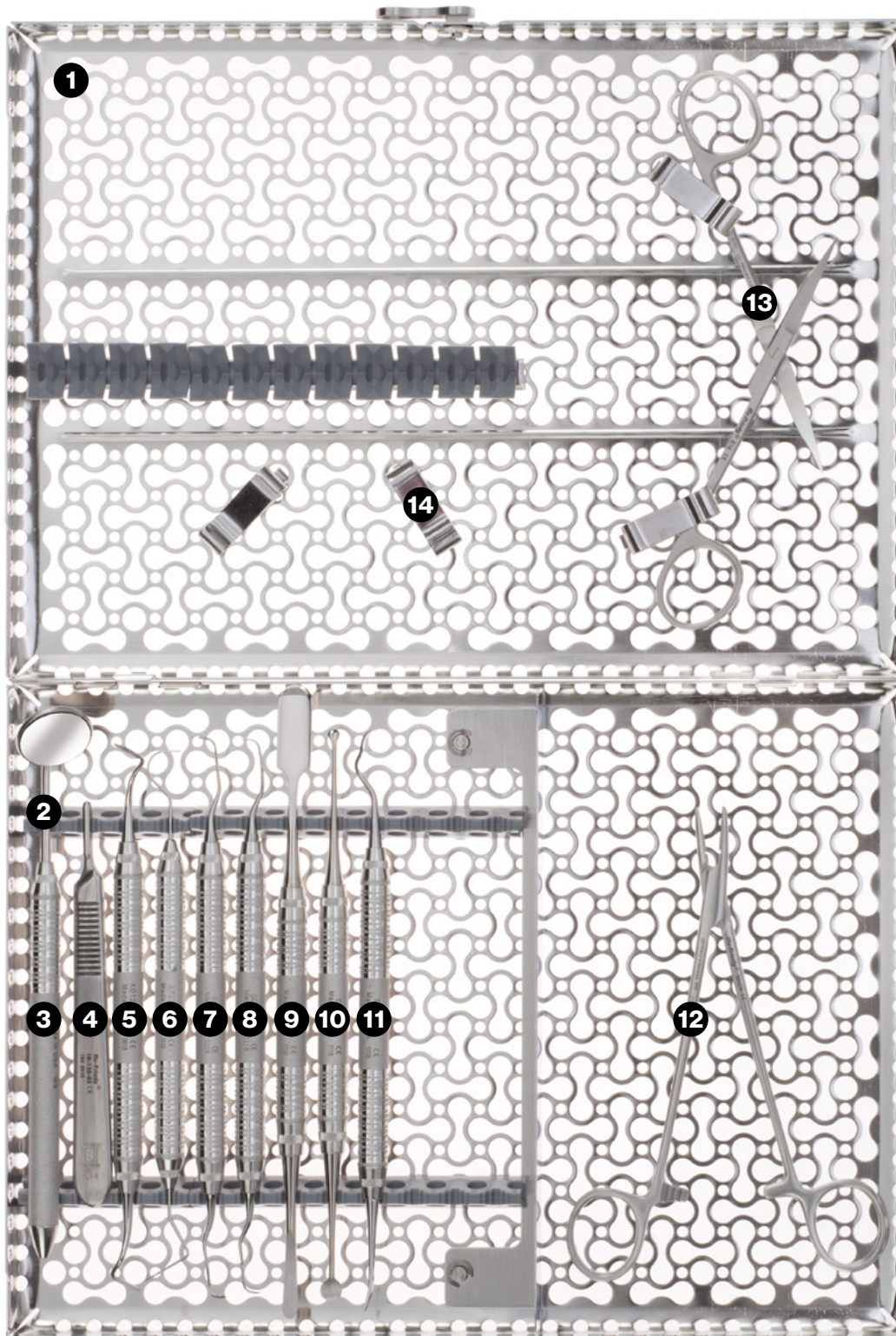

Pictured in IMN8121 - 12 Instrument Infinity Series™ Cassette, Gray

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | [HuFriedyGroup.com](http://HuFriedyGroup.com)

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.  
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-336/1023