



HU-FRIEDY PROCEDURAL SET-UP

RUBBER DAM CLAMP KIT LIST



PART CODE	DESCRIPTION	QUANTITY
RDCMW2	Clamp, Rubber Dam #W2, Satin Steel	
RDCMW3	Clamp, Rubber Dam #W3, Satin Steel	
RDCMW7	Clamp, Rubber Dam #W7, Satin Steel	
RDCM27	Clamp, Rubber Dam #27, Satin Steel	
RDCM212	Rubber Dam Clamp #212, Satin Steel	
RDPA	Rubber Dam Punch-Ainsworth	
RDF	Rubber Dam Forceps	
RDTFDI	Rubber Dam Template, FDI	
IM8128CE2L	MTO 12 IN SIG SERIES,CAS,BLUE,ECR,E2,8LP	
RDAF6	Rubber Dam Frame, Adult	
KIT-RDC2	Complete Rubber Dam Clamp Kit	

RUBBER DAM CLAMP KIT 2 LIST

A RDCMW2

B RDCMW3

C RDCMW7

D RDCM27

E RDCM212


F RDPA

G RDF

H RDTFDI

I IM8128CE2L

J RDAF6

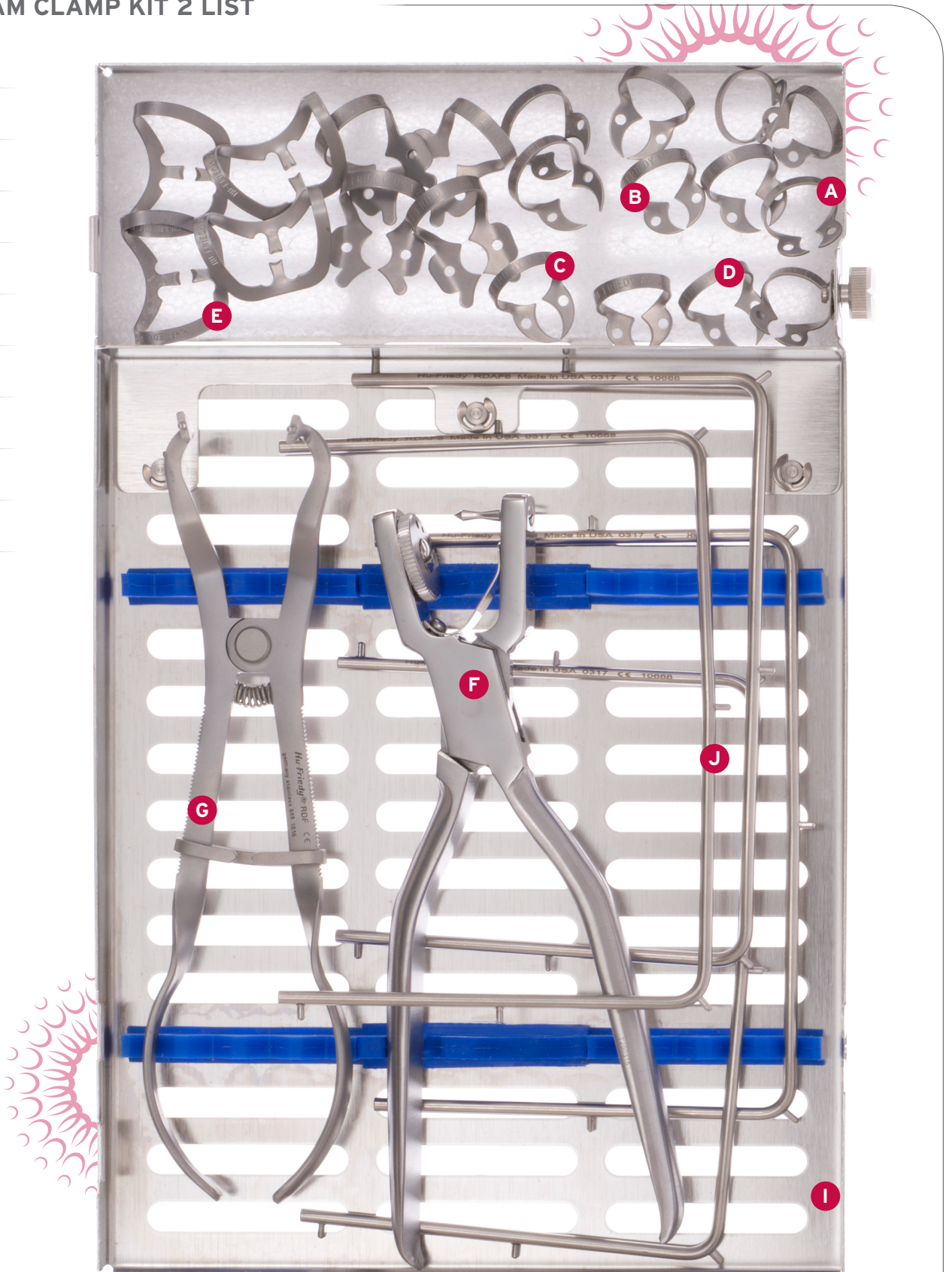

H Rubber Dam Template, FDI [not pictured]

RUBBER DAM CLAMP KIT 2 LIST



- A RDCMW2**
Clamp, Rubber Dam #W2, Satin Steel
Flat jaws for premolars.



- B RDCMW3**
Clamp, Rubber Dam #W3, Satin Steel
Flat narrow jaws for small lower molars.



- C RDCMW7**
Clamp, Rubber Dam #W7, Satin Steel
Flat jaws for lower molars.



- D RDCM27**
Clamp, Rubber Dam #27, Satin Steel
Flat jaws for large molars.




- E RDCM212**
Rubber Dam Clamp #212, Satin Steel
Flat jaws for labial caries on anterior teeth.



- F RDPA**
Rubber Dam Punch-Ainsworth
17 cm (6 2/4) Ainsworth Rubber Dam Punch, essential for efficient treatment. A popular rubber dam clamp accessory, this classic punch is solidly built to for lasting longevity. Hu-Friedy also has a wide variety of rubber dam clamps in winged, wingless and butterfly designs, all made with Satin Steel® to help reduce the glare caused by operator lights for less eye fatigue and improved clinical outcomes.



- G RDF**
Rubber Dam Forceps
17 cm (6 1/2") Rubber Dam Clamp Forceps essential for efficient treatment. A popular rubber dam clamp accessory, this instrument features a useful attachment design which makes it easy to pick up and place rubber dam clamps.

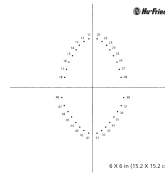

- H RDTFDI**
Rubber Dam Template, FDI

A useful way to determine exactly where to punch the rubber dam hole based on tooth location. Simply place the dental dam on the template, find the appropriate tooth location and create your punch hole. Available in Adult and Pedodontic.



- I IM8128CE2L**
MT0 12 IN SIG
SERIES,CAS,BLUE,ECR,E2,8LP

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- J RDAF6**
Rubber Dam Frame, Adult

One of two stainless steel rubber dam frames, this larger version has a 15cm/6" diameter. Made from high quality stainless steel, its dental dam attachment points provide the durable longevity associated with stainless steel rubber dam frames.