

PWR Pair

Piezoelectric Scaler & Air Polisher



Power
Equipment



 **HuFriedyGroup**
The Best In Practice

Introducing the **PWR Pair**

The power of all-in-one brings versatile air polishing to every clinician treating every patient. The PWR Pair brings you the power of supra- and subgingival calculus and stain removal, biofilm disruption, and implant cleaning.

Learn more about PWR Piezo Scaling & Air Polishing at
HuFriedyGroup.com/PWR-ON



PWR Pair
APPAIR

PIEZO + AIR = PWR PAIR

PWR PAIR PIEZO SCALER

Use the PWR Pair Piezo Scaler to remove calculus and biofilm. The advanced feedback system gives the clinician perfect control over performance to improve results and increase comfort.

PWR PAIR AIR POLISHING

Use PWR Pair Air Polisher to remove biofilm and tough stains supra- and subgingivally. PWR Air Polishing also helps remove trapped food debris and plaque build-up from between your patients' teeth.

What's in the Box

PWR Pair Piezo Electric Scaler & Air Polisher		APPAIR	
PWR Piezo LED handpiece	UPHPLED	PWR Air Powder Chamber Cap	APPCCAP
PWR Foot Pedal	UPFTPEDAL	PWR Air 120° Handpiece	APHP120
PWR Piezo Tip P3	UPP3	PWR Air Perio Handpiece	APHPPERIO
PWR Piezo Tip S1-S	UPS1S	PWR Air Perio Subgingival Tips (40/box)	APPERIOTIPS
PWR Piezo K10 Tip Wrench	UPK10	PWR Air K9 Wrench Tip	APK9
PWR Piezo Bottle Connection O-Ring Set	OPORING	PWR Air Unit Water Supply Hose	APH20HOSE
PWR Plastic Bottle, 500 ml	UPBTL	PWR Air Unit Supply Hose	APAIRHOSE
PWR Bottle Cap, Grey	UPBTLCAPGY	PWR Air O-Ring Set	APORINGS
PWR Piezo Peristaltic Water Pump	UPPUMP	PWR Air Irrigation Circuit Disinfection	APICDKIT
PWR Air Powder Chamber Prophy	WAPPCPROPHY	PWR Air Unit Water Filter	APFILTER
PWR Air Powder Chamber Perio	APPCPERIO	PWR Air Cleaning Needle Kit	APNDLKIT

PWR Pair • Piezoelectric Scaling

Put the PWR in power scaling with the PWR Pair. Our Piezo technology allows for a comfortable experience for patient and clinician. The advanced feedback system gives the clinician perfect control over performance to improve results and increase comfort.

Slim LED Handpiece

The PWR Pair sports a lightweight, slim piezo handpiece with circular LED light, which prevents the “fog” effect due to the cavitation of the irrigation.

500ml Bottle

The PWR Pair 500ml Bottle illuminates to easily identify the level of liquid remaining. The 500ml Bottle also allows you to use different medicaments such as chlorhexidine for perio or sodium hypochlorite for endo.

Three Different Power Zones

The PWR Pair includes Perio, Endo, and Restorative treatment functions. These three power zones provide the ability to use the PWR Pair on more patients.

Pulse Mode

The PWR Pair Pulse Mode (Level 5 of Restorative/ Prosthetics function) transforms the traditional sine wave formed piezoelectric oscillation (standard function) to a new profile, characterized by specific power peaks, which results in optimized performance for prosthetics.

Irrigation System

The PWR Pair includes two irrigation lines, which allow for the flexible use of water (via tap or 500 ml bottle) or medicinal solutions to improve ease of use and reduce clinician downtime.

Soft Mode

Avoid excessive piezoelectric oscillation using Soft Mode, which allows gentle and efficient piezo tip movement. This results in **nearly painless treatment** for patients and optimal comfort for clinicians.



Constant Performance Monitoring

The PWR Pair Feedback System provides constant and optimal tuning of PWR Piezo Tip movement and automatically detects and adjusts the necessary power according to the clinical need. Proper power usage results in the quickest and most effective treatment while reducing treatment time.



“

The PWR Pair device has transformed the way I deliver care. Its precise settings for air polishing make biofilm disruption effortless, improving both patient comfort and clinical outcomes. The versatility between supra- and subgingival treatments streamlines my workflow, saving time without compromising quality. This is a game-changer for any practice focused on elevating patient care and operational efficiency!

AMBER AUGER, MPH, RDH


”




Soft and Pulse Modes

Soft Mode is available when your priority is a gentle touch, or use Pulse Mode when dealing with stubborn calculus.

**Perio
Soft**

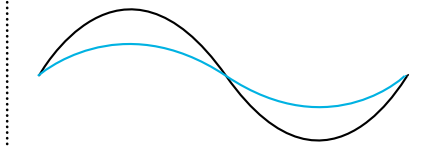

**Restorative
Pulse**

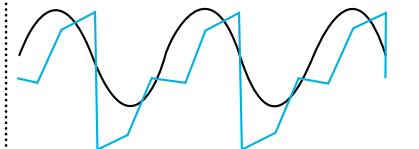

PWR Pair • Piezo Tips & Accessories

Your volume of patients varies every day as do their needs. Not all patients are the same. Put the **PWR** in power scaling with our comprehensive variety of **PWR Piezo Tips***, compatible with our **PWR Piezo** and **PWR Pair** Devices.

SCALING

Our selection of PWR Scaling Piezo Tips includes standard piezo tips for calculus removal and biofilm elimination. With six designs to choose from, there is a PWR Piezo Tip no matter the patient.



S1 PWR Piezo Tip **UPS1**



S2 PWR Piezo Tip **UPS2**



S4 PWR Piezo Tip **UPS4**



S1-S PWR Piezo Tip **UPS1S**



S3 PWR Piezo Tip **UPS3**



S5 PWR Piezo Tip **UPS5**

PERIO

Our selection of PWR Perio Piezo Tips includes 10 universal and anatomic piezo tips for calculus and biofilm removal. Specializing in both supra- and subgingival treatments, there is a PWR Perio Piezo Tip perfect for all of your patients' needs.



P2 PWR Piezo Tip **UPP2**



P10 PWR Piezo Tip **UPP10**



P12 PWR Piezo Tip **UPP12**



P3 PWR Piezo Tip **UPP3**



P11 PWR Piezo Tip **UPP11**



P13 PWR Piezo Tip **UPP13**



P14 PWR Piezo Tip **UPP14**



P16L PWR Piezo Tip **UPP16L**



P16R PWR Piezo Tip **UPP16R**



P15 PWR Piezo Tip **UPP15**

Implant Cleaning



PWR piezo Implant
Cleaning Tip (QTY:5)
UPIC1



PWR Piezo Implant
Cleaning Tip Holder
UPICTH

ENDODONTIC

Discover the PWR of minimally invasive retrograde root canal preparation. PWR Endodontic Retro Piezo Tips assure a Best in Practice solution to the root canal access issue. Their thin and 90° angled structure allows for a conservative cut through a minimum width bone window. These advantages facilitate wound healing and increase patient comfort.

From exposure and removal of root pins to fast and effective removal of calcification in the coronal third of the root and removal of fractured root inserts—the PWR Endodontic Revision Piezo Tips cover the complete spectrum of endodontic revision treatment.



ER1 PWR Piezo Tip **UPER1**



ER2 PWR Piezo Tip **UPER2**



ER3 PWR Piezo Tip **UPER3**



ER4 PWR Piezo Tip **UPER4**



ER5 PWR Piezo Tip **UPER5**



RD3 PWR Piezo Tip,
Diamond Coated **UPRD3**



RD4 PWR Piezo Tip,
Diamond Coated **UPRD4**



RD5 PWR Piezo Tip,
Diamond Coated **UPRD5**

TIP SETS*

PWR Piezo Tip Set Scaling	UPSETS
S1 PWR Piezo Tip	UPS1
S2 PWR Piezo Tip	UPS2
S3 PWR Piezo Tip	UPS3
S4 PWR Piezo Tip	UPS4
S5 PWR Piezo Tip	UPS5
K10 Tip Wrench (QTY: 5)	UPK10

PWR Piezo Tip Set Maintenance	UPSETM
S1 PWR Piezo Tip	UPS1
S2 PWR Piezo Tip	UPS2
S5 PWR Piezo Tip	UPS5
P2 PWR Piezo Tip	UPP2
P3 PWR Piezo Tip	UPP3
K10 Tip Wrench (QTY: 5)	UPK10

PWR Piezo Tip Set Perio Maintenance	UPSETPM
P10 PWR Piezo Tip	UPP10
P11 PWR Piezo Tip	UPP11
P12 PWR Piezo Tip	UPP12
P13 PWR Piezo Tip	UPP13
P14 PWR Piezo Tip	UPP14
K10 Tip Wrench (QTY: 5)	UPK10

PWR Piezo Tip Set Endo Revision	UPSETER
ER1 PWR Piezo Tip	UPER1
ER2 PWR Piezo Tip	UPER2
ER3 PWR Piezo Tip	UPER3
ER4 PWR Piezo Tip	UPER4
ER5 PWR Piezo Tip	UPER5
K10 Tip Wrench (QTY: 5)	UPK10

PWR Piezo Tip Set Perio Treatment	UPSETPT
P2 PWR Piezo Tip	UPP2
P3 PWR Piezo Tip	UPP3
P15 PWR Piezo Tip	UPP15
P16L PWR Piezo Tip	UPP16L
P16R PWR Piezo Tip	UPP16R
K10 Tip Wrench (QTY: 5)	UPK10

PWR Piezo Tip Set Implant Cleaning	UPSETIC
PWR Piezo Tip— Implant Cleaning (QTY:5)	UPIC1
PWR Piezo Implant Cleaning Tip Holder	UPICTH
K10 Tip Wrench	UPK10

PWR Pair • Air Polishing

The PWR Pair utilizes PWR Prophylaxis Powder combined with air and water to supra- and subgingivally remove biofilm, and the tough stains caused by everyday foods and beverages. Air polishing also helps remove interproximal biofilm and food debris.





Air Polishing is more accessible for biofilm removal in areas of pocketing and around implants. By using the PWR Pair Piezo Scaler and Air Polisher, we are able to provide thorough hygiene treatment in areas that are difficult to reach with conventional instrumentation.

DR. EDWARD P. ALLEN, DDS



Slim Stainless Steel Air-Polishing Handpieces

The PWR Pair includes three different slim air polishing handpieces: 90°, 120°, and Perio. Being made from steel and resilient to cracking, you can practice with the power peace of mind: fewer interruptions and less time spent replacing handpieces.

Perio Function

The one-touch Perio Function of the PWR Pair reduces the air pressure for a safe and comfortable subgingival treatment.

Supra- and Subgingival

The PWR Pair allows you to switch easily during treatment between supra- and subgingival biofilm, stain and calculus removal with the simple push of a button, Prophy (Supragingival) or Perio (Subgingival) functions, respectively.

Refill Function

The PWR Pair Refill Function improves ease of use and decreases time to complete maintenance with easy powder chamber removal without switching the device off.

Continuous Flow of Air

PWR Pair Device provides a continuous flow of air which helps to prevent clogging. Practice with the power peace of mind: less interruptions and less time spent trying to unclog lines.



Subgingival Perio Tip

Flexible and soft, the disposable Subgingival Perio Tip (included with the PWR Pair) comes sterile and easily attaches to the PWR Perio Handpiece.

PWR Pair • Air Polishing Handpieces & Powders

The PWR of Three

Achieve Best in Practice access with the PWR Pair 90° and 120° Angled Handpieces for supra- and subgingival use in periodontal pockets up to 5 mm deep. The PWR Pair Perio Handpiece equipped with the disposable sterile PWR Subgingival Perio Tip (included with the PWR Pair) provides optical access in pockets deeper than 5 mm. PWR Air Polishing Handpieces easily attach with just a click and are safely secured, thanks to a dedicated security connection system.

Did You Know?

All PWR Air Polishing Handpieces are resilient to cracking, being made from steel. Practice with the power peace of mind: fewer interruptions and less time spent replacing handpieces.



PWR Air Perio Subgingival Tips
APPERIOTIPS
QTY: 40



PWR Air 90° Handpiece
APHP90



PWR Air 120° Handpiece
APHP120

Included
with the
PWR Pair



PWR Air Perio Handpiece
APPPERIO

Included
with the
PWR Pair

PWR Prophylaxis Powders

PWR Prophylaxis Powders remove supra- and subgingival biofilm, stains and pigmentations along with treating peri-implantitis and cleaning of prosthesis and restorative materials. Choose from PWR Perio Powder (Glycine) and PWR Classic Powder (Sodium Bicarbonate), both coming in mint and neutral flavors. Our specially developed PWR Prophylaxis Powders also come with hydrophobic effect, which helps decrease clogging in the PWR Pair.



Powder	Flavor	Type	SKU
Perio	Mint	Glycine	AP-PM
Perio	Neutral	Glycine	AP-PN
Classic	Mint	Sodium Bicarbonate	AP-CM
Classic	Neutral	Sodium Bicarbonate	AP-CN

Tips for PWRful Procedures

PWR Prophylaxis Powders have you covered no matter the situation.

	Classic	Perio
Supragingival Biofilm Removal	●	●
Removal of Stains & Pigmentations on Enamel	●	
Subgingival Biofilm Removal		●
Periodontal Pocket Detoxification		●
Implant Maintenance		●
Peri-Implantitis Treatment		●
Fixed Orthodontic Appliances	●	●
Cleaning of Prosthesis & Restoration Materials		●
Restorations Maintenance		●

EverCare

SERVICE AND REPAIRS

.....
Let Us Take Care of the Maintenance



Ensure your **PWR devices** are maintained for optimal performance with our **Best in Practice EverCare Service & Repairs PWR Program**.



Annual
maintenance
of your PWR
device



Real-time
online status
tracking



Discounts
on repairs &
rental devices



Two-week
turnaround
time



Complimentary
shipping



Exclusive offers
for members who
register their
device



For More Membership Information, Visit
EverCareService.com

HuFriedyGroup.com

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.

©2025 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-914 /2505

HuFriedyGroup
The Best In Practice