

# PWR Piezo

## Piezoelectric Scaler



Power  
Equipment



 **HuFriedyGroup**  
The Best In Practice

# Introducing **The PWR Piezo**

Each patient provides unique challenges and only you have the power to choose how to apply your wealth of knowledge and skills. Power ON with the **PWR Piezo** Piezoelectric Scaler from HuFriedyGroup.

Choose Best in Practice clinical outcomes with the PWR Piezo:  
Built for your patient, designed for your power.

Learn more about PWR Piezo Scaling & Air Polishing at  
[HuFriedyGroup.com/PWR-ON](http://HuFriedyGroup.com/PWR-ON)





# The **PWR** of Versatility




## **PWR PIEZO SCALING**

Use the PWR Piezo to effectively and efficiently remove your patients' stains and calculus.



## **PWR PIEZO TIPS**

Power scale the way you choose with the 24 PWR Piezo Tips, including scaling, perio and endodontic PWR Piezo Tips for each indication, each treatment, and each patient.

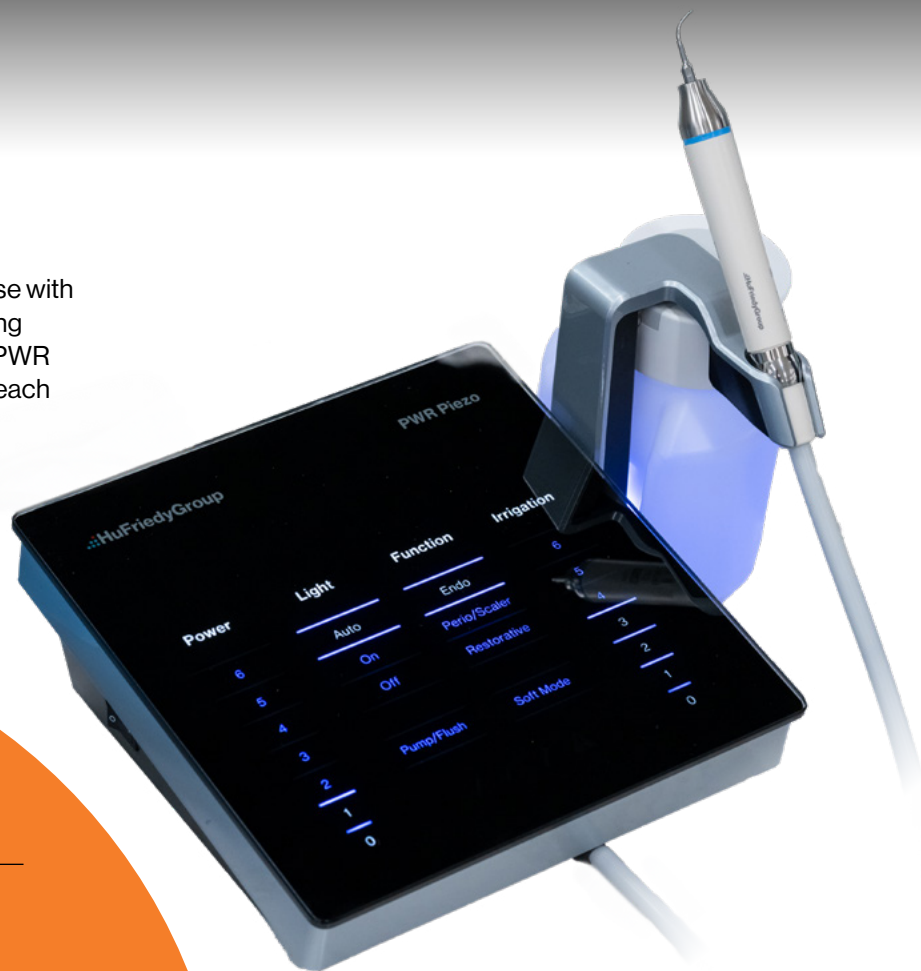

## **What's in the Box**

### **PWR Piezo Piezoelectric Scaler**

PWR Piezo LED handpiece  
PWR Foot Pedal  
PWR Piezo Tip P3  
PWR Piezo Tip S1-S  
PWR Piezo K10 Tip Wrench  
PWR Piezo Bottle Connection O-Ring Set  
PWR Plastic Bottle, 500 ml  
PWR Bottle Cap, Grey  
PWR Piezo Peristaltic Water Pump

### **UPPIEZO**

**UPHPLED**  
**UPFTPEDAL**  
**UPP3**  
**UPS1S**  
**UPK10**  
**UPORING**  
**UPBTL**  
**UPBTLCAPGY**  
**UPPUMP**

# PWR Piezo • Piezoelectric Scaling

Put the PWR in power scaling with the **PWR Piezo**. Our Piezo technology allows for a comfortable experience for patient and clinician. The advanced feedback system gives the clinician perfect control over performance to improve results and increase patient comfort.

## CONSTANT PERFORMANCE MONITORING


The PWR Piezo Feedback System provides constant and optimal tuning of PWR Piezo Tip movement and automatically detects and adjusts the necessary power according to the clinical need. Proper power usage results in the quickest and most effective treatment while reducing treatment time.

### Pulse Mode

The PWR Piezo Pulse Mode (Level 5 of Restorative/Prosthetics function) transforms the traditional sine wave formed piezoelectric oscillation (standard function) to a new profile, characterized by specific power peaks, which results in optimized performance for prosthetics.

### Soft Mode

Avoid excessive piezoelectric oscillation using Soft Mode, which allows gentle and efficient piezo tip movement. This results in **nearly painless treatment** for patients and optimal comfort for clinicians.

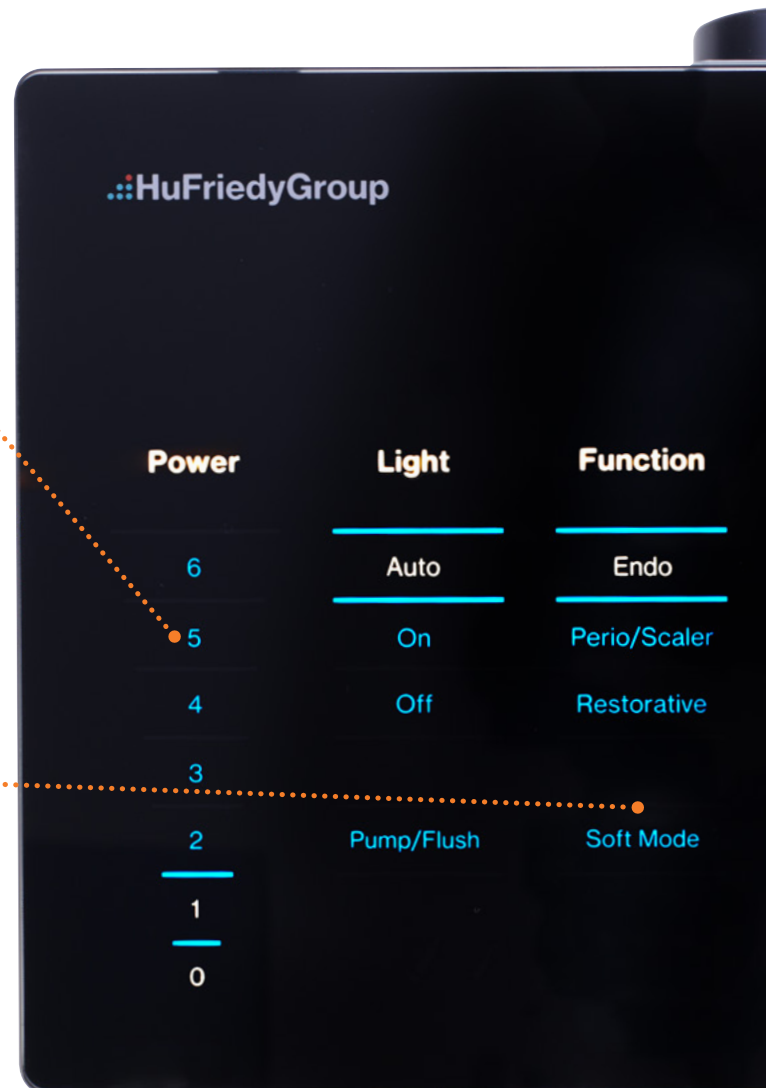

PWR PIEZO UPPIEZO



“

Since integrating PWR Piezo Scaling into our practice, patient comfort has noticeably increased, thanks to the innovative Soft Mode feature of the PWR Piezo device. This mode offers gentle, precise scaling that prioritizes the patient's experience without compromising clinical efficacy. Even patients with high sensitivity report a more comfortable visit, which allows us to focus on thorough treatment without the typical stress on their end.

**AMBER AUGER, MPH, RDH**

”



### Wide Range of Treatments

Power scale the way you choose with the 24 PWR Piezo Tips, including scaling, perio and endodontic PWR Piezo Tips for each indication, each treatment, and each patient.

### Slim LED Handpiece

The PWR Piezo sports a lightweight, slim piezo handpiece with circular LED light, which prevents the “fog” effect due to the cavitation of the irrigation.

### 500ml Bottle & Irrigation System

The PWR Piezo 500ml Bottle illuminates to easily identify the level of liquid remaining. This 500ml Bottle allows you to use water or medicaments such as chlorhexidine for perio or sodium hypochlorite for endo.

# PWR Piezo • Piezo Tips & Accessories

Your volume of patients varies every day as do their needs. Not all patients are the same. Put the **PWR** in power scaling with our comprehensive variety of **PWR Piezo Tips\***, compatible with our **PWR Piezo** and **PWR Pair** Devices.

## SCALING

Our selection of PWR Scaling Piezo Tips includes standard piezo tips for calculus removal and biofilm elimination. With six designs to choose from, there is a PWR Piezo Tip no matter the patient.



S1 PWR Piezo Tip **UPS1**



S2 PWR Piezo Tip **UPS2**



S4 PWR Piezo Tip **UPS4**



S1-S PWR Piezo Tip **UPS1S**



S3 PWR Piezo Tip **UPS3**



S5 PWR Piezo Tip **UPS5**

## PERIO

Our selection of PWR Perio Piezo Tips includes 10 universal and anatomic piezo tips for calculus and biofilm removal. Specializing in both supra- and subgingival treatments, there is a PWR Perio Piezo Tip perfect for all of your patients' needs.



P2 PWR Piezo Tip **UPP2**



P10 PWR Piezo Tip **UPP10**



P12 PWR Piezo Tip **UPP12**



P3 PWR Piezo Tip **UPP3**



P11 PWR Piezo Tip **UPP11**



P13 PWR Piezo Tip **UPP13**



P14 PWR Piezo Tip **UPP14**



P16L PWR Piezo Tip **UPP16L**



P16R PWR Piezo Tip **UPP16R**



P15 PWR Piezo Tip **UPP15**

## Implant Cleaning



PWR piezo Implant  
Cleaning Tip (QTY:5)  
**UPIC1**



PWR Piezo Implant  
Cleaning Tip Holder  
**UPICTH**

## ENDODONTIC

Discover the PWR of minimally invasive retrograde root canal preparation. PWR Endodontic Retro Piezo Tips assure a Best in Practice solution to the root canal access issue. Their thin and 90° angled structure allows for a conservative cut through a minimum width bone window. These advantages facilitate wound healing and increase patient comfort.

From exposure and removal of root pins to fast and effective removal of calcification in the coronal third of the root and removal of fractured root inserts—the PWR Endodontic Revision Piezo Tips cover the complete spectrum of endodontic revision treatment.



ER1 PWR Piezo Tip **UPER1**



ER2 PWR Piezo Tip **UPER2**



ER3 PWR Piezo Tip **UPER3**



ER4 PWR Piezo Tip **UPER4**



ER5 PWR Piezo Tip **UPER5**



RD3 PWR Piezo Tip,  
Diamond Coated **UPRD3**



RD4 PWR Piezo Tip,  
Diamond Coated **UPRD4**



RD5 PWR Piezo Tip,  
Diamond Coated **UPRD5**

## TIP SETS\*

PWR Piezo Tip Set Scaling	UPSETS
S1 PWR Piezo Tip	<b>UPS1</b>
S2 PWR Piezo Tip	<b>UPS2</b>
S3 PWR Piezo Tip	<b>UPS3</b>
S4 PWR Piezo Tip	<b>UPS4</b>
S5 PWR Piezo Tip	<b>UPS5</b>
K10 Tip Wrench (QTY: 5)	<b>UPK10</b>

PWR Piezo Tip Set Maintenance	UPSETM
S1 PWR Piezo Tip	<b>UPS1</b>
S2 PWR Piezo Tip	<b>UPS2</b>
S5 PWR Piezo Tip	<b>UPS5</b>
P2 PWR Piezo Tip	<b>UPP2</b>
P3 PWR Piezo Tip	<b>UPP3</b>
K10 Tip Wrench (QTY: 5)	<b>UPK10</b>

PWR Piezo Tip Set Perio Maintenance	UPSETPM
P10 PWR Piezo Tip	<b>UPP10</b>
P11 PWR Piezo Tip	<b>UPP11</b>
P12 PWR Piezo Tip	<b>UPP12</b>
P13 PWR Piezo Tip	<b>UPP13</b>
P14 PWR Piezo Tip	<b>UPP14</b>
K10 Tip Wrench (QTY: 5)	<b>UPK10</b>

PWR Piezo Tip Set Endo Revision	UPSETER
ER1 PWR Piezo Tip	<b>UPER1</b>
ER2 PWR Piezo Tip	<b>UPER2</b>
ER3 PWR Piezo Tip	<b>UPER3</b>
ER4 PWR Piezo Tip	<b>UPER4</b>
ER5 PWR Piezo Tip	<b>UPER5</b>
K10 Tip Wrench (QTY: 5)	<b>UPK10</b>

PWR Piezo Tip Set Perio Treatment	UPSETPT
P2 PWR Piezo Tip	<b>UPP2</b>
P3 PWR Piezo Tip	<b>UPP3</b>
P15 PWR Piezo Tip	<b>UPP15</b>
P16L PWR Piezo Tip	<b>UPP16L</b>
P16R PWR Piezo Tip	<b>UPP16R</b>
K10 Tip Wrench (QTY: 5)	<b>UPK10</b>

PWR Piezo Tip Set Implant Cleaning	UPSETIC
PWR Piezo Tip— Implant Cleaning (QTY:5)	<b>UPIC1</b>
PWR Piezo Implant Cleaning Tip Holder	<b>UPICTH</b>
K10 Tip Wrench	<b>UPK10</b>

# EverCare

## SERVICE AND REPAIRS

.....

Let Us Take Care of the Maintenance



Ensure your **PWR devices** are maintained for optimal performance with our **Best in Practice EverCare Service & Repairs PWR Program**.



Annual  
maintenance  
of your PWR  
device



Real-time  
online status  
tracking



Discounts  
on repairs &  
rental devices



Two-week  
turnaround  
time



Complimentary  
shipping



Exclusive offers  
for members who  
register their  
device



For More Membership Information, Visit  
**EverCareService.com**

\*PWR Piezo Device only comes with UPS1-S and UPS3 PWR Piezo Tips. All additional PWR Piezo Tips and PWR Piezo Tip Sets sold separately.

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.

©2025 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-915 /2502

**HuFriedyGroup**  
The Best In Practice