

# NEVI™

# Sickle Scaler Collection



Instrumentation and  
Procedural Products



Innovative Designs  
For Effective Deposit  
Removal

 **HuFriedyGroup**  
The Best In Practice

# Nevi™ Sickle Scalers

Superior Interproximal  
access featuring  
**EverEdge™ 2.0 Technology**

The Nevi™ Series of instruments were invented by legendary HuFriedyGroup artisan Neville Hammond with significant input from clinicians. Each of these anterior and posterior instruments were designed for exceptional efficacy, efficiency and comfort. The optimized blade angles provide superior interproximal access and reach to help make hand scaling the best part of the appointment.

Each Nevi™ Instrument comes with our **EverEdge™ 2.0 Blade Technology** for longer-lasting edges and a more efficient, comfortable scaling experience.

Resin Eight™  
Handle

EverEdge™  
2.0  
Stays  
Sharper  
Longer

Harmony™  
Ergonomic  
Handle

Popular #9  
Handle

#G9 Steel Grip  
Handle



# Anterior Sickle Scalers

## NEVI™ 1



## NEVI™ 1/H5




### Two distinct working ends provide solutions for anterior scaling needs

- **Sickle End:** A moderate length, thin scaler with two cutting edges and a blade that tapers to a point. Ideal for removing anterior interproximal calculus and stain
- **Disc End:** A moderate sized, sharp rounded disc. Excellent for removal of anterior lingual stain

The Nevi 1™ Scaler is also available with a longer sickle working end (Nevi1/H5).



*Mesial Left Central Incisor*

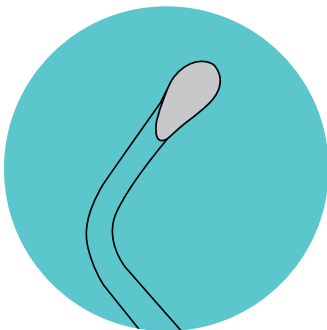

*Nevi™ 1 Disc End*



*Lingual Left Lateral Incisor*

# Posterior Sickie Scalers

## NEVI™ 2



- Perfect for **light to moderate interproximal** calculus removal
- The thinnest dual tapered blades in the Nevi™ Instrument collection remove calculus with minimal tissue trauma
- Deep contra-angles are paired with extra blade length for easy access under contacts



*Mesial Right First Molar*



*Distal Left Second Molar*

## NEVI™ 3



- Used for light to moderate calculus removal in **tight, narrow posterior interproximal areas**
- Ultra slender, dual cutting blades for effective placement
- Contra-angle bends in shank offer superior reach and adaptation to posterior interproximal areas



*Mesial Left First Molar*



*Distal Left Second Premolar*



# Posterior Sickle Scalers

## NEVI™ 4



- Ideal for scaling **heavy posterior interproximal** deposits
- Increased shank size and rigidity generates effective removal of tenacious, heavy and burnished deposits
- Working ends with dual tapered cutting edges allow for precise interproximal access
- Largest blade in the Nevi™ Instrument collection



*Distal Right Premolar*

## NEVI™ 4 Lite




- 10% thinner than the standard Nevi™ 4 blade for enhanced access to tight interproximal spaces



*Mesial Left First Molar*

# Nevi™ Sickle Scalers Featuring EverEdge™ 2.0 Technology

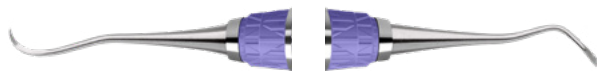


## Nevi™ 1 Anterior Sickle Scaler

|                 |                    |
|-----------------|--------------------|
| Harmony         | <b>SCNEVI1XE2</b>  |
| #9              | <b>SCNEVI19E2</b>  |
| #G9, Steel Grip | <b>SCNEVI1G9E2</b> |
| #C8, Lavender   | <b>SCNEVI1C8E2</b> |




## Nevi™ 1/H5 Anterior Sickle/Discoïd Scaler

|                 |                     |
|-----------------|---------------------|
| Harmony         | <b>SNEV1/H5XE2</b>  |
| #9              | <b>SNEV1/H59E2</b>  |
| #G9, Steel Grip | <b>SNEV1/H5G9E2</b> |
| #C8, Lavender   | <b>SNEV1/H5C8E2</b> |

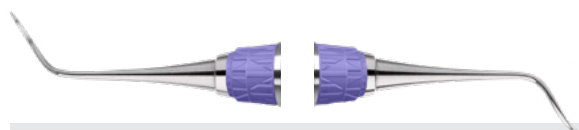


## Nevi™ 2 Posterior Sickle Scaler

|                 |                    |
|-----------------|--------------------|
| Harmony         | <b>SCNEVI2XE2</b>  |
| #9              | <b>SCNEVI29E2</b>  |
| #G9, Steel Grip | <b>SCNEVI2G9E2</b> |
| #C8, Red        | <b>SCNEVI2C8E2</b> |

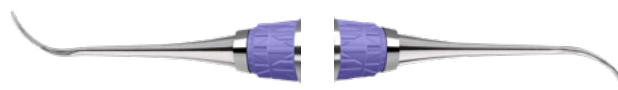

## Nevi™ 3 Posterior Sickle Scaler

|                 |                    |
|-----------------|--------------------|
| Harmony         | <b>SCNEVI3XE2</b>  |
| #9              | <b>SCNEVI39E2</b>  |
| #G9, Steel Grip | <b>SCNEVI3G9E2</b> |
| #C8, Red        | <b>SCNEVI3C8E2</b> |

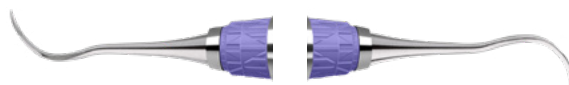

## Nevi™ 4 Posterior Sickle Scaler

|                 |                    |
|-----------------|--------------------|
| Harmony         | <b>SCNEVI4XE2</b>  |
| #9              | <b>SCNEVI49E2</b>  |
| #G9, Steel Grip | <b>SCNEVI4G9E2</b> |
| #C8, Red        | <b>SCNEVI4C8E2</b> |



## Nevi™ 4 Lite Posterior Sickle Scaler

|                 |                     |
|-----------------|---------------------|
| Harmony         | <b>SCNEVI4LXE2</b>  |
| #9              | <b>SCNEVI4L9E2</b>  |
| #G9, Steel Grip | <b>SCNEVI4LG9E2</b> |
| #C8, Red        | <b>SCNEVI4LC8E2</b> |

Several Nevi™ Scaler designs are also available with Immunity Steel™ Handles. Please see the full HuFriedyGroup Product Catalog for more information.

To learn more about our Best in Practice scalers and curettes, visit  
[HuFriedyGroup.com/Scalers](https://www.hufriedygroup.com/scalers)