


COMPLIANT TRANSPORTATION OF REUSABLE CONTAMINATED SHARPS

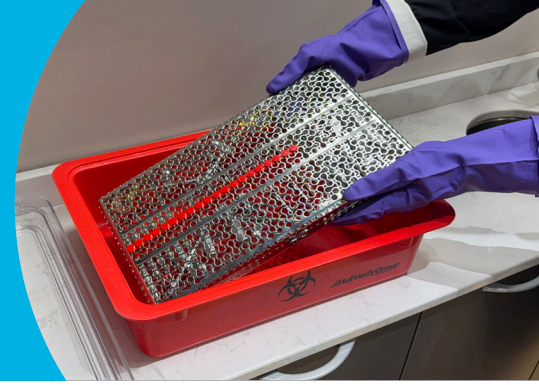

Understand how to comply with the OSHA bloodborne pathogen standard 1910:1030*

To prevent occupational exposure, proper Personal Protective Equipment (PPE) must be worn, which includes utility gloves, a gown or lab coat, a face shield or mask, and eye protection.



Disposable contaminated sharps should not be transported and should be immediately discarded in sharps containers located as close as is feasible to the area of use.

After use, contaminated reusable sharps should be placed in appropriate container as soon as possible until they are reprocessed.

- Containers shall be closed immediately prior to removal or replacement to prevent spillage or protrusion of contents.

Containers used for transportation should be:

- Closeable
- Constructed to contain all contents without exposure
- Leakproof on sides and bottom
- Puncture resistant
- Properly labeled or color-coded

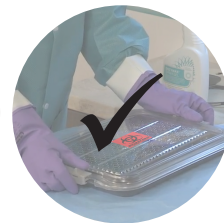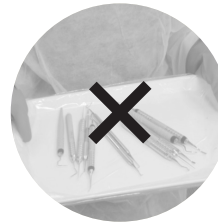

If leakage is possible, a secondary container should be used

Warning labels should be:

- Fluorescent orange or orange-red, or a contrasting color from the container not able to fade, erase, or fall off of the container.
- Excluded if the container or bag is red or if they are uniquely colored in the office and are known to contain contaminated materials.
Note: if contaminated items will leave the facility, proper labeling is required.



Ordering information:
Ref#: 1960
Biohazard Warning Labels
Size 3" x 3"
Quantity 100/roll

HuFriedyGroup has an assortment of products to help you follow the standards established by OSHA, including:

- [IMS Cassettes](#)
- [Lilac Utility Gloves](#)
- [Face Shields](#)
- [Organizational Tubs](#)
- [Face Masks](#)
- [Safety Eyewear](#)
- [Warning Labels](#)

Please note this document is intended as a resource to help with proper transportation of reusable contaminated sharps but it does not cover the entire standard. For full details, you can [review the standard here](#).

*Occupational Safety and Health Standards 1910.1030: bloodborne pathogens.
<https://www.osha.gov/laws-regs/regulations/standardnumber/1910/1910.1030>

Solutions for Transporting Cassettes and Instruments



Large Transport Tub and Cover

- SKU: **IMTRANSPORT**
- Dimensions: 12.6" x 17.4" x 5.5"
- Compatible cassette sizes:
Signa Stat, Large 8x11, Large Thin Ortho, Container Cassette, Full Din, Oral Surgery



Signature Series Tub Only – Red

- SKU: **IMS-1410**
- Dimensions (tub and cover together): 13.75" x 11.1875" x 4"
- Compatible cassette sizes:
Exam, Small 8 or 10, Small Thin Ortho, Double-Deckers, ¼ DIN, ½ DIN



Locking Tub Cover

- SKU: **IMS-1420**
- Compatible with Signature Series Tub (IMS-1410)

For more information, visit [HuFriedyGroup.com](https://www.HuFriedyGroup.com)