

Introducing AIR-FLOW THERAPY

AIR-FLOW Therapy is completely different than other air polishing systems. This advanced technology offers benefits to both the clinician and patient that are difficult to replicate with other procedures. AIR-FLOW Therapy, along with glycine or erythritol powder, is safe and effective for tissue, enamel, root surfaces, porcelain composite restorations, orthodontic appliances, implants and implant abutments.



WHAT IT DOES:

- Polishes teeth
- Removes stain
- Disrupts subgingival biofilm

PATIENT BENEFITS:

- Increased comfort
- Healthy gingiva
- Grit-free polishing
- Faster procedure

CLINICIAN BENEFITS:

- Effective biofilm disruption
- Efficient treatments
- Happier patients

BRING AIR-FLOW TO YOUR PRACTICE WITH A HANDY DEVICE! A great introduction to AIR-FLOW Therapy, the handy 3.0 is a lightweight handheld device that connects to your dental unit. Designed with a team of ergonomics experts, the slim handpiece is easy for the clinician to grip and rotate, while the position of the powder chamber ensures an unrestricted view of the treatment area.

ORDER YOUR HANDY DEVICE!


PRODUCT	PART CODES*
AIR-FLOW handy 3.0 Midwest*	FT-220#HF/A/001
AIR-FLOW handy 3.0 Kavo	FT-220#HF/A/002
AIR-FLOW handy 3.0 PERIO Midwest*	FT-221#HF/A/001
AIR-FLOW handy 3.0 PERIO Kavo	FT-221#HF/A/002
AIR-FLOW handy 3.0 Premium Midwest*	FT-225#HF/A/001
AIR-FLOW handy 3.0 Premium Kavo	FT-225#HF/A/002
handy 3.0 Handpiece	EL-540#HF/A
handy 3.0 PERIO Handpiece	EL-542#HF/A
handy 3.0 PLUS Handpiece	EL-545#HF/A

NOZZLE REFILLS	PART CODES
PERIO-FLOW nozzle refill box (40ct)	DT-113HF
PERIO-FLOW nozzle refill pack (4 boxes of 40)	DT-113HF/HF

POWDERS IN PACKS OF 4	PART CODES
AIR-FLOW Classic Powder, Lemon	DV-1001
AIR-FLOW Classic Powder, Cherry	DV-1006
AIR-FLOW Classic Powder, Mint	DV-1010
AIR-FLOW Classic Powder, Neutral	DV-1009
AIR-FLOW PERIO Powder	DV-110/A/HS
AIR-FLOW PLUS Powder	DV-117/A

*Fits standard, 4-hole connection
*W&H, NSK, Sirona and Bien Air Connections also available

For more information, visit
Hu-Friedy.com/BIOFILM

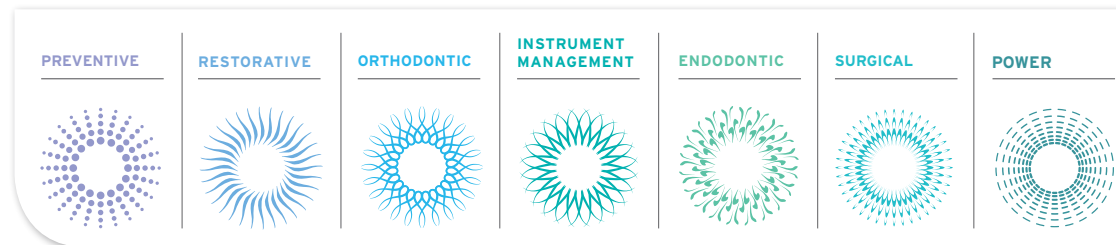

Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA
Telephone: 1-800-Hu-Friedy or 1-773-975-6100 | Website: Hu-Friedy.com

©2018 Hu-Friedy Mfg. Co., LLC. All rights reserved. [557_2] HF-3030/0518
NOTE: Kavo, Midwest, W&H, NSK, Bien Air and Sirona are trademarks or registered trademarks of their respective owners.

Hu-Friedy



AIR-FLOW® HANDY 3.0 DEVICES

THE BETTER WAY FOR
BIOFILM MANAGEMENT

Hu-Friedy  **EMS**
ELECTRO MEDICAL SYSTEMS

POWER



CLINICAL CAPABILITIES

STANDARD HANDPIECE

- Designed for Classic Powder
- Ideal for supragingival stain removal applications
- Included with standard AIR-FLOW handy 3.0 device

PLUS HANDPIECE

- Designed for fine-grain PERIO or PLUS Powder
- Ideal for supragingival biofilm removal, including on ortho and restorative materials, as well as subgingival biofilm removal up to 4mm
- Safe for perio and implant maintenance patients
- Included with AIR-FLOW handy 3.0 PREMIUM device and compatible with AIR-FLOW handy 3.0 PERIO device (sold separately)

PERIO HANDPIECE

- Designed for fine-grain PERIO or PLUS Powder
- Flexible plastic nozzle is inserted directly into pockets to flush out biofilm
- Safe for perio and implant maintenance patients
- Included with AIR-FLOW handy 3.0 PREMIUM and AIR-FLOW handy 3.0 PERIO devices

DEVICE OPTIONS

AIR-FLOW HANDY 3.0 PREMIUM

- For complete biofilm removal
- Includes both PLUS and PERIO handpieces for supragingival and subgingival applications
- Easily switch between handpieces during treatment
- Designed for optimal delivery of PERIO or Plus Powder

DEVICE OPTIONS

AIR-FLOW HANDY 3.0 PERIO

- For direct subgingival biofilm removal
- Includes PERIO handpiece for subgingival applications up to 5mm deep*
- Designed for PERIO or PLUS Powder

*AIR-FLOW handy 3.0 Perio is approved for use in pockets up to 10mm in Canada

AIR-FLOW HANDY 3.0

- For heavy stain removal
- Includes standard handpiece for mainly supragingival applications
- Designed for Classic Powder, but can be used with PERIO or PLUS Powder

AIR-FLOW PLUS POWDER DV-117/A

AIR-FLOW PLUS Powder is the first powder that is safe to use for supra- and subgingival biofilm removal but strong enough to efficiently remove light to moderate stain. The finest grain powder in dentistry, PLUS Powder's erythritol-based, minimally abrasive spray can be used on virtually every patient.

- Erythritol
- 14 µm
- 120g per bottle

AIR-FLOW PERIO POWDER DV-110/A/HS

AIR-FLOW PERIO Powder's fine, low-density grains safely and completely remove harmful biofilm and bacteria both supragingivally and subgingivally. Made from glycine, an amino acid, PERIO Powder is gentle on enamel and restorative surfaces, making the patient experience more comfortable.

- Glycine
- 25 µm
- 120g per bottle

AIR-FLOW CLASSIC POWDER

NEW FORMULA

The Classic Powder offers the well-known cleaning efficiency of sodium bicarbonate with increased patient comfort due to a smaller average grain size than standard sodium bicarbonate powders.

While the Classic Powder is primarily used for supragingival removal of biofilm and colorations, it can also be used for surface preparation before dental treatments.

- Sodium Bicarbonate
- 40 µm
- 300g per bottle

AVAILABLE FLAVORS:

Lemon | DV-1001 Mint | DV-1010
Cherry | DV-1006 Neutral | DV-1009

