



HU-FRIEDY SURGICAL KIT

ATRAUMATIC EXTRACTION KIT




PART CODE	DESCRIPTION
EPTSMS	Luxating Hybrid Elevator, Small Straight
EPTSMC	Luxating Hybrid, Small Curved
EPTSSP	Luxating Hybrid Elevator, Straight Spade
FAF150	Apical Forceps 150
FAF151	Apical Forceps 151
CM2/4S6	2/4 Molt DE Curette, #6 Handle
NHCW	Crile-Wood Needle Holder
S16SC	16 G-F Scissors Super Cut
10-130-05	5 Scalpel Handle, Round Straight
P96	9 Molt Periosteal Elevator
IM1002	6-Pack Signature Series® Instrument Clip
IM4161	16 Instrument Signature Series® Cassette, Gray
IMATRAEXT	Complete Atraumatic Extraction Kit

ATRAUMATIC EXTRACTION KIT

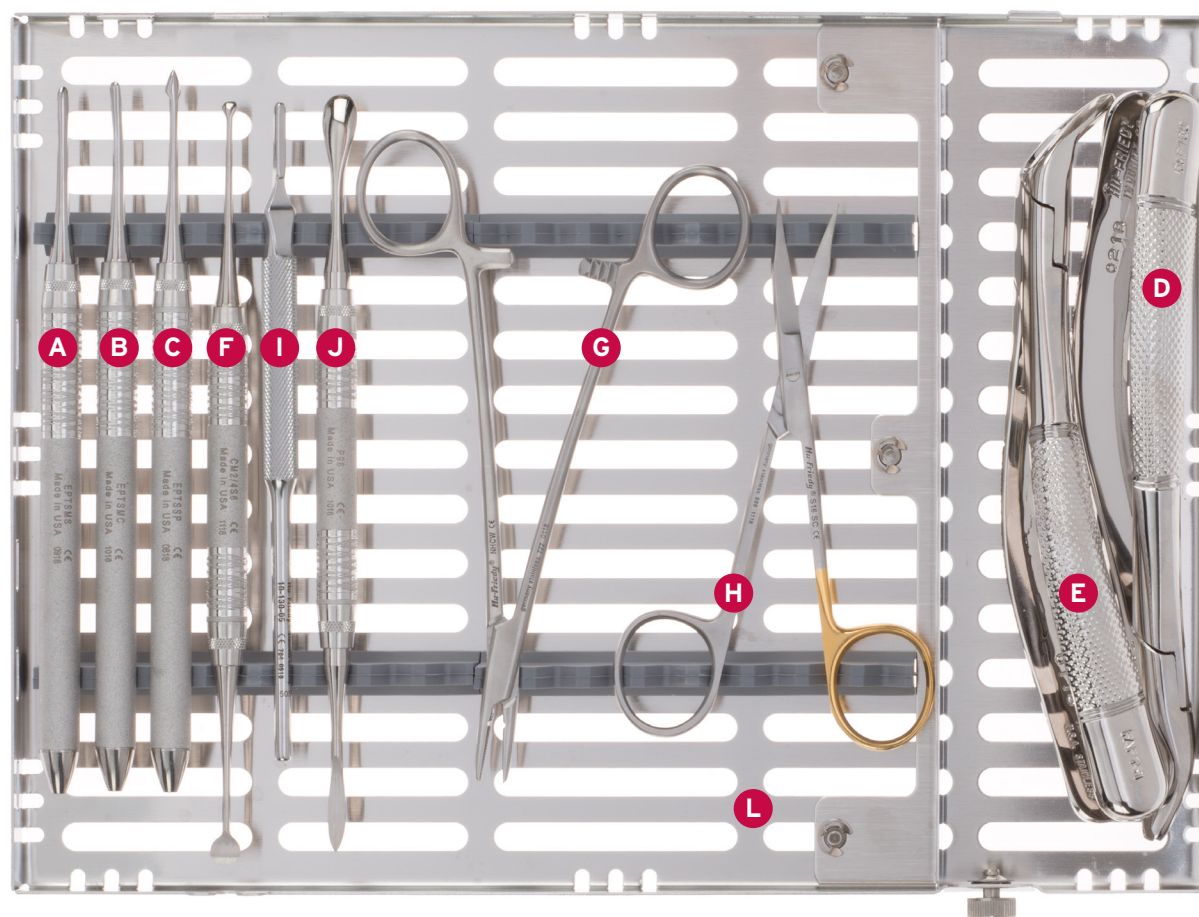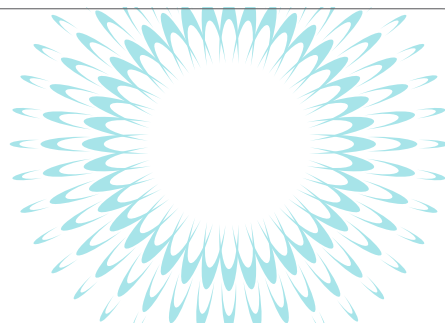

ATRAUMATIC EXTRACTION KIT

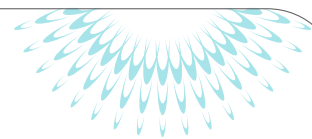

A EPTSMS
Luxating Hybrid Elevator, Small Straight
 A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



B EPTSMC
Luxating Hybrid, Small Curved
 A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



C EPTSSP
Luxating Hybrid Elevator, Straight Spade
 A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



D FAF150
Apical Forceps 150
 Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



E FAF151
Apical Forceps 151
 Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



F CM2/4S6
2/4 Molt DE Curette, #6 Handle
 Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



G NHCW
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



H S16SC
16 Goldman-Fox Scissors, Super Cut
 Used for trimming tissue or cutting sutures. Featuring tapered blade tips, tips that are not sharp on outside and are used as tonotomy scissors. 12.5cm with 1 serrated blade



I 10-130-05
5 Scalpel Handle, Round Straight
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



J P96
9 Molt Periosteal Elevator
 Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



K IM1002
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



L IM4161
16 Instrument Signature Series® Cassette, Gray
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.

