



HU-FRIEDY SURGICAL KIT

EXTRACTION KIT



PART CODE	DESCRIPTION
IM4163E2C4	16 Instrument Signature Series® Cassette, Orange
10-130-05	5 Scalpel Handle, Round Straight
P96	9 Molt Periosteal Elevator
NMH	Mayo-Hegar Needle Holder
S16	16 Goldman-Fox Scissors
CL866	86 Lucas DE Surgical Curette, #6 Handle
E301	301 Apexo Elevator
E34S	34S Seldin Straight Elevator
BF2X6	2X Miller-Colburn Bone File
EHB13/146	13/14 Heidbrink DE Root Pick
RF	Friedman Rongeurs
TP43	43 Adson-Brown Tissue Pliers
EPTSSP	Luxating Hybrid Elevator, Straight Spade
IMEXTKIT	Complete Extraction Kit

EXTRACTION KIT

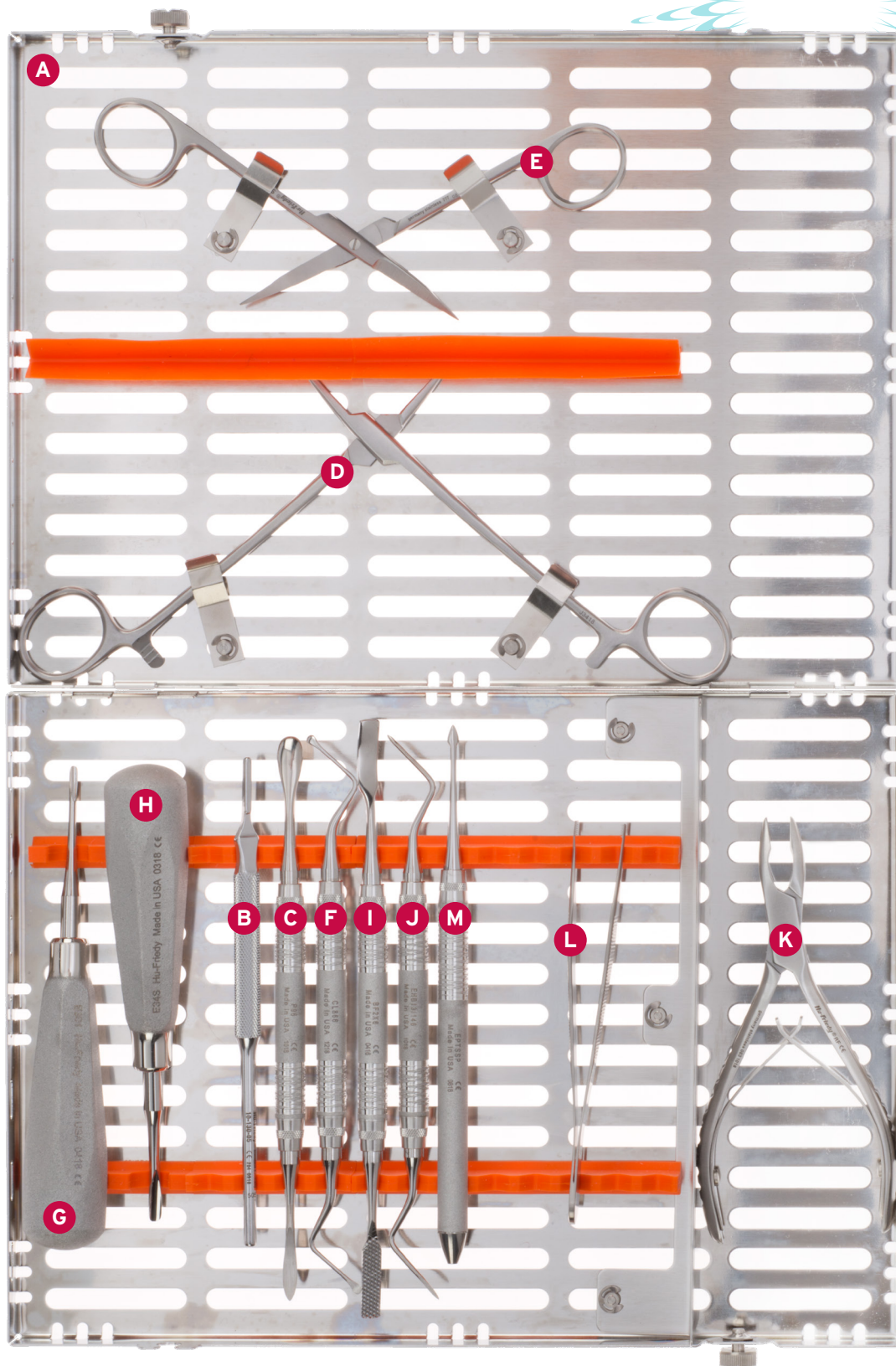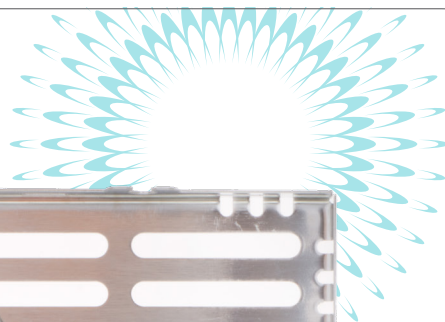

EXTRACTION KIT



A **IM4163E2C4**
16 Instrument Signature Series® Cassette, Orange

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B **10-130-05**
5 Scalpel Handle, Round Straight

Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



C **P96**
9 Molt Periosteal Elevator

Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D **NHMH**
Mayo-Hegar Needle Holder

Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



E **S16**
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



F **CL866**
86 Lucas DE Surgical Curette, #6 Handle

Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



G **E301**
301 Apexo Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction.



H **E34S**
34S Seldin Straight Elevator

Used to loosen the tooth from the periodontal ligament and ease extraction. 4.0mm



I **BF2X6**
2X Miller-Colburn Bone File

Used for final smoothing and shaping of bone with 6.2mm wide working ends.



J **EHB13/146**
13/14 Heidbrink DE Root Pick

Used for the removal of root tips in difficult-to-reach areas. Designed for delicate root teasing procedures. Also available in a "heavy" pattern for more substantial root tip removal procedures.



K **RF**
Friedman Rongeurs

Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



L **TP43**
43 Adson-Brown Tissue Pliers

Used to grasp and stabilize soft tissue during suturing. 12mm



M **EPTSSP**
Luxating Hybrid Elevator, Straight Spade

A unique cross between a periosteal and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.