



HU-FRIEDY SURGICAL KIT

LATERAL SINUS LIFT KIT



PART CODE	DESCRIPTION
IM4162	16 Instrument Signature Series® Cassette, White
IMP65786	Kramer-Nevins Sinus Lift Instr., Obtuse Angles, #6 Handle
IMP6577SC6	Kramer-Nevins Sinus Lift Instr., Acute Angles, #6 Handle
BPACK	Bone Packer DE 6-8mm
SYBG	Syringe Bone Grafting 7mm
BONT	Bone Tamper 7mm
NHCW	Crile-Wood Needle Holder
S16	16 Goldman-Fox Scissors
CRURBAN	Urban Cheek Retractor
IM1002	6-Pack Signature Series® Instrument Clip
IMLATKIT	Complete Lateral Sinus Lift Kit

LATERAL SINUS LIFT KIT

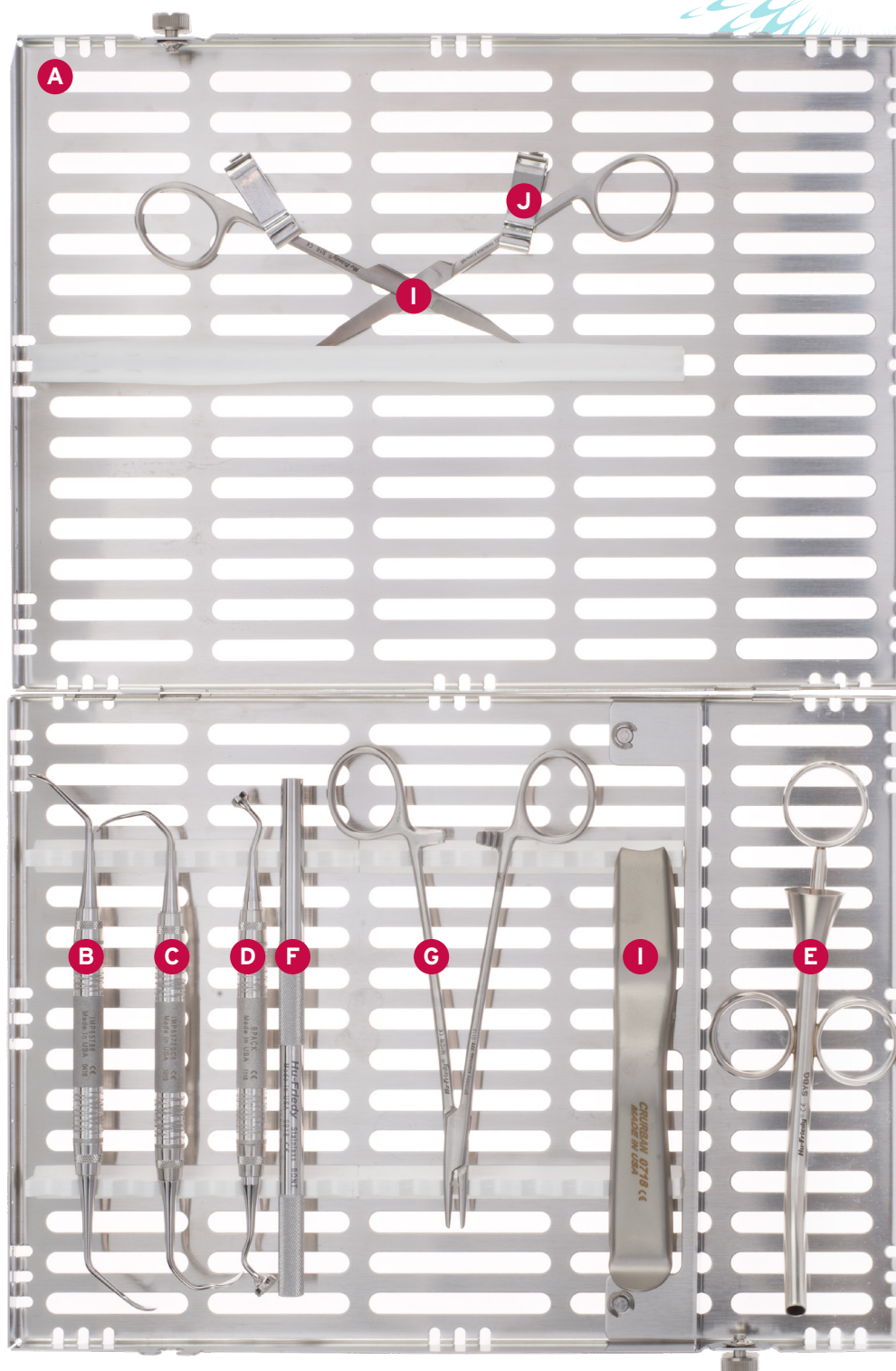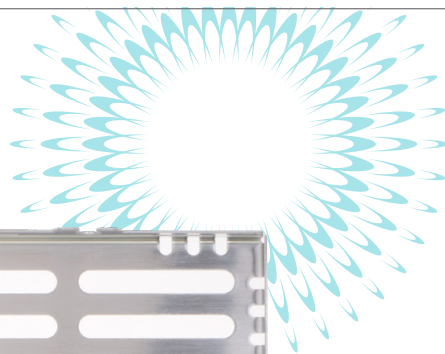

LATERAL SINUS LIFT KIT



A **IM4162**
16 Instrument Signature Series® Cassette, White

IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B **IMP65786**
Kramer-Nevins Sinus Lift Instr., Obtuse Angles, #6 Handle

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With obtuse angles 4.1mm/7.6mm.



C **IMP6577SC6**
Kramer-Nevins Sinus Lift Instr., Acute Angles, #6 Handle

Used during a sinus lift procedure to separate/reflect the Schneiderian membrane from the maxillary bone and to elevate the membrane. With acute angles 4.1mm/7.6mm.



D **BPACK**
Bone Packer DE 6-8mm

Used to condense bone during grafting procedures with working ends 6mm and 8mm wide.



E **SYBG**
Syringe Bone Grafting 7mm

Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.



F **BONT**
Bone Tamper 7mm

Used to separate the Schneiderian membrane from surrounding bone during lateral sinus lift procedures.



G **NHCW**
Crile-Wood Needle Holder

Used to separate the Schneiderian membrane from surrounding bone during lateral sinus lift procedures.




H **S16**
16 Goldman-Fox Scissors

Used for trimming tissue or cutting sutures. 12.5cm with 1 serrated blade.



I **CRURBAN**
Urban Cheek Retractor

Used to hold mucoperiosteal flaps, cheeks, lips and tongue away from the surgical area. 20.6mm

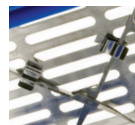

J **IM1002**
6-Pack Signature Series® Instrument Clip

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.