



HU-FRIEDY SURGICAL KIT

BONE GRAFTING KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4128	12 Instrument Infinity Series™ Cassette, Blue
10-130-05E	5 Scalpel Handle
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
SPM20	Micro Suture Pliers
PPBUSERX	Buser Periosteal Elevator, Black Line
P9X	9 Molt Periosteal Elevator, Black Line
PPBUSER13X	PJO/1 & PJO/3 Combined, 13KL, Black Line
PPBUSER23X	PJO/2 & PJO/3 Combined, 13KL, Black Line
CL86X	86 Lucas DE Surgical Curette, Black Line
SPR1/26	PR1/2 Prichard Curette, #6 Handle
MCUPE	Measuring Cup 25ml with Markings
S16SCX	16 Goldman-Fox Scissors, Black Line
NHCW	Crile-Wood Needle Holder
PMP1X	Membrane Placement Instrument, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMNBGKITX	Complete Bone Grafting Kit, Black Line

BONE GRAFTING KIT, BLACK LINE

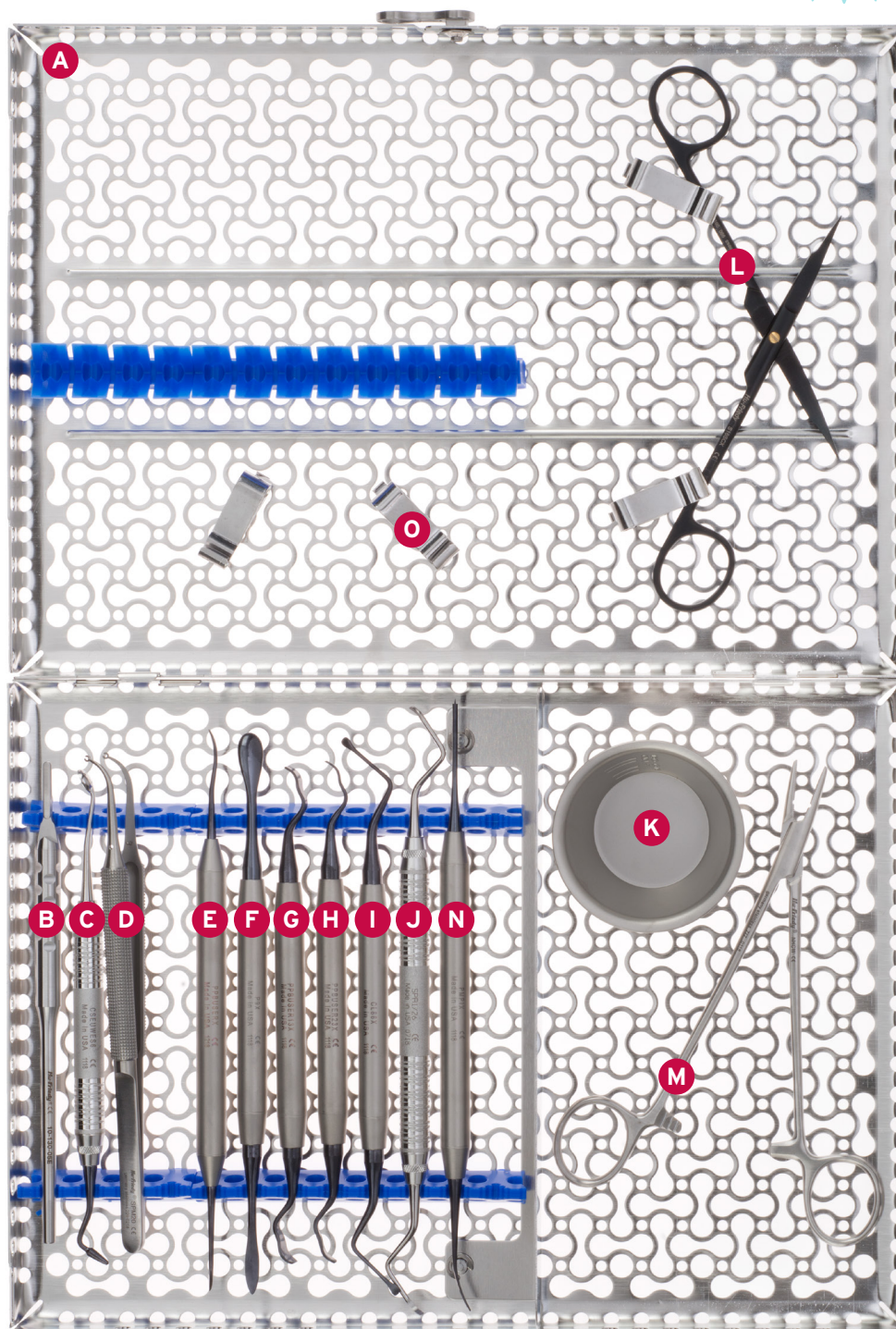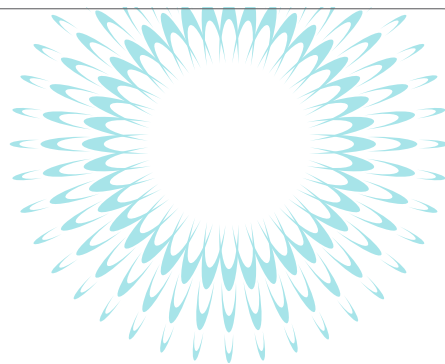

BONE GRAFTING KIT, BLACK LINE



- A IM4128**
12 Instrument Infinity Series™ Cassette, Blue
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- J SPR1/26**
PR1/2 Prichard Curette, #6 Handle
 Larger and heavier curette for the removal of granulation tissue and tenacious subgingival deposits measuring 7.5" long with 1.5mm wide working ends.



- K MCUPE**
Measuring Cup 25ml with Marking




- L S16SCX**
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- M NHCW**
Crile-Wood Needle Holder
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



- N PMPIX**
Membrane Placement Instrument, Black Line
 Micropoint tip for precise placement and manipulation of membranes.

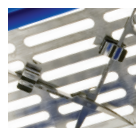

- O IM1002**
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



- 10-130-05E**
B 5 Scalpel Handle
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- CSEUWES6**
C Euwe Bone Scoop/Plugger, #6 Handle
 Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



- SPM20**
D Micro Suture Pliers
 Ideal for use with guided tissue regeneration procedures. 15cm



- PPBUSERX**
E Buser Periosteal Elevator, Black Line
 1/3 Buser Modified Periodontal Chisel 4mm/6mm.



- F P9X**
9 Molt Periosteal Elevator, Black Line
 Featuring a wide sharp blade, non-reflective surface and improved contrast between tissue and elevator. Used for reflecting and retracting the mucoperiosteum.



- G PPBUSER13X**
PJO/1 & PJO/3 Combined 13KL, Black Line
 1/3 Buser Modified Periodontal Chisel 4mm/6mm.



- H PPBUSER23X**
PJO/2 & PJO/3 Combined 13KL, Black Line
 2/3 Buser Modified Periodontal Chisel 5mm/6mm.



- I CL86X**
86 Lucas DE Surgical Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm