



HU-FRIEDY SURGICAL KIT

CRESTAL SINUS LIFT KIT, BLACK LINE



PART CODE	DESCRIPTION
IMN4162	16 Instrument Infinity Series™ Cassette, White
10-130-03	3 Scalpel Handle
P9X	9 Molt Periosteal Elevator, Black Line
CM2/4X	2/4 Molt DE Curette, Black Line
TP5043	Adson-Brown Perma Sharp® Tissue Pliers
CSEUWES6	Euwe Bone Scoop/Plugger, #6 Handle
PCP65006	20mm Probe 8-10-12-14-16-18-20, #6 Handle
OSTPAL1	1 Palti Osteotome 2.7mm
OSTPAL2	2 Palti Osteotome 3.2mm
OSTPAL3	3 Palti Osteotome 3.7mm
OSTPAL4	4 Palti Osteotome 4.2mm
PFIWDS16	1 Woodson Plastic Filling, #6 Handle
IM1002	6-Pack Signature Series® Instrument Clip
NH5038	Crile-Wood Perma Sharp® Needle Holder
IMNCREKITX	Complete Crestal Sinus Lift Kit, Black Line

CRESTAL SINUS LIFT KIT, BLACK LINE

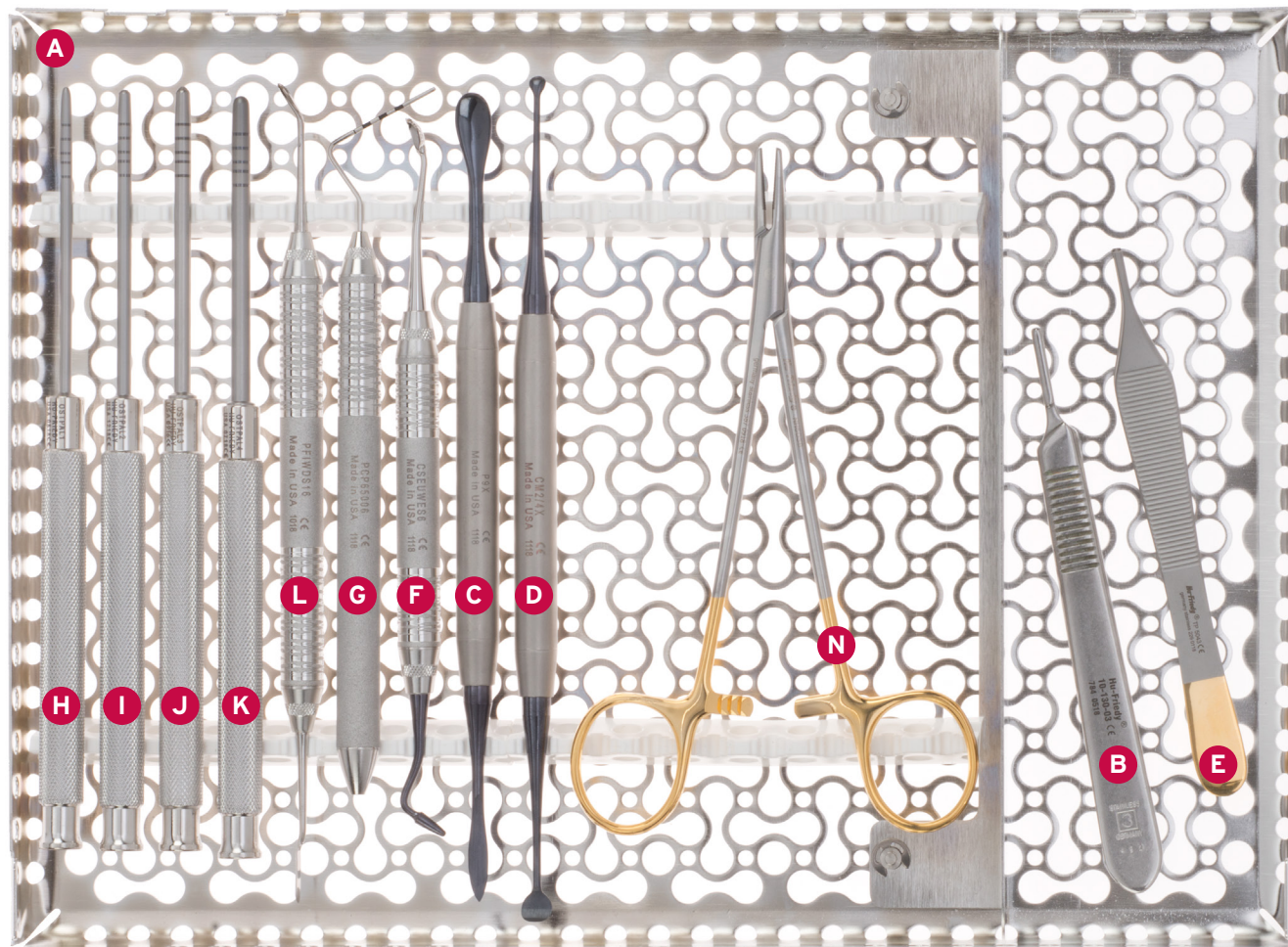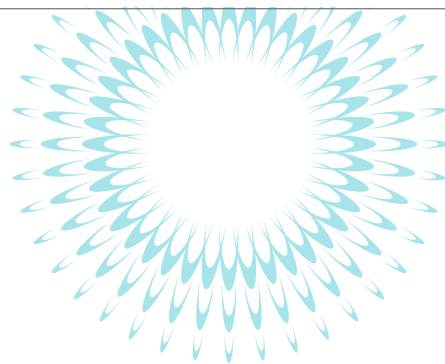

CRESTAL SINUS LIFT KIT, BLACK LINE



A IMN4162
16 Instrument Infinity Series™ Cassette, White
IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-03
3 Scalpel Handle
Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



C P9X
9 Molt Periosteal Elevator, Black Line
Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D CM2/4X
2/4 Molt DE Curette, Black Line
Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



E TP5043
Adson-Brown Perma Sharp® Tissue Pliers
Used to grasp and stabilize soft tissue during suturing. 12cm



F CSEUWES6
Euwe Bone Scoop/Plugger, #6 Handle
Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



G PCP65006
20mm Probe 8-10-12-14-16-18-20, #6 Handle
A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



H OSTPAL1
1 Palti Osteotome 2.7mm
Used for preparing bone. 7-10-13-15-18mm



I OSTPAL2
2 Palti Osteotome 3.2mm
Used for preparing bone. 7-10-13-15-18mm




J OSTPAL3
3 Palti Osteotome 3.7mm
Used for preparing bone. 7-10-13-15-18mm



K OSTPAL4
4 Palti Osteotome 4.2mm
Used for preparing bone. 7-10-13-15-18mm



L PFIWDS16
1 Woodson Plastic Filling, #6 Handle
Designed to reflect and retract the mucoperiosteum after incisions of the gingival tissue with 2.2mm and 3.8mm width.

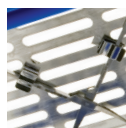

M IM1002
6-Pack Signature Series® Instrument Clip
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



N NH5038
Crile-Wood Perma Sharp® Needle Holder
Used to grasp and guide the needle when suturing. Featuring 3-0 to 5-0 carbide beaks, finger ring style.