



HU-FRIEDY SURGICAL KIT

EXTRACTION KIT, BLACK LINE



PART CODE	DESCRIPTION
IMNOS3	Oral Surgery Infinity Series™ Cassette
10-130-05E	Scalpel Handle #5
P9X	9 Molt Periosteal Elevator, Black Line
NMH	Mayo-Hegar Needle Holder
S16SCX	16 Goldman-Fox Scissors, Black Line
CL86X	86 Lucas DE Surgical Curette, Black Line
E301	301 Apexo Elevator
E21	Mini Cryer Left Elevator
E22	Mini Cryer Right Elevator
BF2X6	2X Miller-Colburn Bone File
EHB13/14X	13/14 Heidbrink DE Root Pick, Black Line
RF	Friedman Rongeurs
TP43	43 Adson-Brown Tissue Pliers
EPTSSPX	Luxating Hybrid Elevator, Straight Spade, Black Line
IM1002	6-Pack Signature Series® Instrument Clip
IMNEXTKITX	Complete Extraction Kit, Black Line

EXTRACTION KIT, BLACK LINE

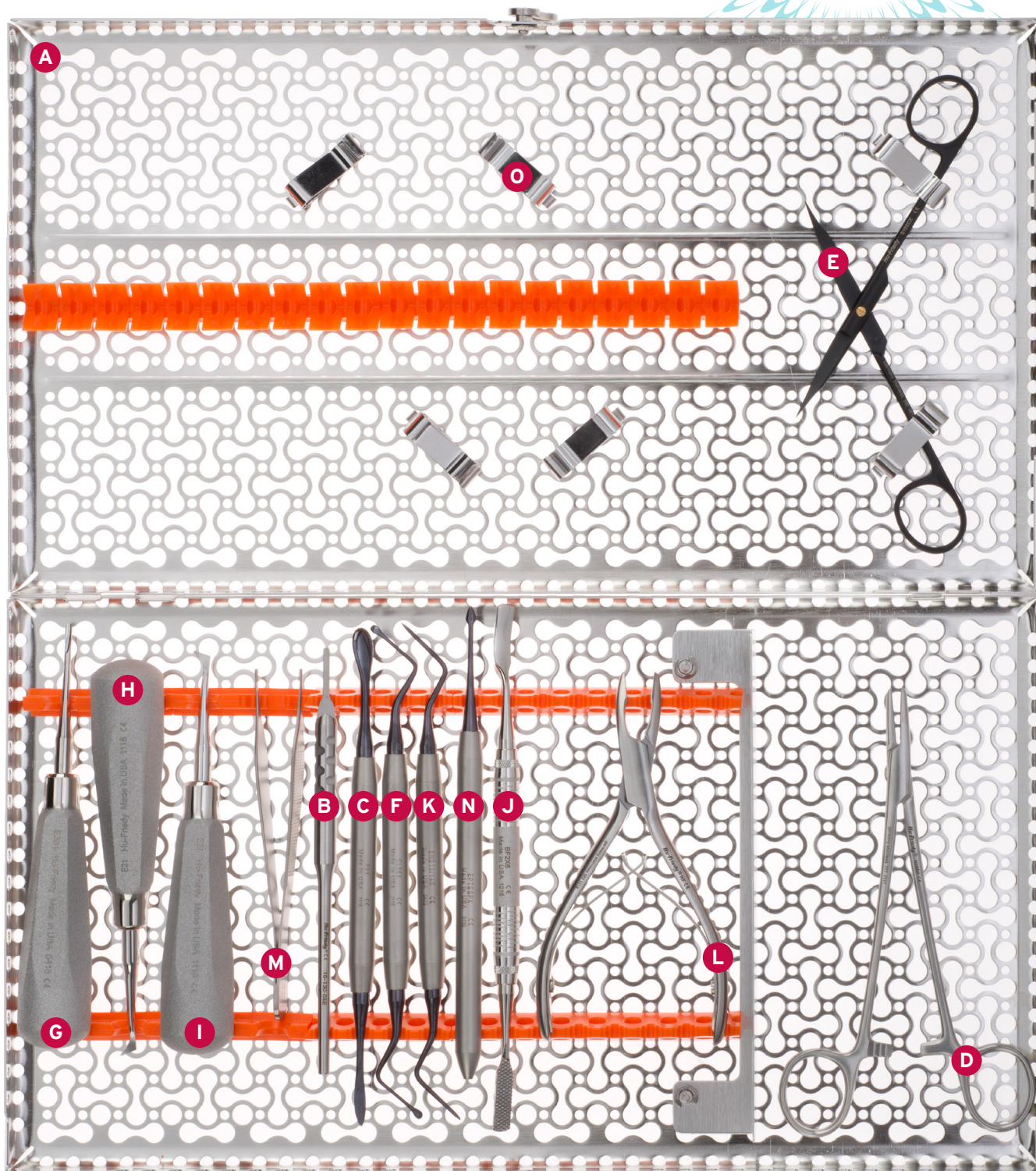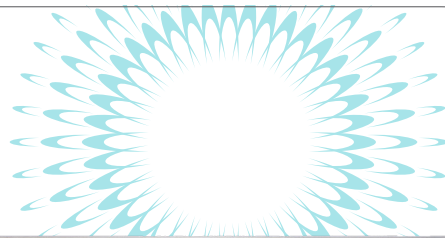

EXTRACTION KIT, BLACK LINE



A IMNOS3
Oral Surgery Infinity Series™ Cassette, Orange
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



B 10-130-05E
Scalpel Handle #5
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



C P9X
9 Molt Periosteal Elevator, Black Line
 Designed for reflecting and retracting the muco-periosteum after incisions of the gingival tissue. 4.3mm/8.4mm



D NHMH
Mayo-Hegar Needle Holder
 Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



E S16SCX
16 Goldman-Fox Scissors, Black Line
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



F CL86X
86 Lucas DE Surgical Curette, Black Line
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



G E301
PJO/2 & PJO/3 Combined 13KL, #6 Handle
 Used to loosen the tooth from the periodontal ligament and ease extraction.



H E21
Mini Cryer Left Elevator
 Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



I E22
Mini Cryer Right Elevator
 Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



J BF2X6
2X Miller-Colburn Bone File
 Used for final smoothing and shaping of bone with 6.2mm wide working ends.



K EHB13/14X
13/14 Heidbrink DE Root Pick, Black Line
 For the removal of root tips in difficult-to-reach areas.



L RF
Friedman Rongeurs
 Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



M TP43
43 Adson-Brown Tissue Pliers
 Used to grasp and stabilize soft tissue during suturing. 12mm



N EPTSSPX
Luxating Hybrid Elevator, Straight Spade, Black Line
 A unique cross between a periosteal and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



O IM1002
6-Pack Signature Series® Instrument Clip
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.