




POWER PLUS INSERTS

REACH FURTHER FOR POCKET DEBRIDEMENT

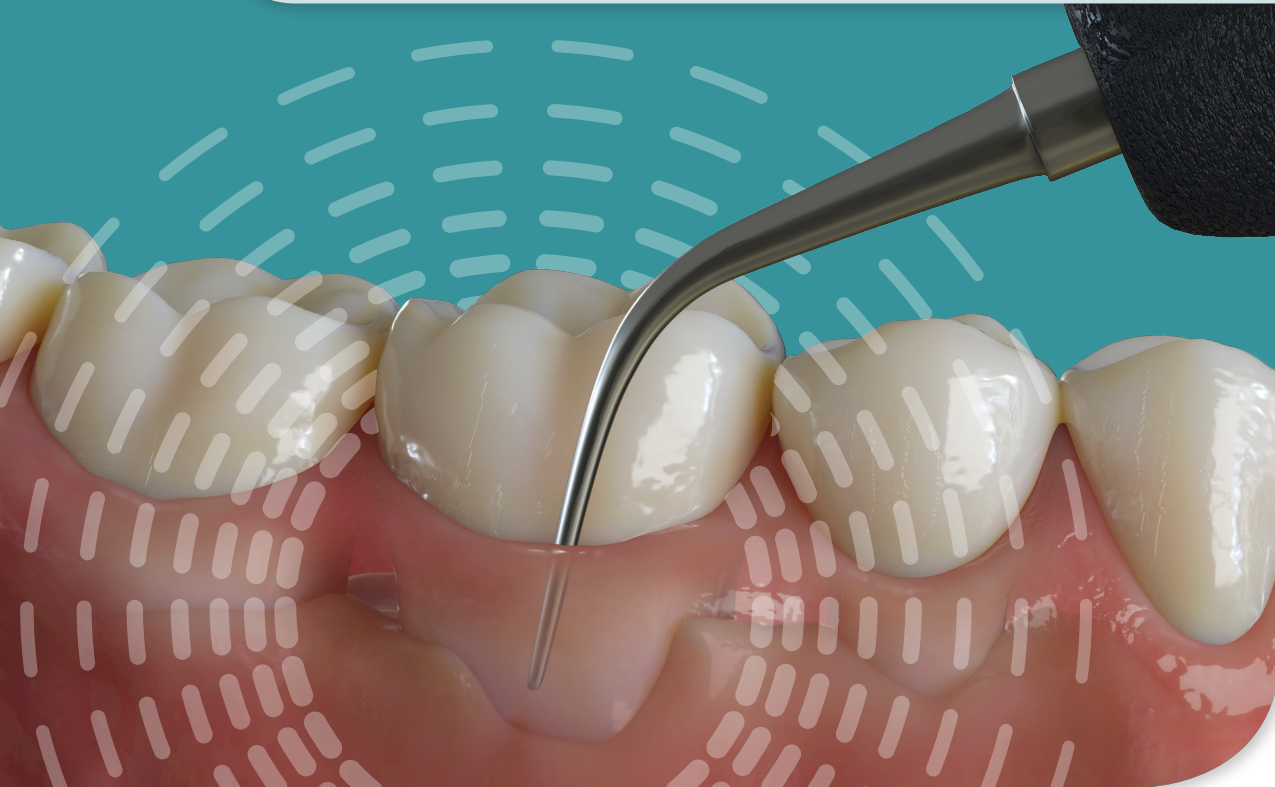

New **POWER PLUS** inserts from Hu-Friedy feature **POWER** for effective removal of tenacious deposits and heavy extrinsic stain **PLUS** extended length to reach deep into periodontal pockets. With a 25% longer tip design,* a variety of tip sizes, and a focus on clinician comfort, Power PLUS inserts offer an advanced solution for top-notch clinical outcomes.

Four Tip Options for improved access to pocket anatomy, tooth anatomy and debris type allowing for optimal debris removal

Conical Tips: Thin and Standard Size

- > Round tip cross-section
- > Allows energy to spread uniformly around the tip
- > Ideal for lighter deposits, new calculus, biofilm removal

Beveled Tips: Thin and Standard Size

- > Diamond shape tip cross section
- > Energy dispersion concentrated on bevels or edges
- > Bevels allow for power boost to fracture mature/tenacious deposits

How the best perform

Hu-Friedy

POWER PLUS INSERTS

Special blend of stainless steel for tip longevity and durability









Popular ergonomic Swivel handle for excellent tip control and maneuverability

Color-coded grip for at-a-glance tip identification

2mm longer tip* for extended reach into periodontal pockets

Single Bend Shank Similar to a periodontal probe provides familiarity, ease of use



Name	Tip Description	Tip Cross-Section	Recommendation	Part Code
Power PLUS Standard Conical	Round Cross-section. Standard diameter shank extending to an extended length tapered tip.	Round 	Moderate to heavy deposits; supragingival; wide, deep pockets	UI30PPSC 
Power PLUS Thin Conical	Round Cross-section. Thin diameter shank extending to an extended length tapered tip.	Round 	Light to moderate deposits; narrow, deep vertical pockets; areas mesial to retromolar pad; crowded areas; general supra use	UI30PPTC 
Power PLUS Standard Beveled	Diamond Cross-section. Standard diameter shank extending to an extended length tapered tip.	Beveled 	Moderate to heavy deposits; tenacious deposits; supragingival; wide, deep pockets	UI30PPSB 
Power PLUS Thin Beveled	Diamond Cross-section. Thin diameter shank extending to an extended length tapered tip.	Beveled 	Light to moderate deposits; tenacious deposits; narrow, deep vertical pockets; crowded areas	UI30PPTB 

Learn more at [Hu-Friedy.com/Inserts](https://www.hu-friedy.com/Inserts)

How the best perform

