

# Surgical Kit

## Atraumatic Extraction Kit (IMATRAEXT)



- 1 IM4161**  
**16 Instrument Signature Series™ Cassette, Gray**  
IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2 EPTSMS**  
**Luxating Hybrid Elevator, Small Straight**  
A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- 3 EPTSMC**  
**Luxating Hybrid, Small Curved**  
A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- 4 EPTSSP**  
**Luxating Hybrid Elevator, Straight Spade**  
A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.



- 5 FAF150**  
**Apical Forceps 150**  
Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- 6 FAF151**  
**Apical Forceps 151**  
Thin, tapered and sharp parallel beaks improve root and crown contact and longitudinal serrations provide additional grip. The design of Apical Forceps facilitates an atraumatic tooth extraction.



- 7 CM2/4S6**  
**2/4 Molt DE Curette, #6 Handle**  
Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



- 8 NHCW**  
**Crile-Wood Needle Holder**  
Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0



- 9 S16SC**  
**16 Goldman-Fox Scissors, Super Cut**  
Used for trimming tissue or cutting sutures. Featuring tapered blade tips, tips that are not sharp on outside and are used as tonotomy scissors. 12.5cm with 1 serrated blade



- 10 10-130-05**  
**5 Scalpel Handle, Round Straight**  
Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.




- 11 P96**  
**9 Molt Periosteal Elevator**  
Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



- IM1002**  
**6-Pack Hinged Instrument Clip**  
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

## Atraumatic Extraction Kit (IMATRAEXT)

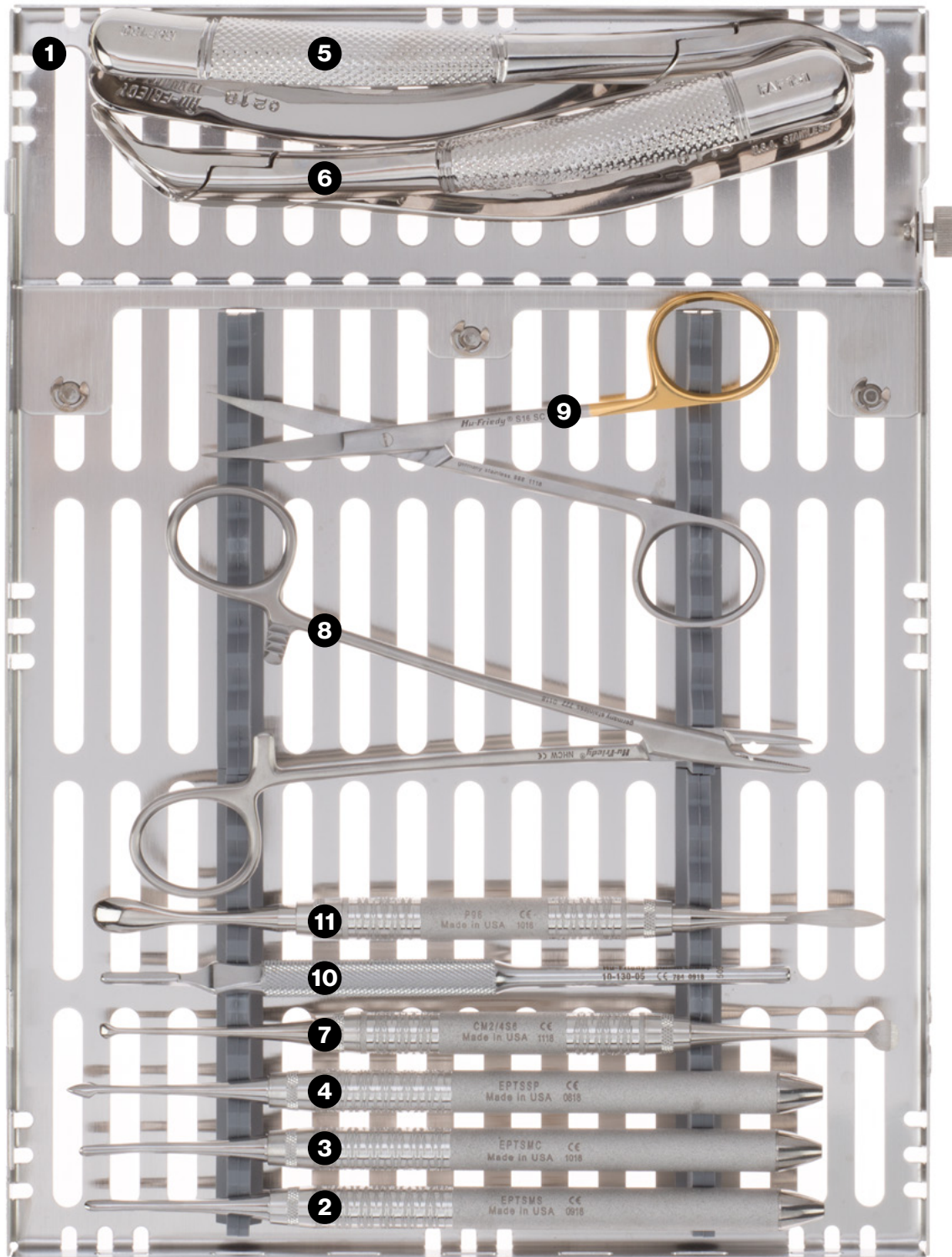

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.  
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-335/1023