

# Surgical Kit

## Bone Grafting Kit, Black Line (IMNBGKITX)



- 1 IMN4128**  
**12 Instrument Infinity Series™ Cassette, Blue\***  
 IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2 10-130-05E**  
**5 Scalpel Handle**  
 Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- 3 CSEUWES6**  
**Euwe Bone Scoop/Plugger, #6 Handle**  
 Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm



- 4 SPM20**  
**Micro Suture Pliers**  
 Ideal for use with guided tissue regeneration procedures. 15cm



- 5 PPBUSERX**  
**Buser Periosteal Elevator, Black Line**  
 1/3 Buser Modified Periodontal Chisel 4mm/6mm



- 6 P9X**  
**9 Molt Periosteal Elevator, Black Line**  
 Featuring a wide sharp blade, non-reflective surface and improved contrast between tissue and elevator. Used for reflecting and retracting the mucoperiosteum.



- 7 PPBUSER13X**  
**PJO/1 & PJO/3 Combined 13KL, Black Line**  
 1/3 Buser Modified Periodontal Chisel 4mm/6mm



- 8 PPBUSER23X**  
**PJO/2 & PJO/3 Combined 13KL, Black Line**  
 2/3 Buser Modified Periodontal Chisel. 5mm/6mm



- 9 CL86X**  
**86 Lucas DE Surgical Curette, Black Line**  
 Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



- 10 SPR1/26**  
**PR1/2 Prichard Curette, #6 Handle**  
 Larger and heavier curette for the removal of granulation tissue and tenacious subgingival deposits measuring 7.5" long with 1.5mm wide working ends.



- 11 PMPIX**  
**Membrane Placement Instrument, Black Line**  
 Micropoint tip for precise placement and manipulation of membranes.



- 12 MCUPE**  
**Measuring Cup 25ml with Marking**



- 13 NHCW**  
**Crile-Wood Needle Holder**  
 Used to grasp and guide the needle when suturing. Recommended suture sizes: 3-0, 4-0, 5-0, 6-0




- 14 S16SCX**  
**16 Goldman-Fox Scissors, Black Line**  
 Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- 15 IM1002**  
**6-Pack Hinged Instrument Clip**  
 Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

\*Pictured in IM4128 - 12 Instrument Signature Series™ Cassette, Blue

**Bone Grafting Kit, Black Line (IMNBGKITX)**

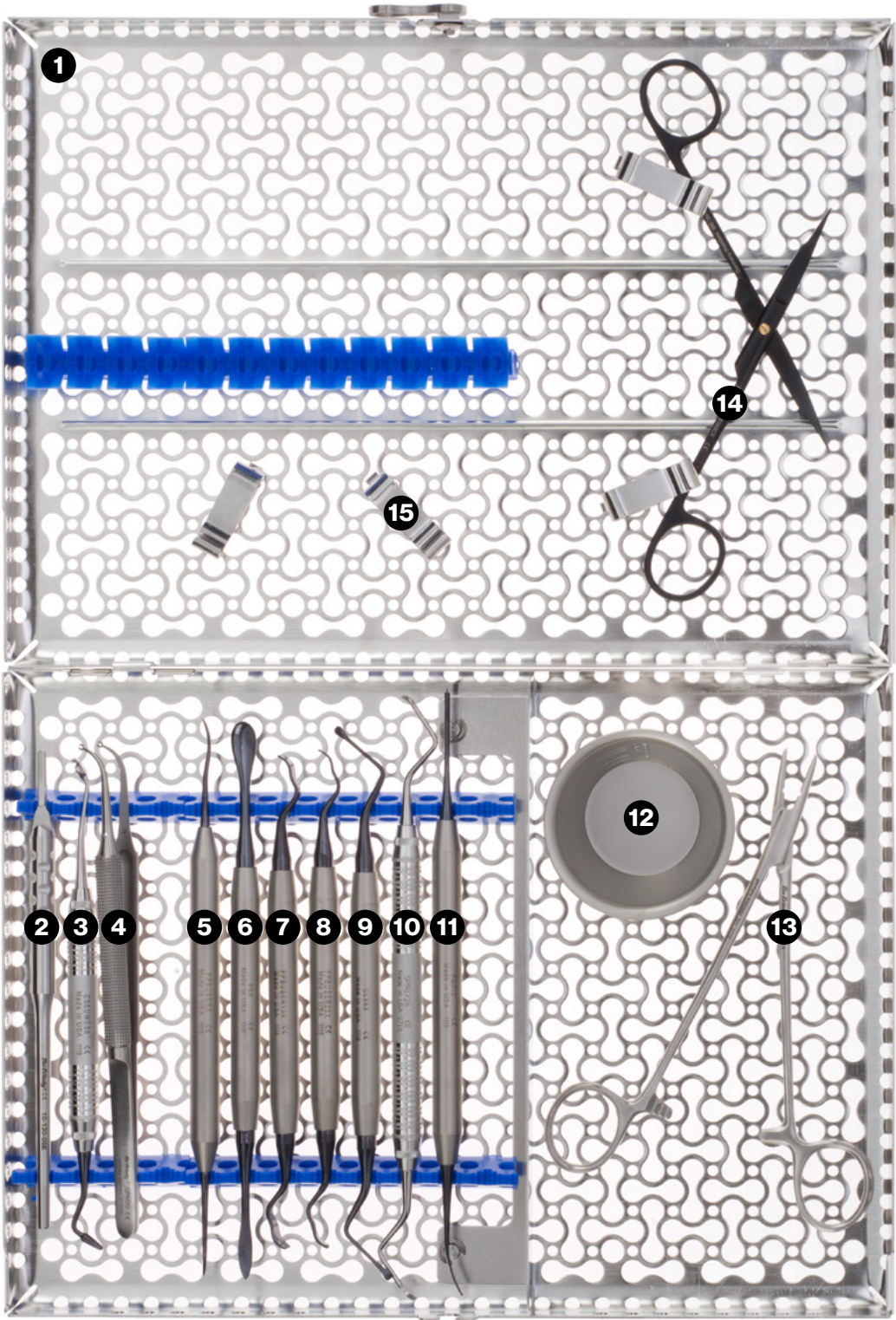

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | [HuFriedyGroup.com](http://HuFriedyGroup.com)  
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.  
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-344 /1023