

# Surgical Kit

## Crestal Sinus Lift Kit, Black Line (IMNCREKITX)



**1** **IMN4162**  
**16 Instrument Infinity Series™**  
**Cassette, White**

IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



**2** **10-130-03**  
**3 Scalpel Handle**

Used for cutting gingival tissue and making surgical incisions. Handle is 4.9" long and 13mm wide at base. Compatible with standard scalpel blade sizes.



**3** **P9X**  
**9 Molt Periosteal Elevator, Black Line**

Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



**4** **CM2/4X**  
**2/4 Molt DE Curette, Black Line**

Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



**5** **TP5043**  
**Adson-Brown Perma Sharp™**  
**Tissue Pliers**

Used to grasp and stabilize soft tissue during suturing. 12cm



**6** **CSEUWES6**  
**Euwe Bone Scoop/Plugger, #6 Handle**

Designed for easy transport and non-stick compaction of bone graft material. 3.0mm/7.0mm.



**7** **PCP65006**  
**20mm Probe 8-10-12-14-16-18-20,**  
**#6 Handle**

A double-ended instrument with an explorer on one end and a probe on the other, excellent for dental examinations and charting.



**8** **OSTPAL1**  
**1 Palti Osteotome 2.7mm**

Used for preparing bone. 7-10-13-15-18mm



**9** **OSTPAL2**  
**2 Palti Osteotome 3.2mm**

Used for preparing bone. 7-10-13-15-18mm



**10** **OSTPAL3**  
**3 Palti Osteotome 3.7mm**

Used for preparing bone. 7-10-13-15-18mm



**11** **OSTPAL4**  
**4 Palti Osteotome 4.2mm**

Used for preparing bone. 7-10-13-15-18mm



**12** **PFIWDS16**  
**1 Woodson Plastic Filling, #6 Handle**


Designed to reflect and retract the mucoperiosteum after incisions of the gingival tissue with 2.2mm and 3.8mm width.



**13** **NH5038**  
**Crile-Wood Perma Sharp™ Needle Holder**

Used to grasp and guide the needle when suturing. Featuring 3-0 to 5-0 carbide beaks, finger ring style.

## Crestal Sinus Lift Kit, Black Line (IMNCREKITX)

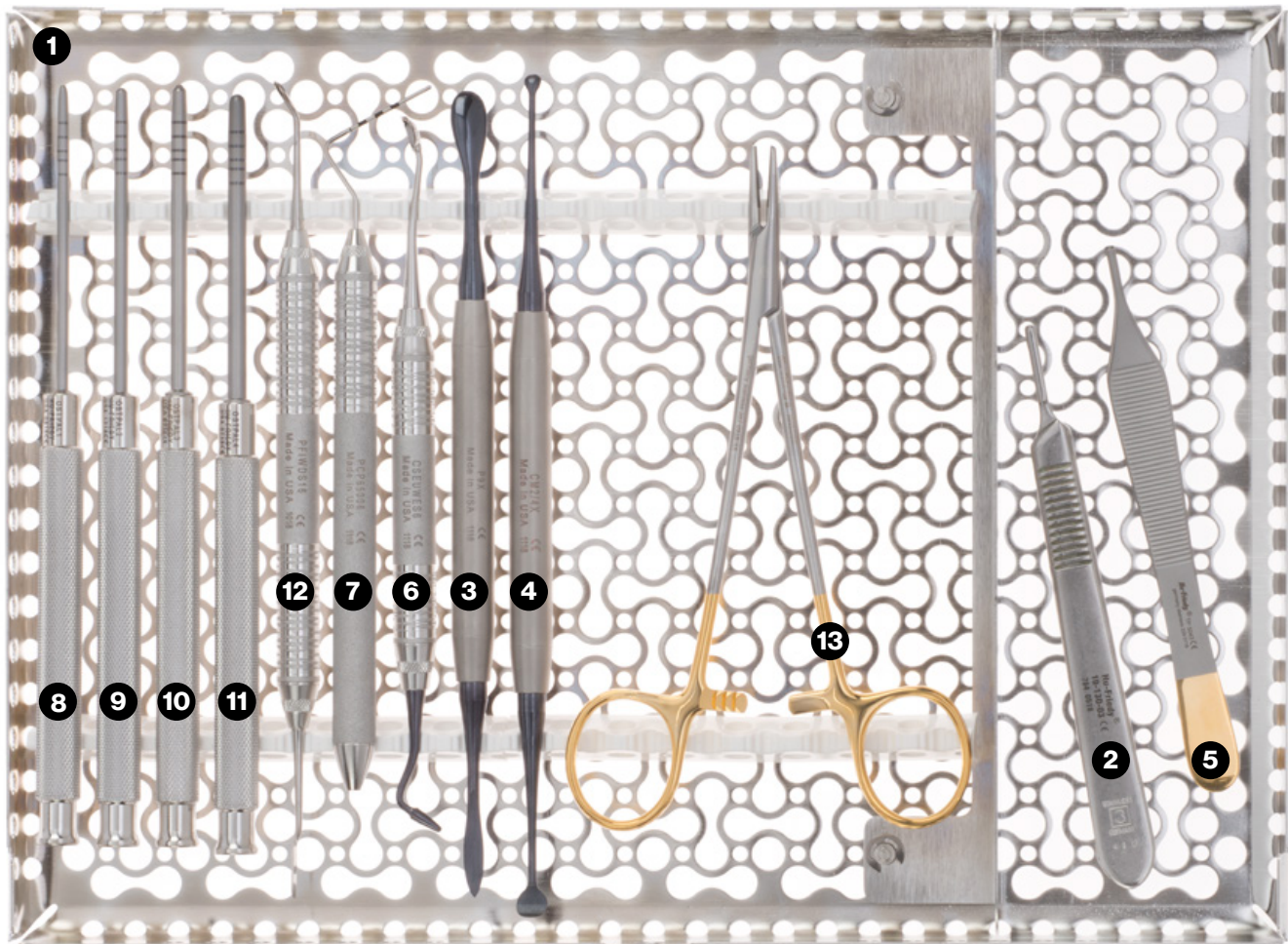

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com  
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.  
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-3458 /1023