

Surgical Kit

Extraction Kit, Black Line (IMNEXTKITX)



- 1** **IMNOS3**
Oral Surgery Infinity Series™ Cassette, Orange
IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- 2** **10-130-05E**
Scalpel Handle #5
Used for cutting gingival tissue and making surgical incisions. Provides greater rotational movement with fingertip control and enhances dexterity. Handle measures 6" long. Compatible with standard scalpel blade sizes.



- 3** **P9X**
9 Molt Periosteal Elevator, Black Line
Designed for reflecting and retracting the mucoperiosteum after incisions of the gingival tissue. 4.3mm/8.4mm



- 4** **NMH**
Mayo-Hegar Needle Holder
Used to grasp and guide the needle when suturing. 16 cm. Recommended suture sizes: 1-0, 2-0, 3-0, 4-0.



- 5** **S16SCX**
16 Goldman-Fox Scissors, Black Line
Used for trimming tissue or cutting sutures. Featuring a non-reflective surface, improved contrast between tissue/scissor. 13.5cm.



- 6** **CL86X**
86 Lucas DE Surgical Curette, Black Line
Designed for curettage, cyst removal and tooth socket debridement. 2.8mm



- 7** **E301**
PJO/2 & PJO/3 Combined 13KL, #6 Handle
Used to loosen the tooth from the periodontal ligament and ease extraction.



- 8** **E21**
Mini Cryer Left Elevator
Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



- 9** **E22**
Mini Cryer Right Elevator
Used to loosen the tooth from the periodontal ligament and ease extraction, 6.5mm wide tip.



- 10** **BF2X6**
2X Miller-Colburn Bone File
Used for final smoothing and shaping of bone with 6.2mm wide working ends.



- 11** **EHB13/14X**
13/14 Heidbrink DE Root Pick, Black Line
For the removal of root tips in difficult-to-reach areas.




- 12** **RF**
Friedman Rongeurs
Used for trimming and recontouring alveolar bone and gross tissue removal. 2.4mm/14cm



- 13** **TP43**
43 Adson-Brown Tissue Pliers
Used to grasp and stabilize soft tissue during suturing. 12mm




- 14** **EPTSSPX**
Luxating Hybrid Elevator, Straight Spade, Black Line
A unique cross between a periosteal elevator and a luxating elevator, used apically within the periodontal ligament space to cut periodontal ligaments prior to extraction.

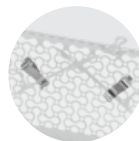

- 15** **IM1002**
6-Pack Hinged Instrument Clip
Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

Extraction Kit, Black Line (IMNEXTKITX)

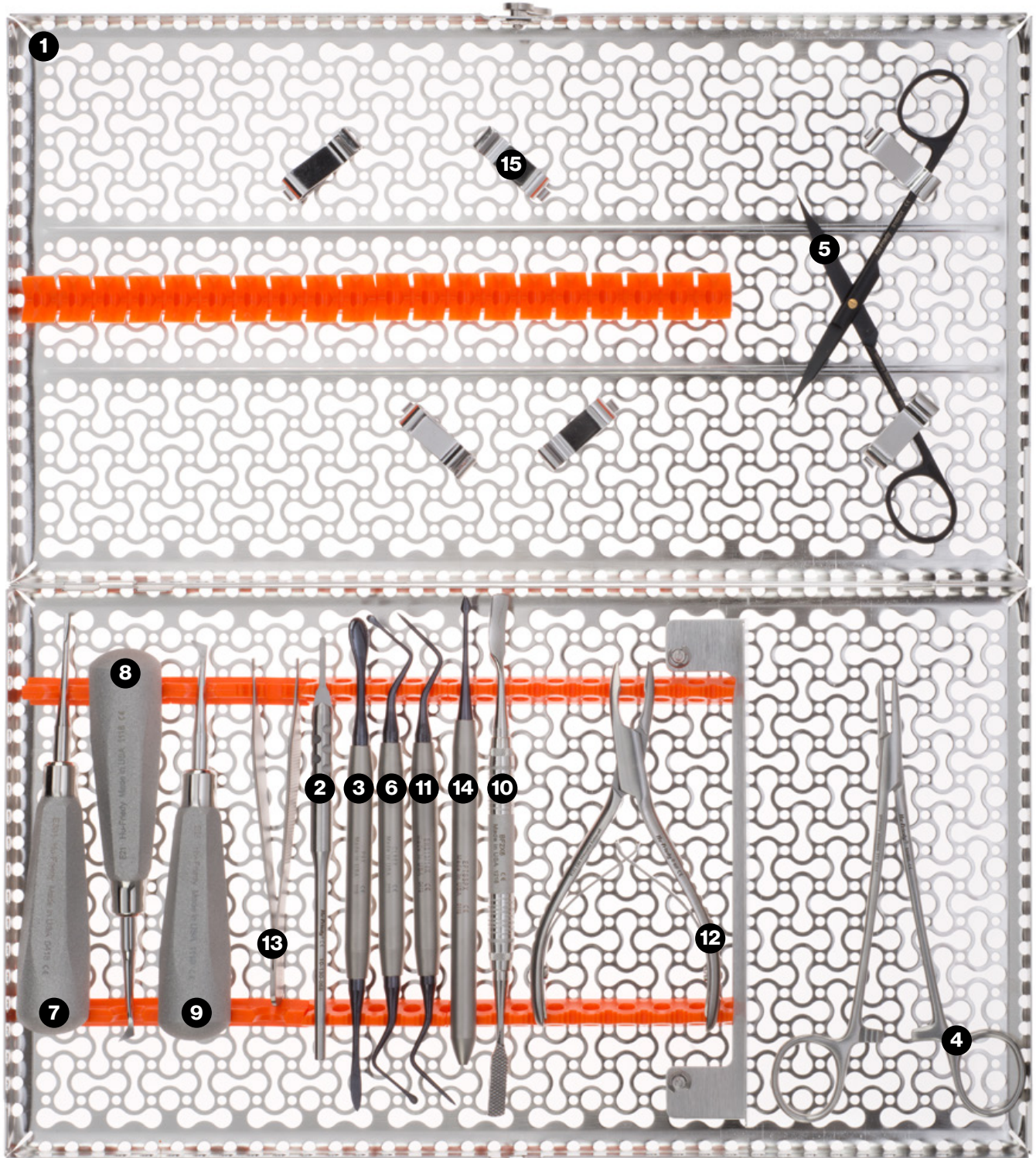

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | HuFriedyGroup.com
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-346/1023