

# Surgical Kit

## Keys Socket Grafting Kit (IMNKEYSGKT)



**1 IMN4164**  
**16 Instrument Infinity Series™ Cassette, Purple**

IMS™ Cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



**2 XP23/86**  
**23/CP-8 Expro**

Double-ended instrument with an explorer on one end and a probe (color-coded) on the other. ergonomic, Satin Steel™ Handle.



**3 H2**  
**2 Kelly Hemostat, Curved**

Used to clamp off blood vessels, remove small root tips and grasp loose objects. Featuring kelly curved beaks.



**4 PLGGR1X**  
**Grisdale Bone Condenser 8-10-13-15-18, Black Line**

Used during grafting procedures with 8-11-13-15-18 mm markings. Featuring a non-reflective surface, improved contrast between tissue/scissor.



**5 CM2/4SX**  
**2/4 Molt DE Curette, Black Line**

Designed for curettage, cyst removal and tooth socket debridement with 3.6mm and 7mm working ends.



**6 10-130-05**  
**5 Scalpel Handle, Round, Straight**

Used for cutting gingival tissue and making surgical incisions.



**7 PPR3X**  
**PR-3 Prichard Periosteal, Black Line**

Used for reflecting and retracting the mucoperiosteum. Featuring a large flat blade and narrow sharp blade.



**8 NH5037**  
**Halsey Perma Sharp™ Needle Holder, Serrated**

Used to grasp and guide the needle when suturing. Featuring carbide inserts, finger ring style. Use on all types of sutures.



**9 S5055**  
**Metzenbaum Perma Sharp™ Scissors, Curved/Blunt**

Used for trimming tissue or cutting sutures. 14.5cm



**10 TP5090**  
**DeBakery Perma Sharp™ Tissue Pliers, 6"**

Used to grasp and stabilize soft tissue during suturing.



**11 DP1**  
**#1 Dressing Pliers, Perry**

Used to grip or transfer materials in and out of the oral cavity. Features a serrated tip to provides a positive grip.



**12 PMPIX**  
**Membrane Placement Instrument, Black Line**

Micropoint tip for precise placement and manipulation of membranes.



**13 MCUPE**  
**Measuring Cup 25 mL w/ Marking**

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.



**14 SYBG**  
**Syringe Bone Grafting 7mm**

Used during grafting procedures. Inner Diameter 5.8 mm. Load from tip to prevent clogging.




**15 CL856**  
**85 Lucas DE Surgical Curette**

Designed for curettage, cyst removal and tooth socket debridement. 2.5mm



**16 PJO/16**  
**S13K/TG Jovanovic 4mm Chisel**


For removing and recontouring bone.

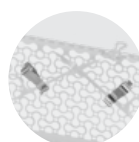

**17 IM1002**  
**6-Pack Hinged Instrument Clip**

Secure and maintain open position of scissors, hemostats, and other hinged instruments in lid of cassette.

## Keys Socket Grafting Kit (IMNKEYSGKT)

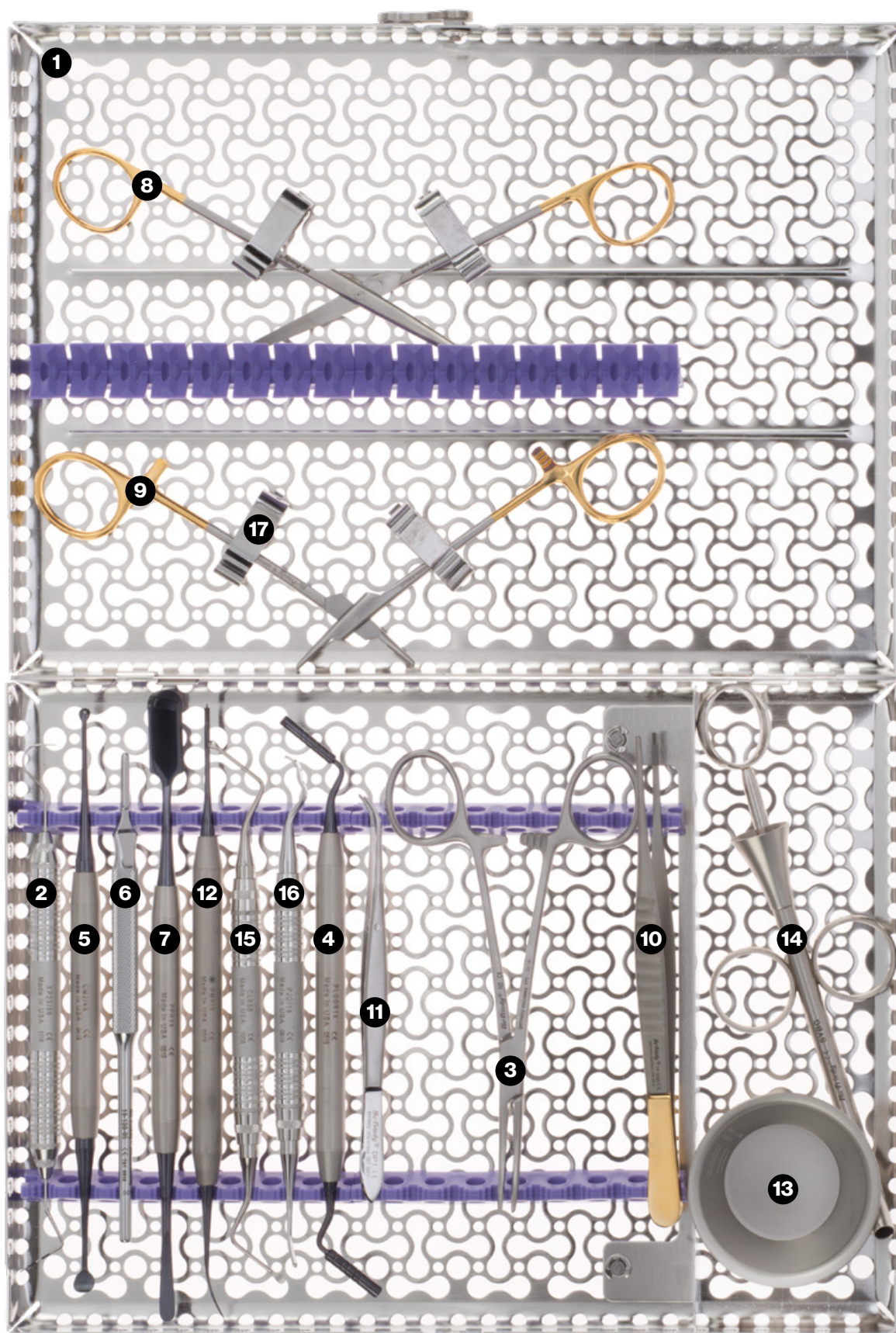

Hu-Friedy Mfg. Co., LLC, 1666 E. Touhy Ave., Des Plaines, IL 60018 | [HuFriedyGroup.com](http://HuFriedyGroup.com)

All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted. Marks not registered in all jurisdictions.

©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. HFL-347/1023