



# HU-FRIEDY PROCEDURAL SET-UP


## RUBBER DAM CLAMP KIT LIST



PART CODE	DESCRIPTION	QUANTITY
<b>RDCMW2</b>	Clamp, Rubber Dam #W2, Satin Steel	
<b>RDCMW3</b>	Clamp, Rubber Dam #W3, Satin Steel	
<b>RDCMW7</b>	Clamp, Rubber Dam #W7, Satin Steel	
<b>RDCM27</b>	Clamp, Rubber Dam #27, Satin Steel	
<b>RDCM212</b>	Rubber Dam Clamp #212, Satin Steel	
<b>RDPA</b>	Rubber Dam Punch-Ainsworth	
<b>RDF</b>	Rubber Dam Forceps	
<b>RDTFDI</b>	Rubber Dam Template, FDI	
<b>IM8128CE2L</b>	MTO 12 IN SIG SERIES,CAS,BLUE,ECR,E2,8LP	
<b>RDAF6</b>	Rubber Dam Frame, Adult	
<b>KIT-RDC2</b>	<b>Complete Rubber Dam Clamp Kit</b>	

## RUBBER DAM CLAMP KIT 2 LIST

- A** RDCMW2
- B** RDCMW3
- C** RDCMW7
- D** RDCM27
- E** RDCM212
- F** RDPA
- G** RDF
- H** RDTFDI
- I** IM8128CE2L
- J** RDAF6

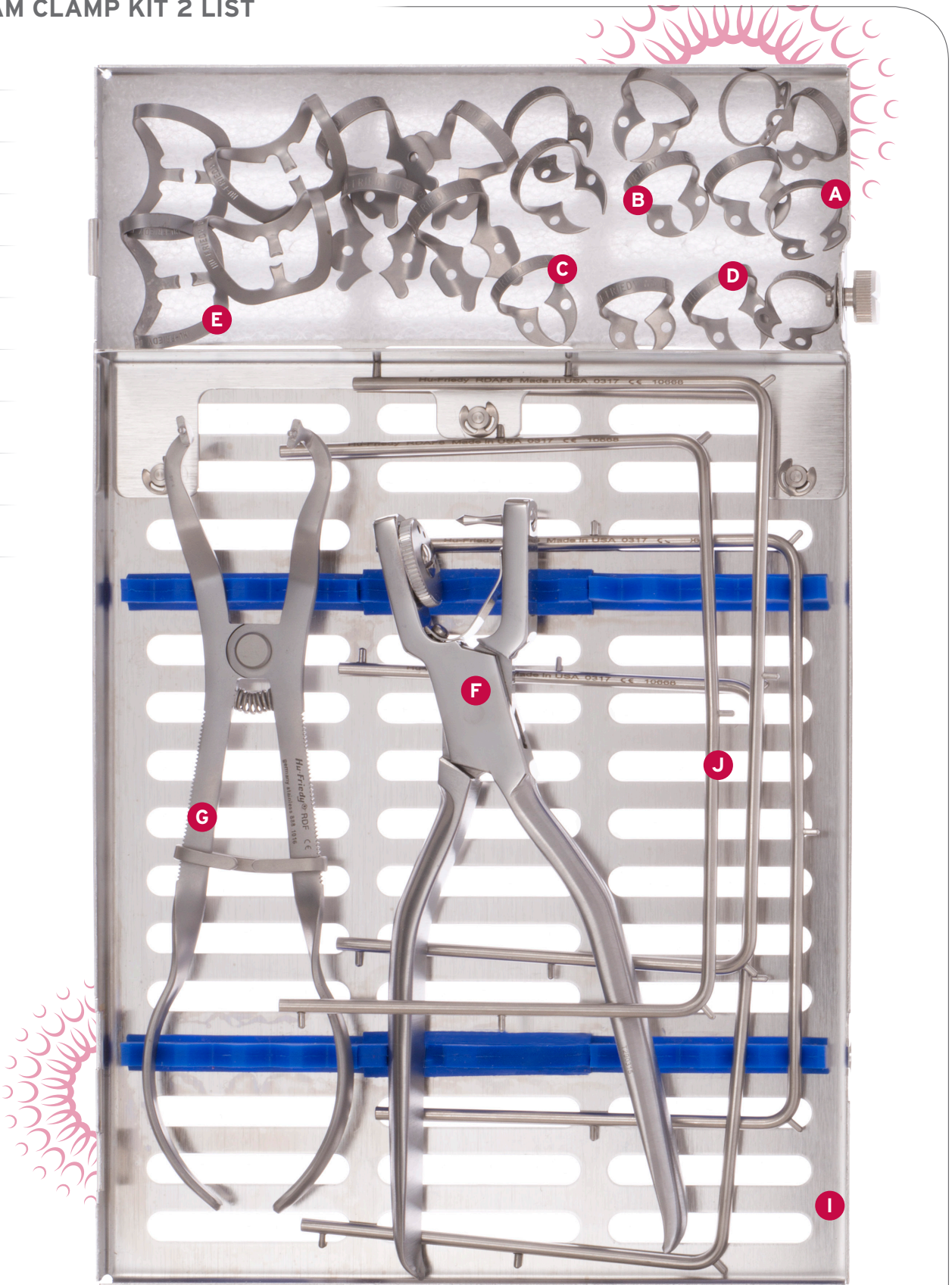

**H** Rubber Dam Template, FDI [not pictured]

## RUBBER DAM CLAMP KIT 2 LIST



- A RDCMW2**  
**Clamp, Rubber Dam #W2, Satin Steel**  
 Flat jaws for premolars.



- B RDCMW3**  
**Clamp, Rubber Dam #W3, Satin Steel**  
 Flat narrow jaws for small lower molars.



- C RDCMW7**  
**Clamp, Rubber Dam #W7, Satin Steel**  
 Flat jaws for lower molars.



- D RDCM27**  
**Clamp, Rubber Dam #27, Satin Steel**  
 Flat jaws for large molars.




- E RDCM212**  
**Rubber Dam Clamp #212, Satin Steel**  
 Flat jaws for labial caries on anterior teeth.



- F RDPA**  
**Rubber Dam Punch-Ainsworth**  
 17 cm (6 2/4) Ainsworth Rubber Dam Punch, essential for efficient treatment. A popular rubber dam clamp accessory, this classic punch is solidly built to for lasting longevity. Hu-Friedy also has a wide variety of rubber dam clamps in winged, wingless and butterfly designs, all made with Satin Steel® to help reduce the glare caused by operator lights for less eye fatigue and improved clinical outcomes.



- G RDF**  
**Rubber Dam Forceps**  
 17 cm (6 1/2") Rubber Dam Clamp Forceps essential for efficient treatment. A popular rubber dam clamp accessory, this instrument features a useful attachment design which makes it easy to pick up and place rubber dam clamps.

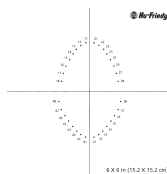

- H RDTFDI**  
**Rubber Dam Template, FDI**  
 A useful way to determine exactly where to punch the rubber dam hole based on tooth location. Simply place the dental dam on the template, find the appropriate tooth location and create your punch hole. Available in Adult and Pedodontic.



- I IM8128CE2L**  
**MTO 12 IN SIG**  
**SERIES, CAS, BLUE, ECR, E2, 8LP**  
 IMS® cassettes are made of long-lasting, low-maintenance stainless steel. Soft, colored silicone rails securely hold instruments throughout the processing cycle. Available in over 100 different sizes, configurations and color options, and can be specialized to meet your needs.



- J RDAF6**  
**Rubber Dam Frame, Adult**  
 One of two stainless steel rubber dam frames, this larger version has a 15cm/6" diameter. Made from high quality stainless steel, its dental dam attachment points provide the durable longevity associated with stainless steel rubber dam frames.